



**TPF**

Fabrikstrasse 8  
Büro A209  
3012 Bern

031 6315513

**GazeToMouse:**  
Usage of the **Tobii Eye Tracker 4C**  
from within a **ztree** program

Tutorial

Simon Maurer - [simon.maurer@humdek.unibe.ch](mailto:simon.maurer@humdek.unibe.ch)  
8th March 2018

# Contents

<b>1</b>	<b>Introduction</b>	<b>3</b>
<b>2</b>	<b>Quick Start</b>	<b>4</b>
2.1	Requesting Access to the Latest ztree Release . . . . .	4
2.2	Requesting a TPF GitLab Account . . . . .	5
<b>3</b>	<b>Use GazeToMouse Toolset in Ztree</b>	<b>6</b>
3.1	Basic Principles of the Tobii Eye Tracker 4C . . . . .	6
3.2	Execute 3rd Party Applications in ztree . . . . .	6
<b>4</b>	<b>Configuration of the GazeToMouse Toolset</b>	<b>8</b>
4.1	Output File . . . . .	10
4.2	Log File . . . . .	11

# Chapter 1

## Introduction

This document describes how to use the **Tobii Eye Tracker 4C** in order to track the gaze of a subject and feed the gaze data to the mouse device. This allows to use the eye tracker to control the mouse pointer position such that the mouse pointer is placed at the screen coordinates where the subject is gazing at. This documentation comes with a set of tools (executable .exe files) that provide this functionality as well as some auxiliary tools that help to calibrate and test the eye tracker.

The project aims at providing a set of executables which allow to use the **Tobii Eye Tracker 4C** in conjunction with **ztree** to perform economic experiments. The following set of executables are provided:

**TobiiCalibrate.exe** This program is a simple wrapper for the Tobii calibration tool. It launches the calibration GUI where the subject is led through the calibration process. The calibration data is stored in the current profile of the eye tracker engine.

**TobiiGuestCalibrate.exe** This program is a simple wrapper for the Tobii guest calibration tool. Similar to **TobiiCalibrate.exe** it launches the calibration tool, however, the calibration data is stored in a guest profile.

**TobiiTest.exe** This program is a simple wrapper for the Tobii eye tracking testing tool. It launches a GUI where the result of the calibration can be verified and a new calibration process can be started if required.

**GazeToMouse.exe** This program uses the **Tobii Core SDK** (default) or the **Tobii Pro SDK** to get the position on the screen where the subject is looking at. The mouse cursor position is updated to this position. As a consequence, the mouse cursor is controlled by the gaze of the subject. This program runs infinitely until it is terminated by an external command. This should **not** be done with a forced kill (e.g. by executing the command `taskkill /F /IM GazeToMouse.exe` or by killing the task with the task manager) because it prevents the program from terminating gracefully. This has several consequences:

- open files are not closed properly and the data stream is cut off. This can lead to corrupt files.
- if the feature of hiding the mouse pointer is used, the mouse will remain hidden.
- memory is not freed properly.

Instead `taskkill /IM GazeToMouse.exe` should be used. This is done in the program **GazeToMouseClose.exe**.

**GazeToMouseClose.exe** This program requests **GazeToMouse.exe** to close gracefully and logs these events to the log file.

**ShowMouse.exe** This program allows to restore the standard mouse pointer. It might be useful if the program **GazeToMouse.exe** crashes or is closed forcefully such that the mouse pointer is not restored after terminating. The subject might end up with a hidden mouse pointer. A good solution for such a case is to install a shortcut to **ShowMouse.exe** on the desktop in order to execute it with the keyboard.

# Chapter 2

## Quick Start

In order to get started with a quick experiment you need the following things (*server* refers to the machine running *ztree* and *client* refers to the machines running *zleaf*):

- a *ztree* installation with a server and one or more clients. To install *ztree* download the latest *ztree* version from the [download section](#) (requires a license and a login, see Section 2.1) and extract the server file *ztree.exe* to the chosen installation path on the server and the client file *zleaf.exe* to the chosen installation path on each client. Throughout this documentation the installation paths of the *ztree* client and server will be referred to as *<zleaf path>* and *<ztree path>*, respectively.
- the *Tobii Eye Tracking Core Software* installed on each client. To install the software [download](#) and execute the installation file. The installation path of this software will be referred to as *<Tobii path>* throughout this documentation.
- the *Tobii Eye Tracker 4C* mounted and connected on each client. When encountering connection problems follow [these](#) instructions which might help to solve the problem.
- the *GazeToMouse toolset* installed on each client. To install the toolset [download](#) the latest version (requires a subscription, see Section 2.2) and extract all files to an installation path of your choosing. The GazeToMouse installation path will be referred to as *<GazeToMouse path>* throughout this document.
- the *ztree* sample file. The file is located in *<GazeToMouse path>\sample\template.ztt*. Make sure to change the *Path* variable in Background of the sample file such that it points to *<GazeToMouse path>*.

The provided sample *ztree* program performs the following stages:

**Calibrate Eye Tracker** A guest calibration is performed where the Tobii calibration tool is started and the calibration data is stored in a guest profile within the Tobii eye tracker software. The calibration is very intuitive and has a game-like feel to it.

**Test Eye Tracker** The Tobii eye tracker test tool is started which allows the subject or experimenter to verify the accuracy of the eye tracker. If the accuracy is not sufficient a recalibration can be started from within the test tool.

**Gaze To Mouse** The gaze of the subject is tracked and transformed to mouse coordinates which allows to control the mouse pointer with the gaze of the subject.

**Terminate** Terminates the gaze-to-mouse functionality and concludes this simple *ztree* program.

The behaviour of the GazeToMouse toolset can be configured with a configuration file (see Chapter 4 for more details). A sample configuration file is provided in *<GazeToMouse path>\sample\config.json* that can be modified to suit different requirements. Copy the configuration file either to *<GazeToMouse path>* or *<zleaf path>* and modify it there.

### 2.1 Requesting Access to the Latest *ztree* Release

To download the latest release of the *ztree* software a login must be requested [here](#). Make a careful note of the terms and conditions.

## 2.2 Requesting a TPF GitLab Account

The necessary code for what is described in this document is contained in a GitLab instance running on the TPF server.

Its address is:

<http://phhum-g111-nns.unibe.ch:10012/>

Please send **Simon Maurer** a mail with a request, specifying your name, the name of the TPF project you are interested in (e.g GazeToMouse) and a user name which you would like to have on GitLab.

## Chapter 3

# Use GazeToMouse Toolset in Ztree

This chapter provides instructions of how to use the Tobii Eye Tracker 4C in a `ztree` program. In order to do this there are two main points to understand:

1. How does the Tobii Eye Tracker 4C work
2. How to execute 3rd party applications from within `ztree`

The first point is addressed in Section 3.1 and the second point in Section 3.2.

### 3.1 Basic Principles of the Tobii Eye Tracker 4C

The eye tracker is mounted on the bottom of the screen such that its cameras are able to capture the eyes of the subject. Several infrared LEDs are visible once the device is connected and properly working (see Figure 3.1). This requires the Tobii software to be installed and running.



Figure 3.1: Picture of the Tobii Eye Tracker 4C in operation.

In order for the eye tracker to be able to track the gaze of the subject it must be calibrated. A user-friendly calibration tool is included in the Tobii software package. For each person different calibration data must be stored. For this reason the Tobii software supports multiple user profiles which can be switched on the fly. Given that in an economic experiment a workstation is used by a subject only once it is suggested to use the *Guest* profile to store the eye tracker calibration data. Once the eye tracker is calibrated for a subject, everything is ready.

### 3.2 Execute 3rd Party Applications in `ztree`

In a first step, it is useful to define a global variable `Path` in a `ztree` program that points to the location of the toolset.

In order to define a path

1. click on the last table in `Background`
2. choose the menu `Treatment → New Program...`
3. define the path of the installation folder of the GazeToMouse toolset, i.e. `<GazeToMouse path>` as `Path` variable (e.g. `Path="C:\\Users\\Max Muster\\Documents\\My Experiment\\GazeToMouse\\";`).  
Note that two backslashes are required for each path delimiter because in `ztree` `\` is used as escape character.

When writing the different stages of a *ztree* program, calls to external applications can be made. This can be achieved as follows:

1. click on the position where you want to include the call to the external program
2. choose the menu *Treatment* → *New External Program...*
3. choose *Run* on *z-Leaf*
4. add the call to the external program to the field *Command Line*  
(e.g a call to the guest calibration tool: `<><Path|-1>TobiiGuestCalibrate.exe`)  
Note that the *Path* variable is used to indicate the location of the program. Figure 3.2 shows an example of calling the Tobii guest calibration tool.

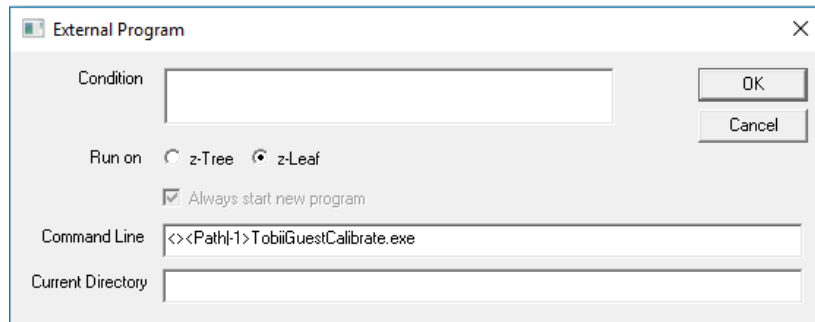


Figure 3.2: External program call definition in *ztree*. The variable *Path* must be globally defined.

## Chapter 4

# Configuration of the GazeToMouse Toolset

The GazeToMouse toolset can be configured to work with different installations. It allows to specify installation paths and gives some control over implemented features (e.g. mouse hiding or gaze filtering). The Listing 4.1 shows the default configuration values and provides a short explanation for each value.

```
1 {
  // Allows to define the order and the delimiters between the different values. The definition is
  // of the form
  // {0}<delim>{1}<delim> ... <delim>{13}
  // where <delim> can be customized (e.g. '\t' for tab, ',' for comma, etc.) and where the numbers are
6  // replaced by the following values
  // - 0: timestamp of the gaze data item (uses ValueFormat.TimeStamp)
  // - 1: x-coordinate of the gaze point of both eyes (pixel value)
  // - 2: x-coordinate of the gaze point of the left eye (pixel value) [SDK Pro only]
  // - 3: x-coordinate of the gaze point of the right eye (pixel value) [SDK Pro only]
11  // - 4: y-coordinate of the gaze point of both eyes (pixel value)
  // - 5: y-coordinate of the gaze point of the left eye (pixel value) [SDK Pro only]
  // - 6: y-coordinate of the gaze point of the right eye (pixel value) [SDK Pro only]
  // - 7: validity of the gaze data of the left eye [SDK Pro only]
  // - 8: validity of the gaze data of the right eye [SDK Pro only]
16  // - 9: average pupil diameter of both eyes (uses ValueFormat.Diameter) [SDK Pro only]
  // - 10: pupil diameter of the left eye (uses ValueFormat.Diameter) [SDK Pro only]
  // - 11: pupil diameter of the right eye (uses ValueFormat.Diameter) [SDK Pro only]
  // - 12: validity of the pupil data of the left eye [SDK Pro only]
  // - 13: validity of the pupul data of the right eye [SDK Pro only]
21  // To log all possible values use:
  // "DataLogColumnOrder": "{0}\t{1}\t{2}\t{3}\t{4}\t{5}\t{6}\t{7}\t{8}\t{9}\t{10}\t{11}\t{12}\t{13}",
  "DataLogColumnOrder": "{0}\t{1}\t{4}",

  // Defines the titles of the data log value columns. A title for all possible columns must be defined.
26  // Titles for values that are removed from the "DataLogColumnOrder" parameter will not be logged but
  // must still be defined here. The index of a title must correspond to the value number of the
  // configuration parameter "DataLogColumnOrder".
  "DataLogColumnTitle": [
    "Timestamp",
31    "coord-x",
    "coord-x-left",
    "coord-x-right",
    "coord-y",
    "coord-y-left",
36    "coord-y-right",
    "coord-valid-left",
    "coord-valid-right",
    "pupil-dia",
    "pupil-dia-left",
41    "pupil-dia-right",
    "pupil-valid-left",
    "pupil-valid-right"
  ],
}
```



```

46 // Number of output data files to allowed the output path. Oldest files are deleted first. When
// setting OutputCount to 0, no files are removed. A value of 1 means that only the output of
// current execution is kept.
"DataLogCount": 200,

51 // Allows to define the format of how the pupil diameter will be logged. Use the .NET syntax to
// specify the format:
// https://docs.microsoft.com/en-us/dotnet/standard/base-types/formatting-types
// Note that the numbers will be represented differently depending in the localisation settings of the
// windows installation (e.g. 123,4 for DE_CH or 123.4 for EN_US).
56 "DataLogFormatDiameter": "0.000",

// Allows to define the format of the timestamp. Use the .NET syntax to specify the format:
// https://docs.microsoft.com/en-us/dotnet/standard/base-types/formatting-types
// Note that special characters (e.g. ':', '.') need to be escaped with '\\'.
61 "DataLogFormatTimeStamp": "hh\\:mm\\:ss\\.fff",

// Defines the location of the output file. It must be the path to a folder (not a file). If empty,
// the output file is produced in the directory of the caller (e.g the directory of zleaf.exe).
"DataLogPath": "",

66 // Defines whether gaze data is written to a log file. If set to false, the configuration items
// "OutputCount", "OutputPath", and "OutputFormat" are ignored.
"DataLogWriteOutput": true,

71 // Defines filter settings for the eye tracker:
// - 0: unfiltered, Tobii: "removal of invalid data points and the averaging of the gaze points from
// both eyes".
// - 1: lightly filtered, Tobii: "an adaptive filter which is weighted based on the age of the gaze
// data points GazePointData and the velocity of the eye movements".
76 // also see https://tobii.github.io/CoreSDK/articles/streams.html
// Note that this parameter only has effect for Tobii SDK Core (see parameter "TobiiSDK").
"GazeFilterCore": 0,

// Defines the location of the license files. It must be the path to a folder (not a file).
81 // This is only required when Tobii SDK Pro is used (see parameter "TobiiSDK")
"LicensePath": "licenses",

// Defines wheter the mouse cursor shall be controlled by the gaze of the subject during the
// experiment. If set to true the mouse cursor will be controlled by the gaze of the subject when
86 // GazeToMouse.exe is executed and control will be released when GazeToMouseClose.exe is executed.
// If set to false the configuration item "HideMouse" is ignored.
"MouseControl": true,

// Defines whether the mouse cursor shall be hidden during the experiemnt. If set to true the
91 // mouse cursor will be hidden when GazeToMouse.exe is executed and restored when
// GazeToMousClose.exe is executed.
"MouseHide": false,

// Defines the Path to the standard mouse pointer icon. This is used to restore the mouse pointer.
96 "MouseStandardIconPath": "C:\\Windows\\Cursors\\aero_arrow.cur",

// Specifies the amount of time (in milliseconds) to wait for the eyetracker to become ready. If the
// eyetracker is not ready withing the specified time the subject will be notified with a popup window.
"ReadyTimer": 5000,

101 // Choose the Tobii SDK (0: Tobii SDK Core, 1: Tobii SDK Pro).
// Note that Tobii SDK Pro requires a license file to work (see parameter "LicesePath").
"TobiiSDK": 0,

106 // Defines the Tobii installation path. It must be the path to a folder (not a file).
"TobiiApplicationPath": "C:\\Program Files (x86)\\Tobii",

// The Tobii EyeX command to run a calibration.
"TobiiCalibrate": "Tobii EyeX Config\\Tobii.EyeX.Configuration.exe",

```

```

111 // The Tobii EyeX parameter to run a calibration.
    "TobiiCalibrateArguments": "--calibrate",

116 // The Tobii EyeX command to run a guest calibration.
    "TobiiGuestCalibrate": "Tobii EyeX Config\\Tobii.EyeX.Configuration.exe",

    // The Tobii EyeX parameter to run a guest calibration.
    "TobiiGuestCalibrateArguments": "--guest-calibration",

121 // The Tobii EyeX command to run calibration test.
    "TobiiTest": "Tobii EyeX Interaction\\Tobii.EyeX.Interaction.TestEyeTracking.exe"
}

```

Listing 4.1: Default configuration values

Note that the configuration file follows the json syntax which must not be violated. If the following points are respected, no problem should arise:

- the configuration values are enclosed in '{' and '}'.
- each configuration line ends with a ',' except for the last line where it is omitted.
- the Windows path delimiter '\' must be escaped (i.e. write '\\' when describing a path)
- json supports standard data types (e.g. integer, boolean, string). Use the same type as the default value.
- everything following a '/' is considered a comment.

Each executable of the toolset uses the same common configuration file. The configuration file must be named `config.json` and is read from the following places with the indicated priority:

1. in the directory of the caller, i.e in `<zleaf path>`
2. in the directory of the executables, i.e. in `<GazeToMouse path>`

If no configuration file can be found, the default values are used.

**Warning:** A word of warning when using the mouse hiding feature. This feature hides the mouse when running `GazeToMouse.exe`. If `GazeToMouse.exe` is forcefully closed or crashes, the mouse pointer stays hidden. For such cases the `ShowMouse.exe` utility can be used.

## 4.1 Output File

When running the program `gazeToMouse.exe` an output file can be generated which holds the gaze data provided by the Tobii engine. The output file is saved in the directory specified by `OutputPath` in `config.json`. The name of the output file follows the form `<yyyyMMddTHHmss>_<hostName>_gaze.txt` where

- `<yyyyMMddTHHmss>` is replaced by the timestamp indicating when the file was created (e.g. 20180129T085521 stands for 29.01.2018 08:55:21).
- `<hostName>` is replaced by the name of the machine

An output file is only generated if the parameter `"DataLogWriteOutput"` in the configuration file is set to `true` (which is the default). The presentation of the gaze data is configurable with the help of the parameters `"DataLogColumnOrder"`, `"DataLogColumnNameTitle"`, `"DataLogFormatDiameter"`, and `"DataLogFormatTimeStamp"` (see Listing 4.1 for more details).

By default the **Tobii Core SDK** is used (parameter `"TobiiSDK"` is set to 0) and only the timestamps, the x and the y coordinates of the gaze point are available. Hence, the parameter `"DataLogColumnOrder"` is set to the value `"{0}\t{1}\t{4}"` in order to not bloat the output file with unnecessary data.

However, when using **Tobii Pro SDK**, much more values are provided and can be logged to the output file. This can be achieved by setting the value of parameter `"DataLogColumnOrder"` as follows:

```
"DataLogColumnOrder": "{0}\t{1}\t{2}\t{3}\t{4}\t{5}\t{6}\t{7}\t{8}\t{9}\t{10}\t{11}\t{12}\t{13}"
```

Note that when using **Tobii Pro SDK**, the values 1 (x coordinate of the gaze point), 4 (y coordinate of the gaze point), and 9 (pupil diameter) are computed by averaging the values of both eyes:

$$\text{val} = \begin{cases} \text{NaN} & \text{if val\_right and val\_left are not valid} \\ \text{val\_left} & \text{if val\_right is not valid} \\ \text{val\_right} & \text{if val\_left is not valid} \\ \frac{\text{val\_left} + \text{val\_right}}{2} & \text{otherwise} \end{cases}$$

When using **Tobii Core SDK**, filter settings can be configured with the parameter "GazeFilterCore" (see Listing 4.1 for more details).

## 4.2 Log File

All executables write continuously to the same log file. This allows to track the eye tracker events that happened throughout a `ztree` session within one log file. The log file is produced at the root directory of the application which is making the calls to the executables (e.g. at the location of `zleaf.exe`: `<zleaf path>`). The name of the log file is of the form `<hostName>_gaze.log` where `<hostName>` is replaced by the name of the machine.