

Gaze Toolset

v3.4.2

Generated by Doxygen 1.8.17

1 Changelog	1
2 Toolset to Control Tobii Eye Tracker	9
3 Sample Files for Experimentation with Eye Tracker Utility	13
4 Namespace Index	15
4.1 Namespace List	15
5 Hierarchical Index	17
5.1 Class Hierarchy	17
6 Class Index	19
6.1 Class List	19
7 Namespace Documentation	23
7.1 CustomCalibrationLibrary Namespace Reference	23
7.2 CustomCalibrationLibrary.Commands Namespace Reference	23
7.3 CustomCalibrationLibrary.Converters Namespace Reference	23
7.4 CustomCalibrationLibrary.Models Namespace Reference	23
7.4.1 Enumeration Type Documentation	24
7.4.1.1 CalibrationEventType	24
7.4.1.2 CalibrationStatus	24
7.5 CustomCalibrationLibrary.ViewModels Namespace Reference	24
7.6 CustomCalibrationLibrary.Views Namespace Reference	25
7.7 GazeControl Namespace Reference	25
7.8 GazeControlLibrary Namespace Reference	25
7.8.1 Enumeration Type Documentation	26
7.8.1.1 LogLevel	26
7.9 GazeToMouse Namespace Reference	26
7.10 GazeUtilityLibrary Namespace Reference	26
7.10.1 Detailed Description	27
7.10.2 Enumeration Type Documentation	27
7.10.2.1 ECalibrationDataError	27
7.10.2.2 EGazeConfigError	28
7.10.2.3 EGazeDataError	28
7.10.2.4 EOutputType	28
7.11 GazeUtilityLibrary.DataStructs Namespace Reference	28
7.11.1 Enumeration Type Documentation	29
7.11.1.1 CalibrationOutputValue	29
7.11.1.2 GazeOutputValue	29
7.11.1.3 ValidationOutputValue	30
7.12 GazeUtilityLibrary.Tracker Namespace Reference	30
7.13 ShowMouse Namespace Reference	30

7.14 Tobii Namespace Reference	30
7.15 Tobii.Research Namespace Reference	30
7.16 Tobii.Research.Addons Namespace Reference	30
7.17 Tobii.Research.Addons.Utility Namespace Reference	31
7.18 TobiiCalibrate Namespace Reference	31
8 Class Documentation	33
8.1 GazeControl.App Class Reference	33
8.1.1 Detailed Description	34
8.2 GazeToMouse.App Class Reference	34
8.2.1 Detailed Description	35
8.2.2 Constructor & Destructor Documentation	35
8.2.2.1 App()	35
8.2.3 Member Function Documentation	36
8.2.3.1 CalibrationValidate()	36
8.2.3.2 CompensateDrift()	36
8.2.3.3 CustomCalibrate()	36
8.2.3.4 GazeRecordingDisable()	36
8.2.3.5 GazeRecordingEnable()	37
8.2.3.6 MouseTrackingDisable()	37
8.2.3.7 MouseTrackingEnable()	37
8.2.3.8 ResetDriftCompensation()	37
8.2.4 Property Documentation	37
8.2.4.1 LastTag	37
8.2.4.2 StartTime	37
8.2.4.3 Tag	38
8.2.4.4 TrialId	38
8.3 ShowMouse.App Class Reference	38
8.3.1 Detailed Description	39
8.4 TobiiCalibrate.App Class Reference	39
8.4.1 Detailed Description	39
8.5 GazeUtilityLibrary.Tracker.BaseTracker Class Reference	40
8.5.1 Detailed Description	43
8.5.2 Member Enumeration Documentation	43
8.5.2.1 DeviceStatus	43
8.5.3 Constructor & Destructor Documentation	43
8.5.3.1 BaseTracker()	43
8.5.4 Member Function Documentation	44
8.5.4.1 ApplyCalibration()	44
8.5.4.2 CollectCalibrationDataAsync()	44
8.5.4.3 CollectValidationDataAsync()	45
8.5.4.4 ComputeValidation()	46

8.5.4.5 Dispose() [1/2]	46
8.5.4.6 Dispose() [2/2]	46
8.5.4.7 DriftCompensationEventHandler()	47
8.5.4.8 FinishCalibration()	47
8.5.4.9 FinishCalibrationAsync()	47
8.5.4.10 FinishValidation()	47
8.5.4.11 GazeDataHandler()	47
8.5.4.12 GetFixationFrameCount()	48
8.5.4.13 GetUnitDirection()	48
8.5.4.14 InitCalibration()	48
8.5.4.15 InitCalibrationAsync()	49
8.5.4.16 InitDriftCompensation()	49
8.5.4.17 InitValidation()	49
8.5.4.18 IsInitialised()	49
8.5.4.19 IsReady()	50
8.5.4.20 OnGazeDataReceived()	50
8.5.4.21 OnPropertyChanged()	50
8.5.4.22 OnTrackerDisabled()	50
8.5.4.23 OnTrackerDisabledTimeout()	51
8.5.4.24 OnTrackerEnabled()	51
8.5.4.25 OnUserPositionDataReceived()	51
8.5.4.26 PatternReplace()	51
8.5.4.27 ResetDriftCompensation()	52
8.5.4.28 StartDriftCompensation()	52
8.5.4.29 UserPositionDataHandler()	52
8.5.5 Member Data Documentation	52
8.5.5.1 config	52
8.5.5.2 DeviceName	53
8.5.5.3 dialogBoxTimer	53
8.5.5.4 driftCompensation	53
8.5.5.5 logger	53
8.5.5.6 screenArea	53
8.5.5.7 trackerMessageBox	53
8.5.6 Property Documentation	54
8.5.6.1 ScreenArea	54
8.5.6.2 State	54
8.5.7 Event Documentation	54
8.5.7.1 DriftCompensationComputed	54
8.5.7.2 GazeDataReceived	54
8.5.7.3 PropertyChanged	54
8.5.7.4 TrackerDisabled	55
8.5.7.5 TrackerEnabled	55

8.5.7.6 UserPositionDataReceived	55
8.6 CustomCalibrationLibrary.Views.Calibration Class Reference	55
8.6.1 Detailed Description	56
8.6.2 Constructor & Destructor Documentation	56
8.6.2.1 Calibration()	56
8.7 CustomCalibrationLibrary.Commands.CalibrationCommand Class Reference	57
8.7.1 Detailed Description	58
8.7.2 Constructor & Destructor Documentation	58
8.7.2.1 CalibrationCommand()	58
8.7.3 Member Function Documentation	58
8.7.3.1 CanExecute()	58
8.7.3.2 Execute()	58
8.7.4 Property Documentation	59
8.7.4.1 CanExecuteChanged	59
8.8 GazeUtilityLibrary.CalibrationDataError Class Reference	59
8.8.1 Detailed Description	60
8.8.2 Member Function Documentation	60
8.8.2.1 GetCalibrationDataErrorString()	60
8.8.3 Property Documentation	61
8.8.3.1 Error	61
8.9 CustomCalibrationLibrary.Views.CalibrationFailed Class Reference	61
8.9.1 Detailed Description	62
8.9.2 Constructor & Destructor Documentation	62
8.9.2.1 CalibrationFailed()	62
8.9.3 Property Documentation	62
8.9.3.1 CalibrationAbortCommand	63
8.9.3.2 CalibrationRestartCommand	63
8.9.3.3 Error	63
8.9.4 Event Documentation	63
8.9.4.1 PropertyChanged	63
8.10 CustomCalibrationLibrary.Views.CalibrationFrame Class Reference	64
8.10.1 Detailed Description	64
8.10.2 Constructor & Destructor Documentation	64
8.10.2.1 CalibrationFrame()	65
8.11 CustomCalibrationLibrary.Models.CalibrationModel Class Reference	65
8.11.1 Detailed Description	67
8.11.2 Constructor & Destructor Documentation	67
8.11.2.1 CalibrationModel()	67
8.11.3 Member Function Documentation	67
8.11.3.1 GazeDataCollected()	68
8.11.3.2 InitCalibration()	68
8.11.3.3 NextCalibrationPoint()	68

8.11.3.4 OnCalibrationEvent()	68
8.11.3.5 RedoCalibrationPoint()	68
8.11.3.6 SetCalibrationResult()	68
8.11.3.7 UpdateGazePoint()	69
8.11.4 Property Documentation	69
8.11.4.1 CalibrationPoints	69
8.11.4.2 CursorType	69
8.11.4.3 Error	69
8.11.4.4 GazePoint	70
8.11.4.5 Index	70
8.11.4.6 LastStatus	70
8.11.4.7 Points	70
8.11.4.8 Status	70
8.11.4.9 UserPositionGuide	70
8.11.4.10 ValidationData	71
8.11.5 Event Documentation	71
8.11.5.1 CalibrationEvent	71
8.11.5.2 GazePointChanged	71
8.11.5.3 PropertyChanged	71
8.11.5.4 UserPositionGuideChanged	71
8.12 GazeUtilityLibrary.DataStructs.CalibrationPoint Class Reference	72
8.12.1 Detailed Description	73
8.12.2 Constructor & Destructor Documentation	73
8.12.2.1 CalibrationPoint()	73
8.12.3 Property Documentation	73
8.12.3.1 GazePositionAverage	73
8.12.3.2 GazePositionLeft	74
8.12.3.3 GazePositionRight	74
8.12.3.4 HasData	74
8.12.3.5 Index	74
8.12.3.6 Position	74
8.12.4 Event Documentation	74
8.12.4.1 PropertyChanged	74
8.13 CustomCalibrationLibrary.Views.CalibrationPoint Class Reference	75
8.13.1 Detailed Description	75
8.13.2 Constructor & Destructor Documentation	75
8.13.2.1 CalibrationPoint()	76
8.14 CustomCalibrationLibrary.ViewModels.CalibrationPointViewModel Class Reference	76
8.14.1 Detailed Description	77
8.14.2 Constructor & Destructor Documentation	77
8.14.2.1 CalibrationPointViewModel() [1/2]	77
8.14.2.2 CalibrationPointViewModel() [2/2]	77

8.15 CustomCalibrationLibrary.Views.CalibrationResult Class Reference	78
8.15.1 Detailed Description	78
8.15.2 Constructor & Destructor Documentation	79
8.15.2.1 CalibrationResult()	79
8.16 CustomCalibrationLibrary.Views.CalibrationResultPoint Class Reference	79
8.16.1 Detailed Description	80
8.16.2 Constructor & Destructor Documentation	80
8.16.2.1 CalibrationResultPoint()	80
8.17 CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel Class Reference	81
8.17.1 Detailed Description	82
8.17.2 Constructor & Destructor Documentation	82
8.17.2.1 CalibrationResultViewModel()	82
8.17.3 Member Function Documentation	82
8.17.3.1 OnGazeToggle()	82
8.17.4 Property Documentation	83
8.17.4.1 CalibrationAcceptCommand	83
8.17.4.2 CalibrationRestartCommand	83
8.17.4.3 GazePoint	83
8.17.4.4 GazeVisibilityCommand	83
8.18 Tobii.Research.Addons.CalibrationValidationPoint Class Reference	83
8.18.1 Detailed Description	84
8.18.2 Member Function Documentation	84
8.18.2.1 ToString()	84
8.18.3 Property Documentation	84
8.18.3.1 AccuracyLeftEye	85
8.18.3.2 AccuracyRightEye	85
8.18.3.3 Coordinates	85
8.18.3.4 GazeData	85
8.18.3.5 PrecisionLeftEye	85
8.18.3.6 PrecisionRightEye	85
8.18.3.7 PrecisionRMSLeftEye	86
8.18.3.8 PrecisionRMSRightEye	86
8.18.3.9 TimedOut	86
8.19 Tobii.Research.Addons.CalibrationValidationResult Class Reference	86
8.19.1 Detailed Description	87
8.19.2 Member Function Documentation	87
8.19.2.1 ToString()	87
8.19.3 Property Documentation	87
8.19.3.1 AverageAccuracyLeftEye	87
8.19.3.2 AverageAccuracyRightEye	87
8.19.3.3 AveragePrecisionLeftEye	87
8.19.3.4 AveragePrecisionRightEye	88

8.19.3.5 AveragePrecisionRMSLeftEye	88
8.19.3.6 AveragePrecisionRMSRightEye	88
8.19.3.7 Points	88
8.20 CustomCalibrationLibrary.ViewModels.CalibrationViewModel Class Reference	89
8.20.1 Detailed Description	90
8.20.2 Constructor & Destructor Documentation	90
8.20.2.1 CalibrationViewModel()	90
8.20.3 Member Data Documentation	90
8.20.3.1 _model	90
8.20.4 Property Documentation	90
8.20.4.1 CalibrationPoints	91
8.21 CustomCalibrationLibrary.Views.CalibrationWindow Class Reference	91
8.21.1 Detailed Description	92
8.22 CustomCalibrationLibrary.Views.Computing Class Reference	92
8.22.1 Detailed Description	93
8.22.2 Constructor & Destructor Documentation	93
8.22.2.1 Computing()	93
8.23 GazeUtilityLibrary.ConfigItem Class Reference	93
8.23.1 Detailed Description	95
8.23.2 Constructor & Destructor Documentation	95
8.23.2.1 ConfigItem()	95
8.23.3 Property Documentation	95
8.23.3.1 CalibrationLogColumnOrder	95
8.23.3.2 CalibrationLogColumnTitle	95
8.23.3.3 CalibrationLogWriteOutput	96
8.23.3.4 CalibrationPoints	96
8.23.3.5 ConfigName	96
8.23.3.6 DataLogColumnOrder	96
8.23.3.7 DataLogColumnTitle	96
8.23.3.8 DataLogCount	96
8.23.3.9 DataLogDisabledOnStartup	97
8.23.3.10 DataLogFormatDiameter	97
8.23.3.11 DataLogFormatNormalizedPoint	97
8.23.3.12 DataLogFormatOrigin	97
8.23.3.13 DataLogFormatTimeStamp	97
8.23.3.14 DataLogFormatTimeStampRelative	97
8.23.3.15 DataLogFormatValidation	98
8.23.3.16 DataLogPath	98
8.23.3.17 DataLogWriteOutput	98
8.23.3.18 DriftCompensationDispersionThreshold	98
8.23.3.19 DriftCompensationDispersionThresholdMax	98
8.23.3.20 DriftCompensationDurationThreshold	98

8.23.3.21 DriftCompensationTimer	99
8.23.3.22 DriftCompensationWindowShow	99
8.23.3.23 LicensePath	99
8.23.3.24 MouseCalibrationHide	99
8.23.3.25 MouseControl	99
8.23.3.26 MouseControlHide	99
8.23.3.27 MouseStandardIconPath	100
8.23.3.28 ReadyTimer	100
8.23.3.29 ScreenArea	100
8.23.3.30 TobiiApplicationPath	100
8.23.3.31 TobiiCalibrate	100
8.23.3.32 TobiiCalibrateArguments	100
8.23.3.33 TrackerDevice	101
8.23.3.34 ValidationDurationThreshold	101
8.23.3.35 ValidationLogColumnOrder	101
8.23.3.36 ValidationLogColumnName	101
8.23.3.37 ValidationLogWriteOutput	101
8.23.3.38 ValidationPoints	101
8.23.3.39 ValidationTimer	102
8.24 GazeUtilityLibrary.ConfigScreenArea Class Reference	102
8.24.1 Detailed Description	102
8.24.2 Constructor & Destructor Documentation	102
8.24.2.1 ConfigScreenArea() [1/2]	103
8.24.2.2 ConfigScreenArea() [2/2]	103
8.24.3 Property Documentation	103
8.24.3.1 BottomLeft	103
8.24.3.2 BottomRight	103
8.24.3.3 Center	103
8.24.3.4 Height	104
8.24.3.5 TopLeft	104
8.24.3.6 TopRight	104
8.24.3.7 Width	104
8.25 CustomCalibrationLibrary.Views.Disconnect Class Reference	105
8.25.1 Detailed Description	105
8.25.2 Constructor & Destructor Documentation	106
8.25.2.1 Disconnect()	106
8.25.3 Property Documentation	106
8.25.3.1 CalibrationAbortCommand	106
8.26 GazeUtilityLibrary.DriftCompensation Class Reference	106
8.26.1 Detailed Description	107
8.26.2 Constructor & Destructor Documentation	107
8.26.2.1 DriftCompensation()	107

8.26.3 Member Function Documentation	107
8.26.3.1 Reset()	107
8.26.3.2 Start()	108
8.26.3.3 Update()	108
8.26.4 Property Documentation	108
8.26.4.1 DeviationAngle	108
8.26.4.2 Dispersion	108
8.26.4.3 Q	109
8.27 GazeUtilityLibrary.DataStructs.DriftCompensationData Class Reference	109
8.27.1 Detailed Description	109
8.27.2 Constructor & Destructor Documentation	109
8.27.2.1 DriftCompensationData()	109
8.27.3 Property Documentation	110
8.27.3.1 Compensation	110
8.27.3.2 GazePosition2d	110
8.27.3.3 GazePosition3d	110
8.28 CustomCalibrationLibrary.ViewModels.DriftCompensationViewModel Class Reference	110
8.28.1 Detailed Description	111
8.28.2 Constructor & Destructor Documentation	111
8.28.2.1 DriftCompensationViewModel()	111
8.28.3 Property Documentation	111
8.28.3.1 FixationPoint	111
8.29 CustomCalibrationLibrary.Views.DriftCompensationWindow Class Reference	112
8.29.1 Detailed Description	112
8.29.2 Constructor & Destructor Documentation	112
8.29.2.1 DriftCompensationWindow()	113
8.30 GazeUtilityLibrary.DataStructs.EyeData Class Reference	113
8.30.1 Detailed Description	113
8.30.2 Constructor & Destructor Documentation	113
8.30.2.1 EyeData()	113
8.30.3 Property Documentation	114
8.30.3.1 IsPupilDiameterValid	114
8.30.3.2 PupilDiameter	114
8.31 GazeUtilityLibrary.Tracker.EyeTrackerPro Class Reference	114
8.31.1 Detailed Description	116
8.31.2 Constructor & Destructor Documentation	116
8.31.2.1 EyeTrackerPro()	116
8.31.3 Member Function Documentation	116
8.31.3.1 ApplyCalibration()	116
8.31.3.2 CollectCalibrationDataAsync()	117
8.31.3.3 CollectValidationDataAsync()	117
8.31.3.4 ComputeValidation()	118

8.31.3.5 FinishCalibration()	118
8.31.3.6 FinishCalibrationAsync()	118
8.31.3.7 FinishValidation()	118
8.31.3.8 GetFixationFrameCount()	118
8.31.3.9 GetUnitDirection()	119
8.31.3.10 InitCalibration()	119
8.31.3.11 InitCalibrationAsync()	119
8.31.3.12 InitDriftCompensation()	120
8.31.3.13 InitValidation()	120
8.31.3.14 IsInitialised()	120
8.31.3.15 IsLicenseOk()	120
8.31.3.16 PatternReplace()	121
8.32 CustomCalibrationLibrary.Views.FixationPoint Class Reference	121
8.32.1 Detailed Description	122
8.32.2 Constructor & Destructor Documentation	122
8.32.2.1 FixationPoint()	122
8.33 GazeUtilityLibrary.DataStructs.GazeCalibrationData Class Reference	122
8.33.1 Detailed Description	123
8.33.2 Constructor & Destructor Documentation	123
8.33.2.1 GazeCalibrationData()	123
8.33.3 Member Function Documentation	124
8.33.3.1 Prepare()	124
8.33.4 Property Documentation	124
8.33.4.1 ValidityLeft	124
8.33.4.2 ValidityRight	124
8.33.4.3 XCoord	124
8.33.4.4 XCoordLeft	125
8.33.4.5 XCoordRight	125
8.33.4.6 YCoord	125
8.33.4.7 YCoordLeft	125
8.33.4.8 YCoordRight	125
8.34 GazeUtilityLibrary.GazeConfigError Class Reference	126
8.34.1 Detailed Description	127
8.34.2 Member Function Documentation	127
8.34.2.1 GetGazeConfigErrorString()	127
8.34.3 Property Documentation	127
8.34.3.1 Error	127
8.35 GazeUtilityLibrary.GazeConfiguration Class Reference	127
8.35.1 Detailed Description	128
8.35.2 Constructor & Destructor Documentation	128
8.35.2.1 GazeConfiguration()	128
8.35.3 Member Function Documentation	129

8.35.3.1 CleanupCalibrationOutputFile()	129
8.35.3.2 CleanupGazeOutputFile()	129
8.35.3.3 CleanupValidationOutputFile()	129
8.35.3.4 DumpCurrentConfigurationFile()	130
8.35.3.5 InitConfig()	130
8.35.3.6 PrepareCalibrationOutputFile()	130
8.35.3.7 PrepareGazeOutputFile()	131
8.35.3.8 PrepareValidationOutputFile()	131
8.35.3.9 WriteToCalibrationOutput()	131
8.35.3.10 WriteToGazeOutput()	132
8.35.3.11 WriteToValidationOutput()	132
8.35.4 Property Documentation	132
8.35.4.1 Config	132
8.36 GazeUtilityLibrary.DataStructs.GazeData Class Reference	133
8.36.1 Detailed Description	133
8.36.2 Constructor & Destructor Documentation	133
8.36.2.1 GazeData() [1/3]	133
8.36.2.2 GazeData() [2/3]	134
8.36.2.3 GazeData() [3/3]	134
8.36.3 Member Function Documentation	135
8.36.3.1 Prepare()	135
8.36.4 Property Documentation	136
8.36.4.1 Combined	136
8.36.4.2 DriftCompensation	136
8.36.4.3 Left	136
8.36.4.4 Right	136
8.36.4.5 Timestamp	136
8.36.4.6 TimestampReceived	137
8.37 GazeUtilityLibrary.DataStructs.GazeData2d Class Reference	137
8.37.1 Detailed Description	137
8.37.2 Constructor & Destructor Documentation	137
8.37.2.1 GazeData2d()	137
8.37.3 Property Documentation	138
8.37.3.1 GazePoint	138
8.37.3.2 IsGazePointValid	138
8.38 GazeUtilityLibrary.DataStructs.GazeData3d Class Reference	138
8.38.1 Detailed Description	139
8.38.2 Constructor & Destructor Documentation	139
8.38.2.1 GazeData3d()	139
8.38.3 Property Documentation	139
8.38.3.1 GazeDirection	139
8.38.3.2 GazeDistance	139

8.38.3.3 GazeOrigin	140
8.38.3.4 GazePoint	140
8.38.3.5 IsGazeOriginValid	140
8.38.3.6 IsGazePointValid	140
8.39 GazeUtilityLibrary.DataStructs.GazeDataCollection Class Reference	140
8.39.1 Detailed Description	141
8.39.2 Constructor & Destructor Documentation	141
8.39.2.1 GazeDataCollection() [1/2]	141
8.39.2.2 GazeDataCollection() [2/2]	141
8.39.3 Property Documentation	142
8.39.3.1 EyeData	142
8.39.3.2 GazeData2d	142
8.39.3.3 GazeData3d	142
8.40 GazeUtilityLibrary.GazeDataError Class Reference	143
8.40.1 Detailed Description	144
8.40.2 Member Function Documentation	144
8.40.2.1 GetGazeDataErrorString()	144
8.40.3 Property Documentation	144
8.40.3.1 Error	144
8.41 GazeUtilityLibrary.GazeError Class Reference	144
8.41.1 Detailed Description	145
8.41.2 Member Function Documentation	145
8.41.2.1 ConvertToBinString()	145
8.42 GazeUtilityLibrary.DataStructs.GazeValidationData Class Reference	145
8.42.1 Detailed Description	146
8.42.2 Constructor & Destructor Documentation	146
8.42.2.1 GazeValidationData() [1/2]	146
8.42.2.2 GazeValidationData() [2/2]	146
8.42.3 Member Function Documentation	147
8.42.3.1 AddPoint()	147
8.42.4 Property Documentation	147
8.42.4.1 AccuracyLeft	148
8.42.4.2 AccuracyRight	148
8.42.4.3 Points	148
8.42.4.4 PrecisionLeft	148
8.42.4.5 PrecisionRight	148
8.42.4.6 PrecisionRmsLeft	148
8.42.4.7 PrecisionRmsRight	149
8.43 GazeUtilityLibrary.DataStructs.GazeValidationPoint Class Reference	149
8.43.1 Detailed Description	149
8.43.2 Constructor & Destructor Documentation	149
8.43.2.1 GazeValidationPoint()	149

8.43.3 Member Function Documentation	150
8.43.3.1 Prepare()	150
8.43.4 Property Documentation	150
8.43.4.1 Point	150
8.43.4.2 Result	150
8.44 CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter Class Reference	151
8.44.1 Detailed Description	151
8.44.2 Member Function Documentation	152
8.44.2.1 Convert()	152
8.44.2.2 ConvertBack()	152
8.45 GazeUtilityLibrary.JsonConfigParser Class Reference	153
8.45.1 Detailed Description	153
8.45.2 Constructor & Destructor Documentation	153
8.45.2.1 JsonConfigParser()	153
8.45.3 Member Function Documentation	154
8.45.3.1 GetDefaultConfig()	154
8.45.3.2 ParseJsonConfig()	154
8.45.3.3 SerializeJsonConfig()	154
8.46 GazeUtilityLibrary.DataStructs.LiveGazePoint Class Reference	155
8.46.1 Detailed Description	156
8.46.2 Property Documentation	156
8.46.2.1 Visibility	156
8.46.2.2 X	156
8.46.2.3 Y	156
8.46.3 Event Documentation	156
8.46.3.1 PropertyChanged	156
8.47 CustomCalibrationLibrary.ViewModels.Monitor Class Reference	157
8.47.1 Detailed Description	157
8.47.2 Constructor & Destructor Documentation	157
8.47.2.1 Monitor()	157
8.47.3 Property Documentation	157
8.47.3.1 Index	157
8.47.3.2 Name	158
8.48 GazeUtilityLibrary.MouseHider Class Reference	158
8.48.1 Detailed Description	158
8.48.2 Constructor & Destructor Documentation	158
8.48.2.1 MouseHider()	158
8.48.3 Member Function Documentation	159
8.48.3.1 HideCursor()	159
8.48.3.2 ShowCursor()	159
8.49 GazeUtilityLibrary.Tracker.MouseTracker Class Reference	159
8.49.1 Detailed Description	161

8.49.2 Constructor & Destructor Documentation	161
8.49.2.1 MouseTracker()	161
8.49.3 Member Function Documentation	162
8.49.3.1 ApplyCalibration()	162
8.49.3.2 CollectCalibrationDataAsync()	162
8.49.3.3 CollectValidationDataAsync()	162
8.49.3.4 ComputeValidation()	163
8.49.3.5 Dispose()	163
8.49.3.6 FinishCalibration()	163
8.49.3.7 FinishCalibrationAsync()	164
8.49.3.8 FinishValidation()	164
8.49.3.9 GetFixationFrameCount()	164
8.49.3.10 GetUnitDirection()	164
8.49.3.11 InitCalibration()	165
8.49.3.12 InitCalibrationAsync()	165
8.49.3.13 InitDriftCompensation()	165
8.49.3.14 InitValidation()	165
8.49.3.15 Start()	166
8.49.3.16 Stop()	166
8.50 GazeControlLibrary.PipeCommand Class Reference	166
8.50.1 Detailed Description	166
8.50.2 Constructor & Destructor Documentation	166
8.50.2.1 PipeCommand()	166
8.50.3 Property Documentation	167
8.50.3.1 Command	167
8.50.3.2 Label	167
8.50.3.3 ResetStartTime	167
8.50.3.4 TrialId	167
8.51 CustomCalibrationLibrary.Converters.PositionConverter Class Reference	168
8.51.1 Detailed Description	169
8.51.2 Member Function Documentation	169
8.51.2.1 Convert()	169
8.51.2.2 ConvertBack()	169
8.51.3 Member Data Documentation	170
8.51.3.1 OffsetProperty	170
8.51.4 Property Documentation	170
8.51.4.1 Offset	170
8.52 CustomCalibrationLibrary.Converters.ProximityColorConverter Class Reference	171
8.52.1 Detailed Description	171
8.52.2 Member Function Documentation	172
8.52.2.1 Convert()	172
8.52.2.2 ConvertBack()	172

8.53 GazeUtilityLibrary.ScreenArea Class Reference	173
8.53.1 Detailed Description	174
8.53.2 Constructor & Destructor Documentation	174
8.53.2.1 ScreenArea()	174
8.53.3 Member Function Documentation	174
8.53.3.1 Dump()	174
8.53.3.2 GetIntersectionPoint()	175
8.53.3.3 GetPoint2d()	175
8.53.3.4 GetPoint2dNormalized()	175
8.53.4 Property Documentation	177
8.53.4.1 BottomLeft	177
8.53.4.2 BottomRight	177
8.53.4.3 Center	177
8.53.4.4 Height	177
8.53.4.5 TopLeft	178
8.53.4.6 TopRight	178
8.53.4.7 Width	178
8.54 Tobii.Research.Addons.ScreenBasedCalibrationValidation Class Reference	178
8.54.1 Detailed Description	180
8.54.2 Member Enumeration Documentation	180
8.54.2.1 ValidationState	180
8.54.3 Constructor & Destructor Documentation	180
8.54.3.1 ScreenBasedCalibrationValidation()	180
8.54.4 Member Function Documentation	181
8.54.4.1 Compute()	181
8.54.4.2 DiscardData()	181
8.54.4.3 Dispose()	181
8.54.4.4 EnterValidationMode()	181
8.54.4.5 LeaveValidationMode()	182
8.54.4.6 StartCollectingData()	182
8.54.4.7 ToString()	182
8.54.5 Property Documentation	182
8.54.5.1 Result	182
8.54.5.2 State	183
8.55 CustomCalibrationLibrary.Views.ScreenSelection Class Reference	183
8.55.1 Detailed Description	184
8.55.2 Constructor & Destructor Documentation	184
8.55.2.1 ScreenSelection()	184
8.56 CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel Class Reference	184
8.56.1 Detailed Description	185
8.56.2 Constructor & Destructor Documentation	185
8.56.2.1 ScreenSelectionViewModel()	185

8.56.3 Property Documentation	185
8.56.3.1 CalibrationAbortCommand	185
8.56.3.2 CalibrationStartCommand	185
8.56.3.3 Monitors	186
8.56.3.4 ScreenSwitchCommand	186
8.57 GazeUtilityLibrary.ScreenTriangle Class Reference	186
8.57.1 Detailed Description	187
8.57.2 Constructor & Destructor Documentation	187
8.57.2.1 ScreenTriangle()	187
8.57.3 Member Function Documentation	187
8.57.3.1 GetIntersectionPoint()	187
8.57.4 Property Documentation	188
8.57.4.1 E1	188
8.57.4.2 E2	188
8.57.4.3 V1	188
8.57.4.4 V2	188
8.57.4.5 V3	188
8.58 GazeUtilityLibrary.TrackerLogger Class Reference	189
8.58.1 Detailed Description	189
8.58.2 Constructor & Destructor Documentation	189
8.58.2.1 TrackerLogger()	189
8.58.3 Member Function Documentation	189
8.58.3.1 Debug()	189
8.58.3.2 DumpFatal()	190
8.58.3.3 Error()	190
8.58.3.4 Info()	190
8.58.3.5 Warning()	191
8.59 GazeUtilityLibrary.TrackerMessageBox Class Reference	191
8.59.1 Detailed Description	192
8.60 GazeUtilityLibrary.DataStructs.UserPositionData Class Reference	192
8.60.1 Detailed Description	193
8.60.2 Constructor & Destructor Documentation	193
8.60.2.1 UserPositionData() [1/2]	193
8.60.2.2 UserPositionData() [2/2]	193
8.60.3 Property Documentation	194
8.60.3.1 XCoordLeft	194
8.60.3.2 XCoordRight	194
8.60.3.3 YCoordLeft	194
8.60.3.4 YCoordRight	194
8.60.3.5 ZCoordLeft	195
8.60.3.6 ZCoordRight	195
8.60.4 Event Documentation	195

8.60.4.1 PropertyChanged	195
8.61 CustomCalibrationLibrary.Views.UserPositionGuide Class Reference	195
8.61.1 Detailed Description	196
8.61.2 Constructor & Destructor Documentation	196
8.61.2.1 UserPositionGuide()	196
8.62 CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel Class Reference	197
8.62.1 Detailed Description	197
8.62.2 Constructor & Destructor Documentation	197
8.62.2.1 UserPositionGuideViewModel()	197
8.62.3 Property Documentation	197
8.62.3.1 CalibrationAbortCommand	198
8.62.3.2 CalibrationStartCommand	198
8.62.3.3 UserPosition	198
8.63 CustomCalibrationLibrary.Views.ValidationResult Class Reference	198
8.63.1 Detailed Description	199
8.63.2 Constructor & Destructor Documentation	199
8.63.2.1 ValidationResult()	199
8.64 CustomCalibrationLibrary.ViewModels.ValidationResultViewModel Class Reference	200
8.64.1 Detailed Description	200
8.64.2 Constructor & Destructor Documentation	200
8.64.2.1 ValidationResultViewModel()	200
8.64.3 Property Documentation	200
8.64.3.1 ValidationCloseCommand	201
8.64.3.2 ValidationData	201
8.64.3.3 ValidationRestartCommand	201
Index	203

Chapter 1

Changelog

v3.4.2

Improvements

- Seperate pipe command handler functionality to `NamedPipeClient` class and move `NamedPipeClient` to seperate project based on .NET 4.72. This allows for a direct integration in python with `pythonnet`.

Bug Fixes

- Log `GazeControl` messages to seperate file to avoid race conditions (#1).
- Fix pipe connection problems on startup (#5).
- Only hide the mouse when calibration dots are shown (#7).
- Fix dispersion calculation in drift compensation (#9).

v3.4.1

Improvements

- Add note on foreground behaviour to readme.

Bug Fixes

- Fix validation when only one screen is connected.

v3.4.0

New Features

- Add configuration option `DriftCompensationWindowShow` to enable or disable the drift compensation window.
- Add configuration option `DriftCompensationDurationThreshold` to configure the required fixation time during drift compensation.
- Add configuration option `DriftCompensationWindowShow` to enable or disable the drift compensation window.
- Add configuration option `DriftCompensationDispersionThresholdMax` which allows to define maximal allowed deviation angle during drift compensation: If the computed compensation angle is larger as the here configured angle the drift compensation is not updated.
- During `Gaze.exe` startup check for already running `Gaze` processes and kill them.

Improvements

- Update default configuration settings to allow starting the application without error.
- Improvements to `GazeControl.exe` and annotations:
 - Introduce arguments `/label` and `/trialId` to pass a label and the trial ID.
 - Make argument `/command` optional to allow setting annotations without a command.
 - Set trialID and label based on capture timestamp instead of system timestamp.
- Improve timestamps: Compute the tracker latency for each sample and use this to dump the system time of data capture and the system time of data reading.
- Update opensesame templates. This might help to reduce taskbar flickering and makes the code more maintainable.

Changes

- Remove argument `/value` in `GazeControl.exe` because it became obsolete with the new arguments `/label` and `trialId`.
- Rename configuration option `DispersionThreshold` to `DriftCompensationDispersionThreshold`.
- Hide window icon in taskbar. This might help to avoid taskbar flickering during experimentation.
- Save data files with `csv` extension instead of `txt`.
- By default, use `,` delimiter instead of `\t` delimiter when dumping values.

Bug Fixes

- Fix dispersion computation during drift compensation.

v3.3.2

New Features

- Allow to select the screen on a multi-screen setup with the keyboard.

Improvements

- Improve opensesame template files and add a template for version 3.3 and 4.0.
- Improve ztree template file
- Add a documentation to the sample folder

Bug Fixes

- Move calibration, validation and drift compensation windows to the foreground.

v3.3.1

Improvements

- Improve performance.
- Dump validation results for each validation point.

Bug Fixes

- Represent the relative timestamp in total milliseconds instead of a timespan.

v3.3.0

New Features

- Add pipe command `SET_TRIAL_ID` to allow to annotate data samples.
- On multi-screen setups start calibration with a screen selection page.
- Add Screen Area coordinates to the dumped configuration file.
- Add gaze validation which can be started through the command `VALIDATE`.

Improvements

- Fix compiler warnings.
- Extend helper scripts.

v3.2.0

New Features

- Add relative timestamp to output data.
- Add annotation tag to output data.
- Add pipe command `SET_TAG` to allow annotate data samples.
- Add pipe command `RESET_START_TIMER` to reset the relative timestamp.
- Add a log entry of the version of the gaze application.
- Add helper scripts to generate shortcuts to `GazeControl.exe`.

v3.1.0

New Features

- Add a custom drift compensation process
- Allow to pass the argument `outputPath` to the application for dynamic output path assignment.

Improvements

- Integrate calibration into `Gaze.exe`
- Remove [Tobii](#) research dependencies from everywhere except the eye tracker device class
- Cleanup and rearrangement of code to improve readability

v3.0.0

New Features

- A custom calibration application is added to the portfolio. This allows to calibrate a device without the need for a 3rd party application.
- Proper shutdown handling of [GazeToMouse](#) through named pipes.
- Allow to enable/disable gaze recording through named pipes.
- Allow to enable/disable mouse tracking through named pipes.
- Allow to pass argument subject to the application.

Improvements

- Update all projects to .NET version 6.0.
- Cleanup code base, split functions into separate libraries.
- Apply MVVM architectural pattern where sensible.

Changes

- Remove [Tobii](#) Interaction Library
- Remove all configuration options for [Tobii](#) Core (only [Tobii](#) Pro SDK is supported)
- Remove [Tobii](#) Core application wrapper (TobiiTest, TobiiGuestCalibrate)
- Use the [Tobii](#) pro eye tracker manager for device calibration instead of the [Tobii](#) Core software.
- Rename [GazeToMouse](#) to `Gaze` and `GazeToMouseClose` to `GazeClose`.

v2.3.0

New Features

- A mouse tracker device can now be used instead of an eyetracker device. The mouse tracker logs the timestamp and the x and y coordinates of the mouse pointer whenever the mouse-move event is fired. The mouse tracker is used when the configuration field 'TrackerDevice' is set to the value 2.

Improvements

- Rename the configuration field 'TobiiSDK' to 'TrackerDevice'.

v2.2.0

New Features

- Configuration file
 - Dump the configurations used for an experiment to a file at the "DataLogPath"
 - Allow to configure an experiment name which is used as a postfix of the dumped configuration file name
 - Consider the config file as invalid if not all required configuration parameters are defined
 - Consider the config file as invalid if unknown parameters are defined
 - Allow to configure whether to log data sets where all data is invalid (eyes closed, no subject in front of the screen, etc)
- Error Handling
 - Attach an error string to the output file, indicating errors that occurred during the run
 - Attach an error string to the dumped configuration file, indicating errors of the configuration

Improvements

- Fall back to Core SDK if the license file cannot be applied to the device

v2.1.0

New Features

- Log eye origin coordinates
 - x, y, z coordinates of the left and the right eye
 - compute distance of the left and right eye to the eyetracker
 - compute the average distance of the two eyes

Improvements

- Check the three format values and the column order individually to produce more specific log entries

v2.0.1

Bug Fix

- with SDK Pro, use system timestamp to cope with disconnected device
- fix the path in the z-tree sample file

v2.0.0

New Features

- Support for [Tobii](#) Pro SDK
 - apply license to eyetracker device at startup
 - logging of pupil diameter
 - logging of individual eye data
- Allow to configure column headers of output file

Improvements

- Improved configuration options for the output file

v1.0.0

New Features

- Notify user with popup if eyetracker is not ready
- Allow to configure time interval for the software to wait for the eyetracker to become ready

Improvements

- Rename default output file for data from `<prefix>_data.txt` to `<prefix>_gaze.txt`

v0.3.2

Improvements

- add header to the data log file.
- change the default value of allowed gaze data files.
- check and wait for ready state of the eye tracker before performing operations with it.

Bug Fix

- create a log file per machine to prevent concurrency conflicts.

v0.3.1

Improvements

- ignore the option "HideMouse" when "ControlMouse" is disabled.

Bug Fix

- remove double log entry of mouse hiding and restoring event.

v0.3.0

New Features

- allow to configure whether the gaze data is logged.
- allow to configure the maximum allowed amount of gaze data files in the output folder. Oldest files are deleted first.

Improvements

- limit the logfile size to 1MB. If the size is exceeded a new file is created. At any time only two log files are allowed, The older file is overwritten once both files exceed 1MB.

v0.2.0

New Features

- allow to configure whether the mouse is controlled by the gaze of the subject or not.
- allow to configure the output format of the gaze data.

v0.1.0

First release of the [GazeToMouse](#) toolset.

The toolset was tested on **Windows 7** in conjunction with **ztree v3.6.7** and **Tobii Eye Tracking Core v2.11.1.6952**.

Chapter 2

Toolset to Control Tobii Eye Tracker

This repository contains the source code for multiple simple tools that allow to control a [Tobii](#) eye tracker from a 3rd party application. Specifically, this project aims at providing a set of executables that can be called from within [ztree](#) to allow eye tracker support for economic experiments.

For more details please refer to the [documentation](#).

Important Configure the task manager to be always in the foreground (In task manager enable "Options->↔ Always on top"). Why: The application `Gaze.exe` may open windows that are put to the foreground in a very aggressive manner. This is done in order to cope with experimentation software that uses this same behaviour (e.g. Opensesame with psychopy or expyriment backend). If something goes wrong with `Gaze.exe` the user could be locked out from the computer because a window keeps blocking access to the system. With the task manager set to "Always on top" there is a way out.

Installation

The complete toolset package can be downloaded from the [releases](#). The package contains the following executables:

- ****Gaze.exe**** This program uses the [Tobii Pro SDK](#) to extract the gaze position on the screen where the subject is looking at. The extracted data is recorded and stored to a file. Optionally, the mouse cursor position is updated to this position such that the mouse cursor is controlled by the gaze of the subject. Instead of using an eye tracker device it is also possible to simply log the mouse coordinates. ****Gaze.exe**** runs infinitely until it is terminated by an external command. This should **not** be done with a forced kill (e.g. by executing the command `taskkill /F /IM Gaze.exe` or by killing the task with the task manager) because it prevents the program from terminating gracefully. This has several consequences:
 - open files are not closed properly and the data stream is cut off. This can lead to corrupt files.
 - if the feature of hiding the mouse pointer is used, the mouse will remain hidden.
 - memory is not freed properly. Instead the program ****GazeControls.exe /command TERMINATE**** should be used.
- ****GazeControl.exe**** This program allows to interact with ****Gaze.exe****. `GazeControl.exe` accepts the following optional arguments:
 - `/reset`: Allows to reset the relative timestamp of the gaze data.
 - `/trialId <ID>`: Sets a trial ID `<ID>` which will be added to each data sample in the output file.**Important:** Make sure that only integer numbers are used as trial ID.

- /label <LABEL>: Sets a custom label <LABEL> which will be added to each data sample in the output file. Any string is accepted here.
- /command <COMMAND>: A command allow to activate/deactivate features of Gaze.exe. The following commands are supported:
 - * CUSTOM_CALIBRATE uses the [Tobii Pro SDK](#) and launches a custom calibration process which allows to calibrate the eye tracker without having to rely on the calibration software provided by [Tobii](#).
 - * VALIDATE uses the [Tobii Pro SDK Addon](#) and launches a validation process.
 - * DRIFT_COMPENSATION launches a custom drift compensation process to compensate gaze drifts that may occur during experimentation.
 - * GAZE_RECORDING_DISABLE requests `**Gaze.exe**` to stop recording gaze data. Gaze.exe will continue to run (and update the mouse pointer if configured accordingly) but no longer store gaze data to the disk.
 - * GAZE_RECORDING_ENABLE requests `**Gaze.exe**` to start recording gaze data.
 - * MOUSE_TRACKING_DISABLE requests `**Gaze.exe**` to stop updating the mouse pointer by the gaze position.
 - * MOUSE_TRACKING_ENABLE requests `**Gaze.exe**` to start updating the mouse pointer by the gaze position.
 - * RESET_DRIFT_COMPENSATION resets the drift compensation computed with the command DRIFT_COMPENSATION.
 - * TERMINATE requests `**Gaze.exe**` to close gracefully and logs these events to the log file.

Multiple arguments can be passed to the application but each argument can only be passed once. Passing an argument to an application can be done in command line or by crating a shortcut to the program. Corresponding shortcuts for all available <COMMAND>s are provided in the release package.

- `**ShowMouse.exe**` This program allows to restore the standard mouse pointer. It might be useful if the program Gaze.exe crashes or is closed forcefully such that the mouse pointer is not restored after terminating. The subject might end up with a hidden mouse pointer. A good solution for such a case is to install a shortcut to ShowMouse.exe on the desktop in order to execute it with the keyboard.

In order to run the executables the following files need to be placed in the same directory as the executables:

- tobii_pro.dll
- tobii_firmware_upgrade.dll
- assets/blank.cur
- config.json

In order to use tje [GazeControlLibrary](#), the following files need to be placed in the same directory as the executables:

- GazeControlLibrary.dll
- Newtonsoft.Json.dll

Further, the [Tobii](#) engine must be running and the eye tracker must be enabled.

Tobii Eye Tracker 4c

To install the driver for the `Tobii Eye Tracker 4c` install `Tobii Experience Driver`.

This will start the following services:

- `Tobii Runtime Service`
- `Tobii Service`

and the following processes:

- `Tobii Interaction Engine`

Tobii Pro Spark

To install the driver for the `Tobii Pro Spark` use the `Tobii Pro Eye Tracker Manager`:

1. Install `Tobii Pro Eye Tracker Manager (ETM)`
2. Connect the `Tobii Pro SPark` device to the computer
3. Install the driver with the ETM

This starts the service `Tobii Pro Spark Runtime`.

Scripts

The folder `scripts` contains two files `CreateShortcut.ps1` and `CreateShortcuts.bat` which allow to create shortcuts to the application `GazeControl.exe` with predefined command arguments. In order to generate the shortcut files perform the following steps:

1. copy the two script files into the installation folder
1. execute the script `CreateShortcuts.bat`

Note that the generated shortcuts are tied to the installation folder. Copying the installation folder to another location will break the links.

3rd Party Applications

This section provides some information on how to run the here provided executables from within 3rd party applications.

ztree

For quick starters, a simple `ztree sample program` is available.

Opensesame

For quick starters, a simple `opensesame sample program` is available.

Release Notes

Information about the releases can be found in the `CHANGELOG`

Chapter 3

Sample Files for Experimentation with Eye Tracker Utility

This folder holds some sample files to use the gaze utility in an experiment management tool.

config.json

A sample configuration file which can be used as a starting point to configure the gaze utility.

config_libgac.json

A sample configuration file which produces minimalistic gaze output data which matches with the example script of `libgac` (a gaze analysis library).

template.osexp

A sample file which demonstrates how to start the gaze utility from openseame. This was tested on opensesame version 3.3.14 and 4.0.5 on Windows.

Note that the application only worked with the PyGame (legacy) backend because otherwise the gaze windows kept being covered by the opensesame fullscreen window.

It might be possible (and potentially a better solution) to manually control the window through python (e.g. with win32gui on Windows or with xdotool on Linux).

template.ztt

A sample file which demonstrates how to start the gaze utility from openseame. The sample file was generated with the ztree version 5.1.11.

Chapter 4

Namespace Index

4.1 Namespace List

Here is a list of all documented namespaces with brief descriptions:

CustomCalibrationLibrary	23
CustomCalibrationLibrary.Commands	23
CustomCalibrationLibrary.Converters	23
CustomCalibrationLibrary.Models	23
CustomCalibrationLibrary.ViewModels	24
CustomCalibrationLibrary.Views	25
GazeControl	25
GazeControlLibrary	25
GazeToMouse	26
GazeUtilityLibrary	
helper class to show and hide the system cursor	26
GazeUtilityLibrary.DataStructs	28
GazeUtilityLibrary.Tracker	30
ShowMouse	30
Tobii	30
Tobii.Research	30
Tobii.Research.Addons	30
Tobii.Research.Addons.Utility	31
TobiiCalibrate	31

Chapter 5

Hierarchical Index

5.1 Class Hierarchy

This inheritance list is sorted roughly, but not completely, alphabetically:

Application	
GazeControl.App	33
GazeToMouse.App	34
ShowMouse.App	38
TobiiCalibrate.App	39
CalibrationPoint	
CustomCalibrationLibrary.ViewModels.CalibrationPointViewModel	76
Tobii.Research.Addons.CalibrationValidationPoint	83
Tobii.Research.Addons.CalibrationValidationResult	86
CustomCalibrationLibrary.ViewModels.CalibrationViewModel	89
CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel	81
GazeUtilityLibrary.ConfigItem	93
GazeUtilityLibrary.ConfigScreenArea	102
DependencyObject	
CustomCalibrationLibrary.Converters.PositionConverter	168
GazeUtilityLibrary.DriftCompensation	106
GazeUtilityLibrary.DataStructs.DriftCompensationData	109
CustomCalibrationLibrary.ViewModels.DriftCompensationViewModel	110
GazeUtilityLibrary.DataStructs.EyeData	113
Frame	
CustomCalibrationLibrary.Views.CalibrationFrame	64
GazeUtilityLibrary.DataStructs.GazeCalibrationData	122
GazeUtilityLibrary.GazeConfiguration	127
GazeUtilityLibrary.DataStructs.GazeData	133
GazeUtilityLibrary.DataStructs.GazeData2d	137
GazeUtilityLibrary.DataStructs.GazeData3d	138
GazeUtilityLibrary.DataStructs.GazeDataCollection	140
GazeUtilityLibrary.GazeError	144
GazeUtilityLibrary.CalibrationDataError	59
GazeUtilityLibrary.GazeConfigError	126
GazeUtilityLibrary.GazeDataError	143
GazeUtilityLibrary.DataStructs.GazeValidationData	145
GazeUtilityLibrary.DataStructs.GazeValidationPoint	149
ICommand	
CustomCalibrationLibrary.Commands.CalibrationCommand	57

IDisposable	
GazeUtilityLibrary.Tracker.BaseTracker	40
GazeUtilityLibrary.Tracker.EyeTrackerPro	114
GazeUtilityLibrary.Tracker.MouseTracker	159
Tobii.Research.Addons.ScreenBasedCalibrationValidation	178
INotifyPropertyChanged	
CustomCalibrationLibrary.Models.CalibrationModel	65
CustomCalibrationLibrary.Views.CalibrationFailed	61
GazeUtilityLibrary.DataStructs.CalibrationPoint	72
GazeUtilityLibrary.DataStructs.LiveGazePoint	155
GazeUtilityLibrary.DataStructs.UserPositionData	192
GazeUtilityLibrary.Tracker.BaseTracker	40
IValueConverter	
CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter	151
CustomCalibrationLibrary.Converters.PositionConverter	168
CustomCalibrationLibrary.Converters.ProximityColorConverter	171
GazeUtilityLibrary.JsonConfigParser	153
CustomCalibrationLibrary.ViewModels.Monitor	157
GazeUtilityLibrary.MouseHider	158
Page	
CustomCalibrationLibrary.Views.Calibration	55
CustomCalibrationLibrary.Views.CalibrationFailed	61
CustomCalibrationLibrary.Views.Computing	92
CustomCalibrationLibrary.Views.Disconnect	105
Page	
CustomCalibrationLibrary.Views.CalibrationResult	78
CustomCalibrationLibrary.Views.ScreenSelection	183
CustomCalibrationLibrary.Views.UserPositionGuide	195
CustomCalibrationLibrary.Views.ValidationResult	198
GazeControlLibrary.PipeCommand	166
GazeUtilityLibrary.ScreenArea	173
CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel	184
GazeUtilityLibrary.ScreenTriangle	186
GazeUtilityLibrary.TrackerLogger	189
UserControl	
CustomCalibrationLibrary.Views.CalibrationPoint	75
CustomCalibrationLibrary.Views.CalibrationResultPoint	79
CustomCalibrationLibrary.Views.FixationPoint	121
CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel	197
CustomCalibrationLibrary.ViewModels.ValidationResultViewModel	200
Window	
CustomCalibrationLibrary.Views.CalibrationWindow	91
CustomCalibrationLibrary.Views.DriftCompensationWindow	112
GazeUtilityLibrary.TrackerMessageBox	191

Chapter 6

Class Index

6.1 Class List

Here are the classes, structs, unions and interfaces with brief descriptions:

GazeControl.App	Interaction logic for App.xaml	33
GazeToMouse.App	Interaction logic for App.xaml	34
ShowMouse.App	Interaction logic for App.xaml	38
TobiiCalibrate.App	Interaction logic for App.xaml	39
GazeUtilityLibrary.Tracker.BaseTracker	The common interface for the Tobii eyetracker Engines Core and Pro	40
CustomCalibrationLibrary.Views.Calibration	Interaction logic for Calibration.xaml	55
CustomCalibrationLibrary.Commands.CalibrationCommand	Comand class to trigger calibration events	57
GazeUtilityLibrary.CalibrationDataError	The calibration data error class to convert error flags to binary strings	59
CustomCalibrationLibrary.Views.CalibrationFailed	Interaction logic for CalibrationFailed.xaml	61
CustomCalibrationLibrary.Views.CalibrationFrame	Interaction logic for CalibrationCollection.xaml	64
CustomCalibrationLibrary.Models.CalibrationModel	The model for the calibration process	65
GazeUtilityLibrary.DataStructs.CalibrationPoint	A calibration point class holding several metrics connected to a calibration point	72
CustomCalibrationLibrary.Views.CalibrationPoint	Interaction logic for CalibrationPoint.xaml	75
CustomCalibrationLibrary.ViewModels.CalibrationPointViewModel	The view model for a calibration point	76
CustomCalibrationLibrary.Views.CalibrationResult	Interaction logic for CalibrationResult.xaml	78
CustomCalibrationLibrary.Views.CalibrationResultPoint	Interaction logic for CalibrationResultPoint.xaml	79
CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel	View model class of the gaze calibration result	81

Tobii.Research.Addons.CalibrationValidationPoint	
Represents a collected point that goes into the calibration validation. It contains calculated values for accuracy and precision as well as the original gaze samples collected for the point	83
Tobii.Research.Addons.CalibrationValidationResult	
Contains the result of the calibration validation	86
CustomCalibrationLibrary.ViewModels.CalibrationViewModel	
The view model class of the calibration view	89
CustomCalibrationLibrary.Views.CalibrationWindow	
Interaction logic for MainWindow.xaml	91
CustomCalibrationLibrary.Views.Computing	
Interaction logic for Computing.xaml	92
GazeUtilityLibrary.ConfigItem	
configuration file class	93
GazeUtilityLibrary.ConfigScreenArea	
The JSON structure of the screen area	102
CustomCalibrationLibrary.Views.Disconnect	
Interaction logic for Disconnect.xaml	105
GazeUtilityLibrary.DriftCompensation	
The class to handle drift compensation	106
GazeUtilityLibrary.DataStructs.DriftCompensationData	
The drift compensation data structure	109
CustomCalibrationLibrary.ViewModels.DriftCompensationViewModel	
The view model class of the drift compensation view	110
CustomCalibrationLibrary.Views.DriftCompensationWindow	
Interaction logic for DriftCompensation.xaml	112
GazeUtilityLibrary.DataStructs.EyeData	
The eye data set, including pupil information	113
GazeUtilityLibrary.Tracker.EyeTrackerPro	
Interface to the Tobii SDK Pro engine	114
CustomCalibrationLibrary.Views.FixationPoint	
Interaction logic for FixationPoint.xaml	121
GazeUtilityLibrary.DataStructs.GazeCalibrationData	
The gaze calibration data structure	122
GazeUtilityLibrary.GazeConfigError	
The gaze config error class to convert error flags to binary strings	126
GazeUtilityLibrary.GazeConfiguration	
The gaze configuration handler	127
GazeUtilityLibrary.DataStructs.GazeData	
The class definition of a gaze data set	133
GazeUtilityLibrary.DataStructs.GazeData2d	
The 2d gaze data set	137
GazeUtilityLibrary.DataStructs.GazeData3d	
The 3d gaze data set	138
GazeUtilityLibrary.DataStructs.GazeDataCollection	
The gaze data set, including 2d and (optionally) 3d gaze data as well as optional eye data . . .	140
GazeUtilityLibrary.GazeDataError	
The gaze data error class to convert error flags to binary strings	143
GazeUtilityLibrary.GazeError	
The base error class to convert error flags to binary strings	144
GazeUtilityLibrary.DataStructs.GazeValidationData	
The gaze validation data structure	145
GazeUtilityLibrary.DataStructs.GazeValidationPoint	
A validation point	149
CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter	
Converts True to Hidden and False to Visible	151
GazeUtilityLibrary.JsonConfigParser	
The config file "config.json" is parsed and its values are attributed to the ConfigItem class . . .	153

GazeUtilityLibrary.DataStructs.LiveGazePoint	
The live gaze point used for verification during the calibration process	155
CustomCalibrationLibrary.ViewModels.Monitor	
A representation of the screen	157
GazeUtilityLibrary.MouseHider	
hide standard mouse pointer and restore it	158
GazeUtilityLibrary.Tracker.MouseTracker	
This class is used to hook into the system mouse events and track the position	159
GazeControlLibrary.PipeCommand	
The JSON structure of a pipe command	166
CustomCalibrationLibrary.Converters.PositionConverter	
Converter class to convert a normalized coordinate to a pixel coordinate	168
CustomCalibrationLibrary.Converters.ProximityColorConverter	
Converter class to convert the proximity of a normalized coordinate to the center point (0.5) into colors	171
GazeUtilityLibrary.ScreenArea	
The class describing the Screen area in 3d and 2d space	173
Tobii.Research.Addons.ScreenBasedCalibrationValidation	
Provides methods and properties for managing calibration validation for screen based eye trackers	178
CustomCalibrationLibrary.Views.ScreenSelection	
Interaction logic for ScreenSelection.xaml	183
CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel	
The view model class for the screen selection view	184
GazeUtilityLibrary.ScreenTriangle	
A class to describe a triangle. This was supposed to be used to construct the ScreenArea but it turned out that it is simpler to work with the screen plane and use the normalised intersection points to check whether the gaze point is outside the screen area	186
GazeUtilityLibrary.TrackerLogger	
Simple logger class	189
GazeUtilityLibrary.TrackerMessageBox	
Interaction logic for TrackerMessageBox.xaml	191
GazeUtilityLibrary.DataStructs.UserPositionData	
The user position to be rendered on the screen	192
CustomCalibrationLibrary.Views.UserPositionGuide	
Interaction logic for UserPositionGuide.xaml	195
CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel	
The view model class for the user position guide view	197
CustomCalibrationLibrary.Views.ValidationResult	
Interaction logic for ValidationResult.xaml	198
CustomCalibrationLibrary.ViewModels.ValidationResultViewModel	
View model class of the gaze validation result	200

Chapter 7

Namespace Documentation

7.1 CustomCalibrationLibrary Namespace Reference

7.2 CustomCalibrationLibrary.Commands Namespace Reference

Classes

- class [CalibrationCommand](#)
Comand class to trigger calibration events.
- class **GazeVisibilityCommand**
Command class to change the gaze visibility
- class **ScreenSwitchCommand**

7.3 CustomCalibrationLibrary.Converters Namespace Reference

Classes

- class [HasDataToVisibilityConverter](#)
Converts True to Hidden and False to Visible
- class [PositionConverter](#)
Converter class to convert a normalized coordinate to a pixel coordinate.
- class [ProximityColorConverter](#)
Converter class to convert the proximito of a normailized coordinate to the center point (0.5) into colors.

7.4 CustomCalibrationLibrary.Models Namespace Reference

Classes

- class [CalibrationModel](#)
The model for the calibration process.

Enumerations

- enum [CalibrationEventType](#) {
 Init, Start, Accept, Restart,
 Abort }
 Events to trigger changes in the calibration process.
- enum [CalibrationStatus](#) {
 ScreenSelection, HeadPosition, DataCollection, Computing,
 CalibrationResult, ValidationResult, Error, Disconnect }
 The status of the calibarion process.

7.4.1 Enumeration Type Documentation

7.4.1.1 CalibrationEventType

enum [CustomCalibrationLibrary.Models.CalibrationEventType](#) [strong]

Events to trigger changes in the calibration process.

7.4.1.2 CalibrationStatus

enum [CustomCalibrationLibrary.Models.CalibrationStatus](#) [strong]

The status of the calibarion process.

7.5 CustomCalibrationLibrary.ViewModels Namespace Reference

Classes

- class [CalibrationPointViewModel](#)
 The view model for a calibration point.
- class [CalibrationResultViewModel](#)
 View model class of the gaze calibration result.
- class [CalibrationViewModel](#)
 The view model class of the calibration view
- class [DriftCompensationViewModel](#)
 The view model class of the drift compensation view.
- class [Monitor](#)
 A representation of the screen.
- class [ScreenSelectionViewModel](#)
 The view model class for the screen selection view.
- class [UserPositionGuideViewModel](#)
 The view model class for the user position guide view.
- class [ValidationResultViewModel](#)
 View model class of the gaze validation result.

7.6 CustomCalibrationLibrary.Views Namespace Reference

Classes

- class [Calibration](#)
Interaction logic for Calibration.xaml
- class [CalibrationFailed](#)
Interaction logic for CalibrationFailed.xaml
- class [CalibrationFrame](#)
Interaction logic for CalibrationCollection.xaml
- class [CalibrationPoint](#)
Interaction logic for CalibrationPoint.xaml
- class [CalibrationResult](#)
Interaction logic for CalibrationResult.xaml
- class [CalibrationResultPoint](#)
Interaction logic for CalibrationResultPoint.xaml
- class [CalibrationWindow](#)
Interaction logic for MainWindow.xaml
- class [Computing](#)
Interaction logic for Computing.xaml
- class [Disconnect](#)
Interaction logic for Disconnect.xaml
- class [DriftCompensationWindow](#)
Interaction logic for DriftCompensation.xaml
- class [FixationPoint](#)
Interaction logic for FixationPoint.xaml
- class [ScreenSelection](#)
Interaction logic for ScreenSelection.xaml
- class [UserPositionGuide](#)
Interaction logic for UserPositionGuide.xaml
- class [ValidationResult](#)
Interaction logic for ValidationResult.xaml

7.7 GazeControl Namespace Reference

Classes

- class [App](#)
Interaction logic for App.xaml

7.8 GazeControlLibrary Namespace Reference

Classes

- class **NamedPipeClient**
The named pipe client handler.
- class [PipeCommand](#)
The JSON structure of a pipe command.

Enumerations

- enum [LogLevel](#) { **Debug**, **Info**, **Warning**, **Error** }
Defines the available log levels.

7.8.1 Enumeration Type Documentation

7.8.1.1 LogLevel

```
enum GazeControllLibrary.LogLevel [strong]
```

Defines the available log levels.

7.9 GazeToMouse Namespace Reference

Classes

- class [App](#)
Interaction logic for App.xaml

7.10 GazeUtilityLibrary Namespace Reference

helper class to show and hide the system cursor

Classes

- class [CalibrationDataError](#)
The calibration data error class to convert error flags to binary strings.
- class **ConfigChecker**
Helper class to check for the validity of configuration options.
- class [ConfigItem](#)
configuration file class
- class [ConfigScreenArea](#)
The JSON structure of the screen area.
- class [DriftCompensation](#)
The class to handle drift compensation.
- class [GazeConfigError](#)
The gaze config error class to convert error flags to binary strings.
- class [GazeConfiguration](#)
The gaze configuration handler.
- class [GazeDataError](#)
The gaze data error class to convert error flags to binary strings.
- class [GazeError](#)

- The base error class to convert error flags to binary strings.*

 - class [JsonConfigParser](#)

The config file "config.json" is parsed and its values are attributed to the [ConfigItem](#) class.
 - class [MouseHider](#)

hide standard mouse pointer and restore it
 - class [ScreenArea](#)

The class describing the Screen area in 3d and 2d space.
 - class [ScreenTriangle](#)

A class to describe a triangle. This was supposed to be used to construct the [ScreenArea](#) but it turned out that it is simpler to work with the screen plane and use the normalised intersection points to check wheter the gaze point is outside the screen area.
 - class [TrackerLogger](#)

Simple logger class.
 - class [TrackerMessageBox](#)

Interaction logic for TrackerMessageBox.xaml

Enumerations

- enum [EOutputType](#) { **gaze**, **calibration**, **validation**, **control** }

A list of output files.
- enum [EGazeConfigError](#) {
FallbackToDefaultConfigName = 0x001, **FallbackToCurrentOutputDir** = 0x002, **FallbackToDefaultConfig** = 0x004, **FallbackToDefaultDiameterFormat** = 0x008, **FallbackToDefaultOriginFormat** = 0x010, **FallbackToDefaultTimestampFormat** = 0x020, **OmitColumnTitles** = 0x040, **FallbackToDefaultColumnOrder** = 0x080, **FallbackToDefaultNormalizedFormat** = 0x100 }

Error values of the configuration
- enum [EGazeDataError](#) { **FallbackToMouse** = 0x01, **DeviceInterrupt** = 0x02 }

Error values of the gaze output data
- enum [ECalibrationDataError](#) { **DeviceNotSupported** = 0x01, **DeviceInterrupt** = 0x02 }

Error values of the gaze output data

7.10.1 Detailed Description

helper class to show and hide the system cursor

7.10.2 Enumeration Type Documentation

7.10.2.1 ECalibrationDataError

```
enum GazeUtilityLibrary.ECalibrationDataError [strong]
```

Error values of the gaze output data

7.10.2.2 EGazeConfigError

enum [GazeUtilityLibrary.EGazeConfigError](#) [strong]

Error values of the configuration

7.10.2.3 EGazeDataError

enum [GazeUtilityLibrary.EGazeDataError](#) [strong]

Error values of the gaze output data

7.10.2.4 EOutputType

enum [GazeUtilityLibrary.EOutputType](#) [strong]

A list of output files.

7.11 GazeUtilityLibrary.DataStructs Namespace Reference

Classes

- class [CalibrationPoint](#)
A calibration point class holding several metrics connected to a calibration point.
- class [DriftCompensationData](#)
The drift compensation data structure
- class [EyeData](#)
The eye data set, including pupil information.
- class [GazeCalibrationData](#)
The gaze calibration data structure
- class [GazeData](#)
The class definition of a gaze data set
- class [GazeData2d](#)
The 2d gaze data set.
- class [GazeData3d](#)
The 3d gaze data set.
- class [GazeDataCollection](#)
The gaze data set, including 2d and (optionally) 3d gaze data as well as optional eye data.
- class [GazeDataConverter](#)
Convert values to strings according to a format.
- class [GazeValidationData](#)
The gaze validation data structure
- class [GazeValidationPoint](#)
A validation point.
- class [LiveGazePoint](#)
The live gaze point used for verification during the calibration process.
- class [UserPositionData](#)
The user position to be rendered on the screen.

Enumerations

- enum [GazeOutputValue](#) {
DateTimeStamp = 0, DateTimeStampReceived, DateTimeStampRelative, TrialId,
Tag, CombinedGazePoint2dCompensatedX, CombinedGazePoint2dCompensatedY, Combined↵
GazePoint2dX,
CombinedGazePoint2dY, CombinedGazePoint2dIsValid, CombinedGazePoint3dCompensatedX,
CombinedGazePoint3dCompensatedY,
CombinedGazePoint3dCompensatedZ, CombinedGazePoint3dX, CombinedGazePoint3dY, Combined↵
GazePoint3dZ,
CombinedGazePoint3dIsValid, CombinedGazeOrigin3dX, CombinedGazeOrigin3dY, Combined↵
GazeOrigin3dZ,
CombinedGazeOrigin3dIsValid, CombinedGazeDistance, CombinedPupilDiameter, CombinedPupil↵
DiameterIsValid,
LeftGazePoint2dX, LeftGazePoint2dY, LeftGazePoint2dIsValid, LeftGazePoint3dX,
LeftGazePoint3dY, LeftGazePoint3dZ, LeftGazePoint3dIsValid, LeftGazeOrigin3dX,
LeftGazeOrigin3dY, LeftGazeOrigin3dZ, LeftGazeOrigin3dIsValid, LeftGazeDistance,
LeftPupilDiameter, LeftPupilDiameterIsValid, RightGazePoint2dX, RightGazePoint2dY,
RightGazePoint2dIsValid, RightGazePoint3dX, RightGazePoint3dY, RightGazePoint3dZ,
RightGazePoint3dIsValid, RightGazeOrigin3dX, RightGazeOrigin3dY, RightGazeOrigin3dZ,
RightGazeOrigin3dIsValid, RightGazeDistance, RightPupilDiameter, RightPupilDiameterIsValid }
enumerates output values produced by the eyetracker
- enum [CalibrationOutputValue](#) {
Point2dX, Point2dY, LeftGazePoint2dX, LeftGazePoint2dY,
LeftGazePoint2dIsValid, RightGazePoint2dX, RightGazePoint2dY, RightGazePoint2dIsValid }
enumerates output values produced by the eyetracker
- enum [ValidationOutputValue](#) {
Point2dX, Point2dY, LeftAccuracy, LeftPrecision,
LeftPrecisionRMS, RightAccuracy, RightPrecision, RightPrecisionRMS }
enumerates output values produced by the eyetracker

7.11.1 Enumeration Type Documentation

7.11.1.1 CalibrationOutputValue

enum [GazeUtilityLibrary.DataStructs.CalibrationOutputValue](#) [strong]

enumerates output values produced by the eyetracker

7.11.1.2 GazeOutputValue

enum [GazeUtilityLibrary.DataStructs.GazeOutputValue](#) [strong]

enumerates output values produced by the eyetracker

7.11.1.3 ValidationOutputValue

enum [GazeUtilityLibrary.DataStructs.ValidationOutputValue](#) [strong]

enumerates output values produced by the eyetracker

7.12 GazeUtilityLibrary.Tracker Namespace Reference

Classes

- class [BaseTracker](#)
The common interface for the [Tobii](#) eyetracker Engines Core and Pro
- class [EyeTrackerPro](#)
Interface to the [Tobii](#) SDK Pro engine
- class [MouseTracker](#)
This class is used to hook into the system mouse events and track the position

7.13 ShowMouse Namespace Reference

Classes

- class [App](#)
Interaction logic for App.xaml

7.14 Tobii Namespace Reference

7.15 Tobii.Research Namespace Reference

7.16 Tobii.Research.Addons Namespace Reference

Classes

- class [CalibrationValidationPoint](#)
Represents a collected point that goes into the calibration validation. It contains calculated values for accuracy and precision as well as the original gaze samples collected for the point.
- class [CalibrationValidationResult](#)
Contains the result of the calibration validation.
- class [ScreenBasedCalibrationValidation](#)
Provides methods and properties for managing calibration validation for screen based eye trackers.

7.17 Tobii.Research.Addons.Utility Namespace Reference

Classes

- class **Extensions**

Extensions with some operations on Point3D and NormalizedPoint2D among other things.

- class **TimeKeeper**

7.18 TobiiCalibrate Namespace Reference

Classes

- class [App](#)

Interaction logic for App.xaml


Chapter 8

Class Documentation


8.1 GazeControl.App Class Reference

Interaction logic for App.xaml

Inheritance diagram for GazeControl.App:

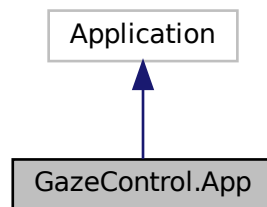

Collaboration diagram for GazeControl.App:

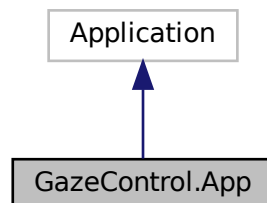

8.1.1 Detailed Description

Interaction logic for App.xaml


The documentation for this class was generated from the following file:

- `source/GazeControl/App.xaml.cs`


8.2 GazeToMouse.App Class Reference

Interaction logic for App.xaml

Inheritance diagram for GazeToMouse.App:

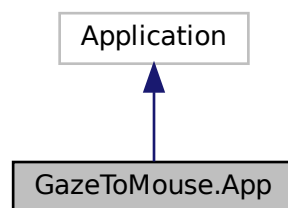

Collaboration diagram for GazeToMouse.App:

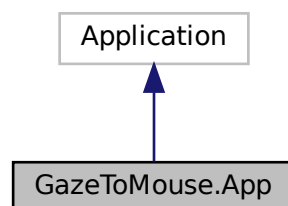

Public Member Functions

- void [GazeRecordingEnable](#) ()
Enable gaze recordings to disk.
- void [GazeRecordingDisable](#) ()
Disable gaze recordings.
- void [MouseTrackingEnable](#) ()
Enable mouse tracking which updates the mouse position to the current gaze point.
- void [MouseTrackingDisable](#) ()
Disable mouse tracking.
- void [ResetDriftCompensation](#) ()
Reset the current drift compensation offset to zero.
- async Task< bool > [CompensateDrift](#) ()
Start the drift compensation process
- async Task< bool > [CustomCalibrate](#) ()
Start the gaze calibration process
- async Task< bool > [CalibrationValidate](#) ()
Start the gaze calibration process
- [App](#) ()
Constructor: initialised logger, gaze configuration, pipe server, and calibration model

Properties

- TimeSpan [StartTime](#) [get, set]
The start time of the application.
- string [LastTag](#) [get, set]
The last tag to annotate gaze data.
- string [Tag](#) [get, set]
An arbitrary tag to annotate gaze data.
- int [TrialId](#) [get, set]
The trial ID to annotate gaze data.

8.2.1 Detailed Description

Interaction logic for App.xaml

8.2.2 Constructor & Destructor Documentation

8.2.2.1 App()

```
GazeToMouse.App.App ( ) [inline]
```

Constructor: initialised logger, gaze configuration, pipe server, and calibration model

8.2.3 Member Function Documentation

8.2.3.1 CalibrationValidate()

```
async Task<bool> GazeToMouse.App.CalibrationValidate ( ) [inline]
```

Start the gaze calibration process

Returns

True on success, false on failure

8.2.3.2 CompensateDrift()

```
async Task<bool> GazeToMouse.App.CompensateDrift ( ) [inline]
```

Start the drift compensation process

Returns

True on success, false on failure

8.2.3.3 CustomCalibrate()

```
async Task<bool> GazeToMouse.App.CustomCalibrate ( ) [inline]
```

Start the gaze calibration process

Returns

True on success, false on failure

8.2.3.4 GazeRecordingDisable()

```
void GazeToMouse.App.GazeRecordingDisable ( ) [inline]
```

Disable gaze recordings.

8.2.3.5 GazeRecordingEnable()

```
void GazeToMouse.App.GazeRecordingEnable ( ) [inline]
```

Enable gaze recordings to disk.

8.2.3.6 MouseTrackingDisable()

```
void GazeToMouse.App.MouseTrackingDisable ( ) [inline]
```

Disable mouse tracking.

8.2.3.7 MouseTrackingEnable()

```
void GazeToMouse.App.MouseTrackingEnable ( ) [inline]
```

Enable mouse tracking which updates the mouse position to the current gaze point.

8.2.3.8 ResetDriftCompensation()

```
void GazeToMouse.App.ResetDriftCompensation ( ) [inline]
```

Reset the current drift compensation offset to zero.

8.2.4 Property Documentation

8.2.4.1 LastTag

```
string GazeToMouse.App.LastTag [get], [set]
```

The last tag to annotate gaze data.

8.2.4.2 StartTime

```
TimeSpan GazeToMouse.App.StartTime [get], [set]
```

The start time of the application.

8.2.4.3 Tag

```
string GazeToMouse.App.Tag [get], [set]
```

An arbitrary tag to annotate gaze data.

8.2.4.4 TrialId

```
int GazeToMouse.App.TrialId [get], [set]
```

The trial ID to annotate gaze data.


The documentation for this class was generated from the following file:

- source/GazeToMouse/App.xaml.cs


8.3 ShowMouse.App Class Reference

Interaction logic for App.xaml

Inheritance diagram for ShowMouse.App:

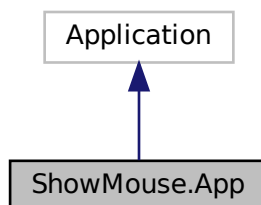

Collaboration diagram for ShowMouse.App:

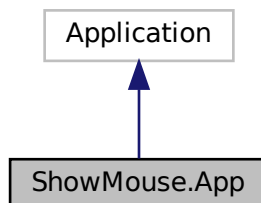

8.3.1 Detailed Description

Interaction logic for App.xaml


The documentation for this class was generated from the following file:

- source/ShowMouse/App.xaml.cs


8.4 TobiiCalibrate.App Class Reference

Interaction logic for App.xaml

Inheritance diagram for TobiiCalibrate.App:

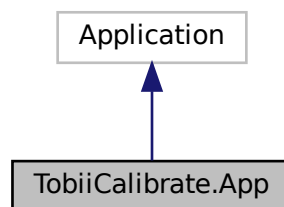

Collaboration diagram for TobiiCalibrate.App:

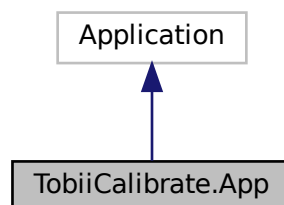

8.4.1 Detailed Description

Interaction logic for App.xaml


The documentation for this class was generated from the following file:

- source/TobiiCalibrate/App.xaml.cs


8.5 GazeUtilityLibrary.Tracker.BaseTracker Class Reference

The common interface for the [Tobii](#) eyetracker Engines Core and Pro

Inheritance diagram for GazeUtilityLibrary.Tracker.BaseTracker:

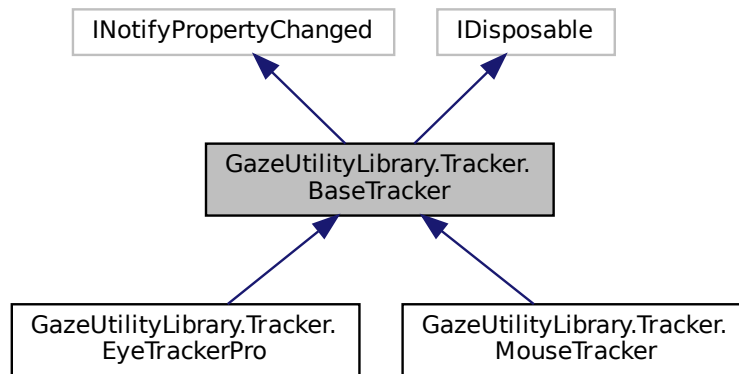

Collaboration diagram for GazeUtilityLibrary.Tracker.BaseTracker:

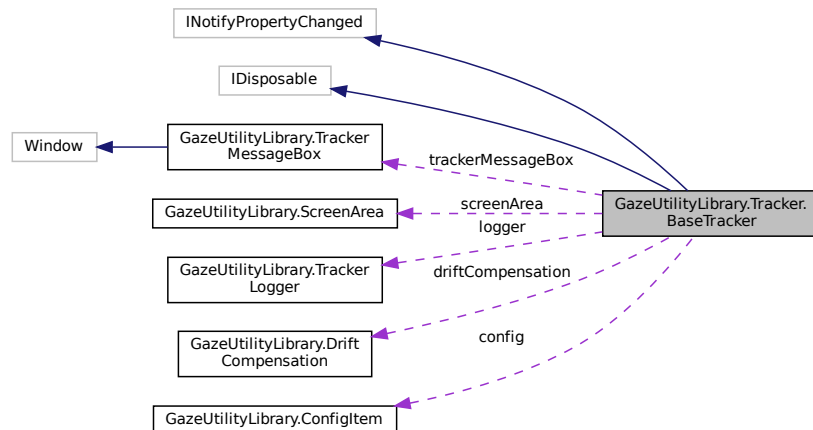

Public Types

- enum `DeviceStatus` {
Configuring, Initializing, InvalidConfiguration, DeviceNotConnected, Tracking }

The tracker device status

Public Member Functions

- delegate void [GazeDataHandler](#) (object sender, [GazeData](#) gazeData)
Event handler for gaze data events of the eyetracker
- delegate void [DriftCompensationEventHandler](#) (object sender, Quaternion [driftCompensation](#))
Event handler for drift compensation events
- delegate void [UserPositionDataHandler](#) (object sender, [UserPositionData](#) e)
Event handler for user position data events of the eyetracker
- [BaseTracker](#) ([TrackerLogger](#) logger, [ConfigItem](#) config, string deviceName)
Initializes a new instance of the EyeTrackerHandler class.
- void [Dispose](#) ()
Performs application-defined tasks associated with freeing, releasing, or resetting unmanaged resources.
- virtual string [PatternReplace](#) (string pattern)
Replaces a patten string with information from the eye tracker. This is device specific and may be overwritten by the device class.
- abstract Task [InitCalibrationAsync](#) ()
Initialise the async calibartion process. This is device specific and must be overwritten by the device class.
- abstract void [InitCalibration](#) ()
Initialise the calibartion process. This is device specific and must be overwritten by the device class.
- abstract void [InitValidation](#) ()
Initialise the validation process. This is device specific and must be overwritten by the device class.
- abstract Task [FinishCalibrationAsync](#) ()
Finish the async calibartion process. This is device specific and must be overwritten by the device class.
- abstract void [FinishCalibration](#) ()
Finish the calibartion process. This is device specific and must be overwritten by the device class.
- abstract void [FinishValidation](#) ()
Finish the validation process. This is device specific and must be overwritten by the device class.
- abstract Task< List< [GazeCalibrationData](#) > > [ApplyCalibration](#) ()
Apply the calibration data. This is device specific and must be overwritten by the device class.
- abstract ? [GazeValidationData](#) [ComputeValidation](#) ()
Apply the validation data. This is device specific and must be overwritten by the device class.
- abstract Task< bool > [CollectCalibrationDataAsync](#) (Point point)
Collect calibration data on a calibration point. This is device specific and must be overwritten by the device class.
- abstract Task< bool > [CollectValidationDataAsync](#) (Point point)
Collect validation data on a validation point. This is device specific and must be overwritten by the device class.
- void [StartDriftCompensation](#) ()
Start the drift compensation process.
- void [ResetDriftCompensation](#) ()
Reset the drift compensation value
- virtual bool [IsInitialised](#) ()
Checks wheter the device is connected and initialised. This is device specific and may be overwritten. Otherwise true is always returned.

Public Attributes

- readonly string [DeviceName](#)
The name of the tracker device

Protected Member Functions

- abstract void [InitDriftCompensation](#) ()
Initialise the drift compensation. This is device specific and must be overwritten by the device class.
- abstract int [GetFixationFrameCount](#) (int durationThreshold)
Get the number of required gaze samples to compute a fixation. This is device specific and must be overwritten by the device because the duration of fixation point detection depends on the frame rate of the device.
- abstract Vector3 [GetUnitDirection](#) ()
Get the unit vector pointing in the direction of the gaze vector. This is device specific as the gaze data are represented in a coordinate system as defined by the device.
- virtual void [Dispose](#) (bool disposing)
Releases unmanaged and - optionally - managed resources.
- bool [IsReady](#) ()
Determines whether this eye tracker is ready.
- virtual void [OnGazeDataReceived](#) ([GazeData](#) gazeData)
Called when [gaze data received].
- virtual void [OnUserPositionDataReceived](#) ([UserPositionData](#) e)
Called when [user position data received].
- virtual void [OnPropertyChanged](#) (string property_name)
Called when when the state property of EyeTracker is changing.
- virtual void [OnTrackerDisabled](#) (EventArgs e)
Raises the E:TrackerDisabled event.
- void [OnTrackerDisabledTimeout](#) (object? source, ElapsedEventArgs e)
Called after a specified amount of time of the eyetracker not being ready.
- virtual void [OnTrackerEnabled](#) (EventArgs e)
Raises the E:TrackerEnabled event.

Protected Attributes

- Timer? [dialogBoxTimer](#)
Timer to control the apperance of the dialog box
- [TrackerLogger](#) logger
The logger
- [TrackerMessageBox?](#) [trackerMessageBox](#)
The dialog box that is controlled by the dialogBoxTimer
- [DriftCompensation?](#) [driftCompensation](#)
drift compensation handler
- [ScreenArea?](#) [screenArea](#) = null
The screen area structure holding the metrics of the screen in 3d space.
- [ConfigItem](#) config
The gaze configuration item

Properties

- [ScreenArea?](#) [ScreenArea](#) [get]
The screen area structure holding the metrics of the screen in 3d space.
- [DeviceStatus](#) [State](#) [get, set]
Gets or sets the state of the eyetracker device.

Events

- EventHandler? [TrackerEnabled](#)
Occurs when [tracker enabled].
- EventHandler? [TrackerDisabled](#)
Occurs when [tracker disabled].
- PropertyChangedEventHandler? [PropertyChanged](#)
Occurs when a property value changes.
- [GazeDataHandler](#)? [GazeDataReceived](#)
Occurs when [gaze data received].
- [DriftCompensationEventHandler](#)? [DriftCompensationComputed](#)
Occurs when drift compensation was computed.
- [UserPositionDataHandler](#)? [UserPositionDataReceived](#)
Occurs when [user position data received].

8.5.1 Detailed Description

The common interface for the [Tobii](#) eyetracker Engines Core and Pro

See also

[INotifyPropertyChanged](#), [IDisposable](#)

8.5.2 Member Enumeration Documentation

8.5.2.1 DeviceStatus

```
enum GazeUtilityLibrary.Tracker.BaseTracker.DeviceStatus [strong]
```

The tracker device status

8.5.3 Constructor & Destructor Documentation

8.5.3.1 BaseTracker()

```
GazeUtilityLibrary.Tracker.BaseTracker.BaseTracker (
    TrackerLogger logger,
    ConfigItem config,
    string deviceName ) [inline]
```

Initializes a new instance of the EyeTrackerHandler class.

Parameters

<i>logger</i>	The logger.
<i>config</i>	The configuration object.
<i>deviceName</i>	Name of the device.

8.5.4 Member Function Documentation

8.5.4.1 ApplyCalibration()

```
abstract Task<List<GazeCalibrationData> > GazeUtilityLibrary.Tracker.BaseTracker.ApplyCalibration ( ) [pure virtual]
```

Apply the calibration data. This is device specific and must be overwritten by the device class.

Returns

The calibration data result wrapped by an async handler.

Implemented in [GazeUtilityLibrary.Tracker.EyeTrackerPro](#), and [GazeUtilityLibrary.Tracker.MouseTracker](#).

8.5.4.2 CollectCalibrationDataAsync()

```
abstract Task<bool> GazeUtilityLibrary.Tracker.BaseTracker.CollectCalibrationDataAsync ( Point point ) [pure virtual]
```

Collect calibration data on a calibration point. This is device specific and must be overwritten by the device class.

Parameters

<i>point</i>	The calibration point for which to collect data
--------------	---

Returns

True on success, false on failure, wrapped by an async handler.

Implemented in [GazeUtilityLibrary.Tracker.MouseTracker](#), and [GazeUtilityLibrary.Tracker.EyeTrackerPro](#).

8.5.4.3 CollectValidationDataAsync()

```
abstract Task<bool> GazeUtilityLibrary.Tracker.BaseTracker.CollectValidationDataAsync (
    Point point ) [pure virtual]
```

Collect validation data on a validation point. This is device specific and must be overwritten by the device class.

Parameters

<i>point</i>	The calibration point for which to collect data
--------------	---

Returns

True on success, false on failure, wrapped by an async handler.

Implemented in [GazeUtilityLibrary.Tracker.MouseTracker](#), and [GazeUtilityLibrary.Tracker.EyeTrackerPro](#).

8.5.4.4 ComputeValidation()

```
abstract ? GazeValidationData GazeUtilityLibrary.Tracker.BaseTracker.ComputeValidation ( )
[pure virtual]
```

Apply the validation data. This is device specific and must be overwritten by the device class.

Returns

The validation data result.

Implemented in [GazeUtilityLibrary.Tracker.EyeTrackerPro](#), and [GazeUtilityLibrary.Tracker.MouseTracker](#).

8.5.4.5 Dispose() [1/2]

```
void GazeUtilityLibrary.Tracker.BaseTracker.Dispose ( ) [inline]
```

Performs application-defined tasks associated with freeing, releasing, or resetting unmanaged resources.

8.5.4.6 Dispose() [2/2]

```
virtual void GazeUtilityLibrary.Tracker.BaseTracker.Dispose (
    bool disposing ) [inline], [protected], [virtual]
```

Releases unmanaged and - optionally - managed resources.

Parameters

<i>disposing</i>	true to release both managed and unmanaged resources; false to release only unmanaged resources.
------------------	--

Reimplemented in [GazeUtilityLibrary.Tracker.MouseTracker](#).

8.5.4.7 DriftCompensationEventHandler()

```
delegate void GazeUtilityLibrary.Tracker.BaseTracker.DriftCompensationEventHandler (
    object sender,
    Quaternion driftCompensation )
```

Event handler for drift compensation events

Parameters

<i>sender</i>	The sender.
<i>driftCompensation</i>	The drift compensation quaternion

8.5.4.8 FinishCalibration()

```
abstract void GazeUtilityLibrary.Tracker.BaseTracker.FinishCalibration ( ) [pure virtual]
```

Finish the calibration process. This is device specific and must be overwritten by the device class.

Implemented in [GazeUtilityLibrary.Tracker.MouseTracker](#), and [GazeUtilityLibrary.Tracker.EyeTrackerPro](#).

8.5.4.9 FinishCalibrationAsync()

```
abstract Task GazeUtilityLibrary.Tracker.BaseTracker.FinishCalibrationAsync ( ) [pure virtual]
```

Finish the async calibration process. This is device specific and must be overwritten by the device class.

Returns

An async handler

Implemented in [GazeUtilityLibrary.Tracker.EyeTrackerPro](#), and [GazeUtilityLibrary.Tracker.MouseTracker](#).

8.5.4.10 FinishValidation()

```
abstract void GazeUtilityLibrary.Tracker.BaseTracker.FinishValidation ( ) [pure virtual]
```

Finish the validation process. This is device specific and must be overwritten by the device class.

Implemented in [GazeUtilityLibrary.Tracker.EyeTrackerPro](#), and [GazeUtilityLibrary.Tracker.MouseTracker](#).

8.5.4.11 GazeDataHandler()

```
delegate void GazeUtilityLibrary.Tracker.BaseTracker.GazeDataHandler (
    object sender,
    GazeData gazeData )
```

Event handler for gaze data events of the eyetracker

Parameters

<i>sender</i>	The sender.
<i>gazeData</i>	The e.

8.5.4.12 GetFixationFrameCount()

```
abstract int GazeUtilityLibrary.Tracker.BaseTracker.GetFixationFrameCount (
    int durationThreshold ) [protected], [pure virtual]
```

Get the number of required gaze samples to compute a fixation. This is device specific and must be overwritten by the device because the duration of fixation point detection depends on the frame rate of the device.

Parameters

<i>durationThreshold</i>	The required fixation duration in milliseconds.
--------------------------	---

Returns

The number of gaze samples to require for fixation detection.

Implemented in [GazeUtilityLibrary.Tracker.EyeTrackerPro](#), and [GazeUtilityLibrary.Tracker.MouseTracker](#).

8.5.4.13 GetUnitDirection()

```
abstract Vector3 GazeUtilityLibrary.Tracker.BaseTracker.GetUnitDirection ( ) [protected],
[pure virtual]
```

Get the unit vector pointing in the direction of the gaze vector. This is device specific as the gaze data are represented in a coordinate system as defined by the device.

Returns

The unit vector

Implemented in [GazeUtilityLibrary.Tracker.EyeTrackerPro](#), and [GazeUtilityLibrary.Tracker.MouseTracker](#).

8.5.4.14 InitCalibration()

```
abstract void GazeUtilityLibrary.Tracker.BaseTracker.InitCalibration ( ) [pure virtual]
```

Initialise the calibration process. This is device specific and must be overwritten by the device class.

Implemented in [GazeUtilityLibrary.Tracker.MouseTracker](#), and [GazeUtilityLibrary.Tracker.EyeTrackerPro](#).

8.5.4.15 InitCalibrationAsync()

```
abstract Task GazeUtilityLibrary.Tracker.BaseTracker.InitCalibrationAsync ( ) [pure virtual]
```

Initialise the async calibration process. This is device specific and must be overwritten by the device class.

Returns

An async handler

Implemented in [GazeUtilityLibrary.Tracker.MouseTracker](#), and [GazeUtilityLibrary.Tracker.EyeTrackerPro](#).

8.5.4.16 InitDriftCompensation()

```
abstract void GazeUtilityLibrary.Tracker.BaseTracker.InitDriftCompensation ( ) [protected],  
[pure virtual]
```

Initialise the drift compensation. This is device specific and must be overwritten by the device class.

Implemented in [GazeUtilityLibrary.Tracker.EyeTrackerPro](#), and [GazeUtilityLibrary.Tracker.MouseTracker](#).

8.5.4.17 InitValidation()

```
abstract void GazeUtilityLibrary.Tracker.BaseTracker.InitValidation ( ) [pure virtual]
```

Initialise the validation process. This is device specific and must be overwritten by the device class.

Implemented in [GazeUtilityLibrary.Tracker.MouseTracker](#), and [GazeUtilityLibrary.Tracker.EyeTrackerPro](#).

8.5.4.18 IsInitialised()

```
virtual bool GazeUtilityLibrary.Tracker.BaseTracker.IsInitialised ( ) [inline], [virtual]
```

Checks wheter the device is connected and initialised. This is device specific and may be overwritten. Otherwise true is always returned.

Returns

True

Reimplemented in [GazeUtilityLibrary.Tracker.EyeTrackerPro](#).

8.5.4.19 IsReady()

```
bool GazeUtilityLibrary.Tracker.BaseTracker.IsReady ( ) [inline], [protected]
```

Determines whether this eye tracker is ready.

Returns

true if this instance is ready; otherwise, false.

8.5.4.20 OnGazeDataReceived()

```
virtual void GazeUtilityLibrary.Tracker.BaseTracker.OnGazeDataReceived (
    GazeData gazeData ) [inline], [protected], [virtual]
```

Called when [gaze data received].

Parameters

<i>gazeData</i>	The gaze data event data.
-----------------	---------------------------

8.5.4.21 OnPropertyChanged()

```
virtual void GazeUtilityLibrary.Tracker.BaseTracker.OnPropertyChanged (
    string property_name ) [inline], [protected], [virtual]
```

Called when when the state property of EyeTracker is changing.

Parameters

<i>property_name</i>	Name of the property in WPF.
----------------------	------------------------------

8.5.4.22 OnTrackerDisabled()

```
virtual void GazeUtilityLibrary.Tracker.BaseTracker.OnTrackerDisabled (
    EventArgs e ) [inline], [protected], [virtual]
```

Raises the E:TrackerDisabled event.

Parameters

<i>e</i>	The EventArgs instance containing the event data.
----------	---

8.5.4.23 OnTrackerDisabledTimeout()

```
void GazeUtilityLibrary.Tracker.BaseTracker.OnTrackerDisabledTimeout (
    object? source,
    ElapsedEventArgs e ) [inline], [protected]
```

Called after a specified amount of time of the eyetracker not being ready.

Parameters

<i>source</i>	The source.
<i>e</i>	The ElapsedEventArgs instance containing the event data.

8.5.4.24 OnTrackerEnabled()

```
virtual void GazeUtilityLibrary.Tracker.BaseTracker.OnTrackerEnabled (
    EventArgs e ) [inline], [protected], [virtual]
```

Raises the E:TrackerEnabled event.

Parameters

<i>e</i>	The EventArgs instance containing the event data.
----------	---

8.5.4.25 OnUserPositionDataReceived()

```
virtual void GazeUtilityLibrary.Tracker.BaseTracker.OnUserPositionDataReceived (
    UserPositionData e ) [inline], [protected], [virtual]
```

Called when [user position data received].

Parameters

<i>e</i>	The gaze data event data.
----------	---------------------------

8.5.4.26 PatternReplace()

```
virtual string GazeUtilityLibrary.Tracker.BaseTracker.PatternReplace (
    string pattern ) [inline], [virtual]
```

Replaces a patten string with information from the eye tracker. This is device specific and may be overwritten by the device class.

Returns

The string where patterns were replaced.

Reimplemented in [GazeUtilityLibrary.Tracker.EyeTrackerPro](#).

8.5.4.27 ResetDriftCompensation()

```
void GazeUtilityLibrary.Tracker.BaseTracker.ResetDriftCompensation ( ) [inline]
```

Reset the drift compensation value

8.5.4.28 StartDriftCompensation()

```
void GazeUtilityLibrary.Tracker.BaseTracker.StartDriftCompensation ( ) [inline]
```

Start the drift compensation process.

8.5.4.29 UserPositionDataHandler()

```
delegate void GazeUtilityLibrary.Tracker.BaseTracker.UserPositionDataHandler (
    object sender,
    UserPositionData e )
```

Event handler for user position data events of the eyetracker

Parameters

<i>sender</i>	The sender.
<i>e</i>	The e.

8.5.5 Member Data Documentation

8.5.5.1 config

```
ConfigItem GazeUtilityLibrary.Tracker.BaseTracker.config [protected]
```


The gaze configuration item

8.5.5.2 DeviceName

```
readonly string GazeUtilityLibrary.Tracker.BaseTracker.DeviceName
```

The name of the tracker device

8.5.5.3 dialogBoxTimer

```
Timer? GazeUtilityLibrary.Tracker.BaseTracker.dialogBoxTimer [protected]
```

Timer to control the apperance of the dialog box

8.5.5.4 driftCompensation

```
DriftCompensation? GazeUtilityLibrary.Tracker.BaseTracker.driftCompensation [protected]
```

drift compensation handler

8.5.5.5 logger

```
TrackerLogger GazeUtilityLibrary.Tracker.BaseTracker.logger [protected]
```

The logger

8.5.5.6 screenArea

```
ScreenArea? GazeUtilityLibrary.Tracker.BaseTracker.screenArea = null [protected]
```

The screen area structure holding the metrics of the screen in 3d space.

8.5.5.7 trackerMessageBox

```
TrackerMessageBox? GazeUtilityLibrary.Tracker.BaseTracker.trackerMessageBox [protected]
```

The dialog box that is controlled by the dialogBoxTimer

8.5.6 Property Documentation

8.5.6.1 ScreenArea

`ScreenArea?` `GazeUtilityLibrary.Tracker.BaseTracker.ScreenArea` `[get]`

The screen area structure holding the metrics of the screen in 3d space.

8.5.6.2 State

`DeviceStatus` `GazeUtilityLibrary.Tracker.BaseTracker.State` `[get]`, `[set]`

Gets or sets the state of the eyetracker device.

The state.

8.5.7 Event Documentation

8.5.7.1 DriftCompensationComputed

`DriftCompensationEventHandler?` `GazeUtilityLibrary.Tracker.BaseTracker.DriftCompensation↔`
`Computed`

Occurs when drift compensation was computed.

8.5.7.2 GazeDataReceived

`GazeDataHandler?` `GazeUtilityLibrary.Tracker.BaseTracker.GazeDataReceived`

Occurs when [gaze data received].

8.5.7.3 PropertyChanged

`PropertyChangedEventHandler?` `GazeUtilityLibrary.Tracker.BaseTracker.PropertyChanged`

Occurs when a property value changes.

8.5.7.4 TrackerDisabled

EventHandler? GazeUtilityLibrary.Tracker.BaseTracker.TrackerDisabled

Occurs when [tracker disabled].

8.5.7.5 TrackerEnabled

EventHandler? GazeUtilityLibrary.Tracker.BaseTracker.TrackerEnabled

Occurs when [tracker enabled].

8.5.7.6 UserPositionDataReceived

UserPositionDataHandler? GazeUtilityLibrary.Tracker.BaseTracker.UserPositionDataReceived

Occurs when [user position data received].


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/Tracker/BaseTracker.cs


8.6 CustomCalibrationLibrary.Views.Calibration Class Reference

Interaction logic for Calibration.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.Calibration:

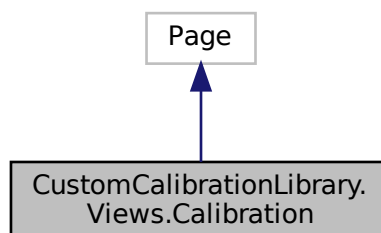

Collaboration diagram for CustomCalibrationLibrary.Views.Calibration:

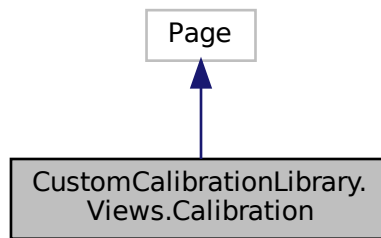

Public Member Functions

- [Calibration](#) (CalibrationModel model)
Initializes a new instance of the [Calibration](#) class.

8.6.1 Detailed Description

Interaction logic for Calibration.xaml

8.6.2 Constructor & Destructor Documentation

8.6.2.1 Calibration()

```
CustomCalibrationLibrary.Views.Calibration.Calibration (
    CalibrationModel model ) [inline]
```

Initializes a new instance of the [Calibration](#) class.

Parameters

<i>model</i>	The calibration model
--------------	-----------------------


The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Views/Calibration.xaml.cs


8.7 CustomCalibrationLibrary.Commands.CalibrationCommand Class Reference

Command class to trigger calibration events.

Inheritance diagram for CustomCalibrationLibrary.Commands.CalibrationCommand:

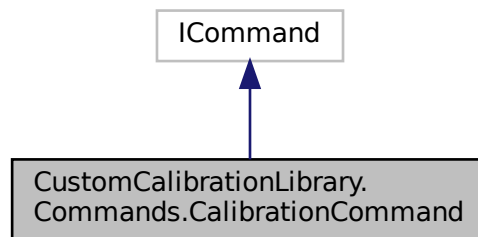

Collaboration diagram for CustomCalibrationLibrary.Commands.CalibrationCommand:

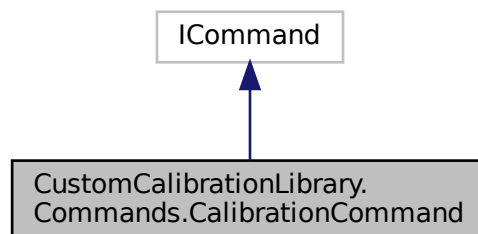

Public Member Functions

- [CalibrationCommand](#) (CalibrationModel model, CalibrationEventType eventType)
Initializes a new instance of the [CalibrationCommand](#) class.
- bool [CanExecute](#) (object? parameter)
Returns whether command can be executed or not.
- void [Execute](#) (object? parameter)
Send calibration event.

Properties

- EventHandler? [CanExecuteChanged](#)
Event handler on can executed flag change.

8.7.1 Detailed Description

Comand class to trigger calibration events.

8.7.2 Constructor & Destructor Documentation

8.7.2.1 CalibrationCommand()

```
CustomCalibrationLibrary.Commands.CalibrationCommand.CalibrationCommand (
    CalibrationModel model,
    CalibrationEventType eventType ) [inline]
```

Initializes a new instance of the [CalibrationCommand](#) class.

Parameters

<i>model</i>	The calibration model
<i>eventType</i>	The type of the calibration event.

8.7.3 Member Function Documentation

8.7.3.1 CanExecute()

```
bool CustomCalibrationLibrary.Commands.CalibrationCommand.CanExecute (
    object? parameter ) [inline]
```

Returns whether command can be executed or not.

Parameters

<i>parameter</i>	The command parameter
------------------	-----------------------

Returns

True

8.7.3.2 Execute()

```
void CustomCalibrationLibrary.Commands.CalibrationCommand.Execute (
    object? parameter ) [inline]
```

Send calibration event.

Parameters

<i>parameter</i>	The command parameter
------------------	-----------------------

8.7.4 Property Documentation

8.7.4.1 CanExecuteChanged

```
EventHandler? CustomCalibrationLibrary.Commands.CalibrationCommand.CanExecuteChanged [add],  
[remove]
```

Event handler on can executed flag change.


The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Commands/CalibrationCommand.cs


8.8 GazeUtilityLibrary.CalibrationDataError Class Reference

The calibration data error class to convert error flags to binary strings.

Inheritance diagram for GazeUtilityLibrary.CalibrationDataError:

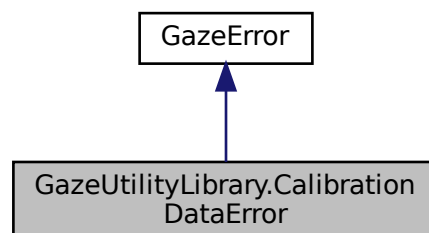

Collaboration diagram for GazeUtilityLibrary.CalibrationDataError:

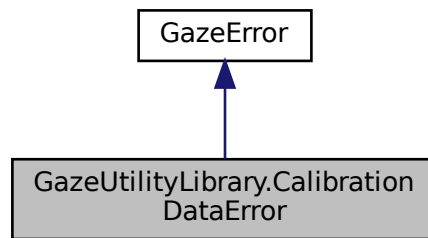

Public Member Functions

- string [GetCalibrationDataErrorString](#) ()
Gets the gaze error string.

Properties

- [ECalibrationDataError Error](#) [set]
The error flags.

Additional Inherited Members

8.8.1 Detailed Description

The calibration data error class to convert error flags to binary strings.

8.8.2 Member Function Documentation

8.8.2.1 GetCalibrationDataErrorString()

```
string GazeUtilityLibrary.CalibrationDataError.GetCalibrationDataErrorString ( ) [inline]
```

Gets the gaze error string.

Returns

the error string with binary error values if errors occurred, the empty string otherwise

8.8.3 Property Documentation

8.8.3.1 Error

`ECalibrationDataError` GazeUtilityLibrary.CalibrationDataError.Error [set]

The error flags.


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/GazeError.cs


8.9 CustomCalibrationLibrary.Views.CalibrationFailed Class Reference

Interaction logic for CalibrationFailed.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.CalibrationFailed:

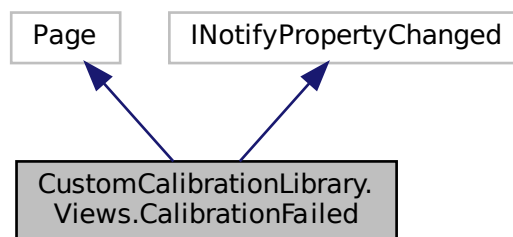

Collaboration diagram for CustomCalibrationLibrary.Views.CalibrationFailed:

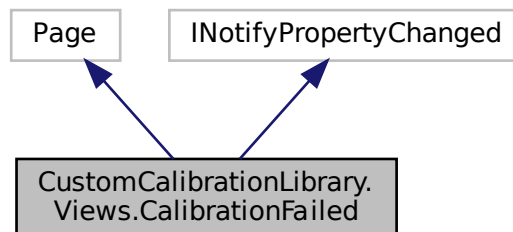

Public Member Functions

- [CalibrationFailed](#) ([CalibrationModel](#) model)
Constructor

Properties

- ICommand [CalibrationRestartCommand](#) [get]
Command to restart the calibration
- ICommand [CalibrationAbortCommand](#) [get]
Command to abort the calibration
- string [Error](#) [get, set]
The error message to be updated on the view.

Events

- PropertyChangedEventHandler? [PropertyChanged](#)
The property change event to update the view.

8.9.1 Detailed Description

Interaction logic for CalibrationFailed.xaml

8.9.2 Constructor & Destructor Documentation

8.9.2.1 CalibrationFailed()

```
CustomCalibrationLibrary.Views.CalibrationFailed.CalibrationFailed (  
    CalibrationModel model ) [inline]
```

Constructor

Parameters

<i>model</i>	The claibration model
--------------	-----------------------

8.9.3 Property Documentation

8.9.3.1 CalibrationAbortCommand

`ICommand CustomCalibrationLibrary.Views.CalibrationFailed.CalibrationAbortCommand [get]`

Command to abort the calibration

8.9.3.2 CalibrationRestartCommand

`ICommand CustomCalibrationLibrary.Views.CalibrationFailed.CalibrationRestartCommand [get]`

Command to restart the calibration

8.9.3.3 Error

`string CustomCalibrationLibrary.Views.CalibrationFailed.Error [get], [set]`

The error message to be updated on the view.

8.9.4 Event Documentation

8.9.4.1 PropertyChanged

`PropertyChangedEventHandler? CustomCalibrationLibrary.Views.CalibrationFailed.PropertyChanged`

The property change event to update the view.


The documentation for this class was generated from the following file:

- `source/CustomCalibrationLibrary/Views/CalibrationFailed.xaml.cs`


8.10 CustomCalibrationLibrary.Views.CalibrationFrame Class Reference

Interaction logic for CalibrationCollection.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.CalibrationFrame:

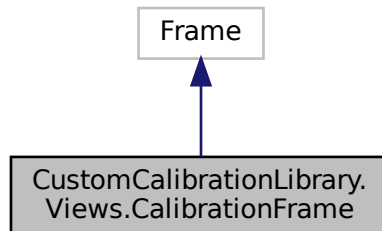

Collaboration diagram for CustomCalibrationLibrary.Views.CalibrationFrame:

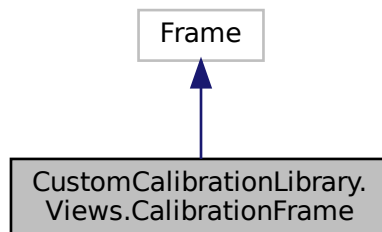

Public Member Functions

- [CalibrationFrame](#) (CalibrationModel model, Window window)
Initializes a new instance of the [CalibrationFrame](#) class.

8.10.1 Detailed Description

Interaction logic for CalibrationCollection.xaml

8.10.2 Constructor & Destructor Documentation

8.10.2.1 CalibrationFrame()

```
CustomCalibrationLibrary.Views.CalibrationFrame.CalibrationFrame (
    CalibrationModel model,
    Window window ) [inline]
```

Initializes a new instance of the [CalibrationFrame](#) class.

Parameters

<i>model</i>	The calibration model.
<i>window</i>	The target window.


The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Views/CalibrationFrame.xaml.cs


8.11 CustomCalibrationLibrary.Models.CalibrationModel Class Reference

The model for the calibration process.

Inheritance diagram for CustomCalibrationLibrary.Models.CalibrationModel:

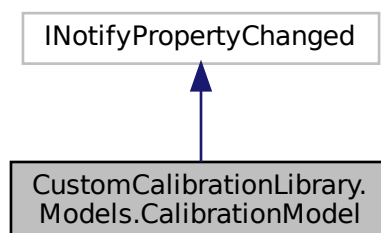

Collaboration diagram for CustomCalibrationLibrary.Models.CalibrationModel:

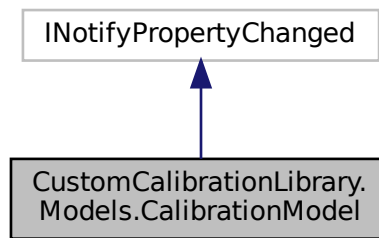

Public Member Functions

- void `OnCalibrationEvent` (`CalibrationEventType` type)
The calibraion event change handler.
- `CalibrationModel` (`TrackerLogger` logger, `double[][]` points)
Initializes a new instance of the `CalibrationModel` class.
- void `UpdateGazePoint` (`double` x, `double` y)
Update the normalized gaze point on the screen.
- void `InitCalibration` ()
Initialise the calibration.
- void `NextCalibrationPoint` ()
Trigger the next calibration point.
- void `RedoCalibrationPoint` ()
Remove and re-add the current calibration point
- void `GazeDataCollected` ()
Trigger the data collected events.
- void `SetCalibrationResult` (`List`< `GazeCalibrationData` > points)
Updates the calibration results on the screen.

Properties

- string `Error` [get, set]
The error message of the calibration process.
- `GazeValidationData ValidationData` [get, set]
The data returned by a successful validation process.
- Cursor `CursorType` [get, set]
The data returned by a successful validation process.
- `CalibrationStatus Status` [get, set]
The status of the calibarion process.
- `CalibrationStatus LastStatus` [get]
The calibration status before an error occured.
- `Point[] Points` [get]
All calibration points.

- `ObservableCollection< CalibrationPoint > CalibrationPoints` [get]
The calibration points to be added during the calibration process.
- `Point GazePoint` [get]
The gaze point position.
- `UserPositionData UserPositionGuide` [get, set]
The user position giude values.
- `int Index` [get]
The index of the current calibration point

Events

- `EventHandler< CalibrationEventType >? CalibrationEvent`
Event to trigger changes in the calibration process.
- `PropertyChangedEventHandler? PropertyChanged`
Event to trigger property changes in this class.
- `EventHandler< Point >? GazePointChanged`
Event to trigger gaze point changes.
- `EventHandler< UserPositionData >? UserPositionGuideChanged`
Event to trigger user position guide changes.

8.11.1 Detailed Description

The model for the calibration process.

8.11.2 Constructor & Destructor Documentation

8.11.2.1 CalibrationModel()

```
CustomCalibrationLibrary.Models.CalibrationModel.CalibrationModel (
    TrackerLogger logger,
    double points[][] ) [inline]
```

Initializes a new instance of the [CalibrationModel](#) class.

Parameters

<i>logger</i>	The log handler.
<i>points</i>	Calibration points.

8.11.3 Member Function Documentation

8.11.3.1 GazeDataCollected()

```
void CustomCalibrationLibrary.Models.CalibrationModel.GazeDataCollected ( ) [inline]
```

Trigger the data collected events.

8.11.3.2 InitCalibration()

```
void CustomCalibrationLibrary.Models.CalibrationModel.InitCalibration ( ) [inline]
```

Initialise the calibration.

8.11.3.3 NextCalibrationPoint()

```
void CustomCalibrationLibrary.Models.CalibrationModel.NextCalibrationPoint ( ) [inline]
```

Trigger the next calibration point.

8.11.3.4 OnCalibrationEvent()

```
void CustomCalibrationLibrary.Models.CalibrationModel.OnCalibrationEvent (
    CalibrationEventType type ) [inline]
```

The calibraion event change handler.

Parameters

<i>type</i>	
-------------	--

8.11.3.5 RedoCalibrationPoint()

```
void CustomCalibrationLibrary.Models.CalibrationModel.RedoCalibrationPoint ( ) [inline]
```

Remove and re-add the current calibration point

8.11.3.6 SetCalibrationResult()

```
void CustomCalibrationLibrary.Models.CalibrationModel.SetCalibrationResult (
    List< GazeCalibrationData > points ) [inline]
```

Updates the calibration results on the screen.

Parameters

<i>points</i>	
---------------	--

8.11.3.7 UpdateGazePoint()

```
void CustomCalibrationLibrary.Models.CalibrationModel.UpdateGazePoint (
    double x,
    double y ) [inline]
```

Update the normalized gaze point on the screen.

Parameters

<i>x</i>	The x coordinate
<i>y</i>	The y coordinate

8.11.4 Property Documentation

8.11.4.1 CalibrationPoints

```
ObservableCollection<CalibrationPoint> CustomCalibrationLibrary.Models.CalibrationModel.↔
CalibrationPoints [get]
```

The calibration points to be added during the calibration process.

8.11.4.2 CursorType

```
Cursor CustomCalibrationLibrary.Models.CalibrationModel.CursorType [get], [set]
```

The data returned by a successful validation process.

8.11.4.3 Error

```
string CustomCalibrationLibrary.Models.CalibrationModel.Error [get], [set]
```

The error message of the calibration process.

8.11.4.4 GazePoint

`Point CustomCalibrationLibrary.Models.CalibrationModel.GazePoint [get]`

The gaze point position.

8.11.4.5 Index

`int CustomCalibrationLibrary.Models.CalibrationModel.Index [get]`

The index of the current calibration point

8.11.4.6 LastStatus

`CalibrationStatus CustomCalibrationLibrary.Models.CalibrationModel.LastStatus [get]`

The calibration status before an error occurred.

8.11.4.7 Points

`Point [] CustomCalibrationLibrary.Models.CalibrationModel.Points [get]`

All calibration points.

8.11.4.8 Status

`CalibrationStatus CustomCalibrationLibrary.Models.CalibrationModel.Status [get], [set]`

The status of the calibration process.

8.11.4.9 UserPositionGuide

`UserPositionData CustomCalibrationLibrary.Models.CalibrationModel.UserPositionGuide [get], [set]`

The user position guide values.

8.11.4.10 ValidationData

```
GazeValidationData CustomCalibrationLibrary.Models.CalibrationModel.ValidationData [get],  
[set]
```

The data returned by a successful validation process.

8.11.5 Event Documentation

8.11.5.1 CalibrationEvent

```
EventHandler<CalibrationEventType>? CustomCalibrationLibrary.Models.CalibrationModel.Calibration↔  
Event
```

Event to trigger changes in the calibration process.

8.11.5.2 GazePointChanged

```
EventHandler<Point>? CustomCalibrationLibrary.Models.CalibrationModel.GazePointChanged
```

Event to trigger gaze point changes.

8.11.5.3 PropertyChanged

```
PropertyChangedEventHandler? CustomCalibrationLibrary.Models.CalibrationModel.PropertyChanged
```

Event to trigger property changes in this class.

8.11.5.4 UserPositionGuideChanged

```
EventHandler<UserPositionData>? CustomCalibrationLibrary.Models.CalibrationModel.UserPosition↔  
GuideChanged
```

Event to trigger user position guide changes.


The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Models/CalibrationModel.cs


8.12 GazeUtilityLibrary.DataStructs.CalibrationPoint Class Reference

A calibration point class holding several metrics connected to a calibration point.

Inheritance diagram for GazeUtilityLibrary.DataStructs.CalibrationPoint:

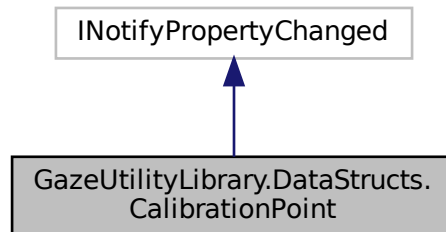

Collaboration diagram for GazeUtilityLibrary.DataStructs.CalibrationPoint:

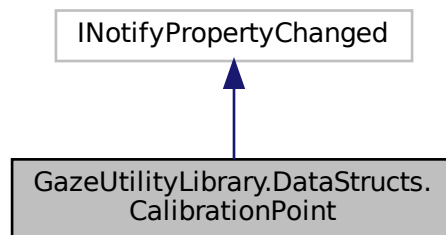

Public Member Functions

- [CalibrationPoint](#) (Point position, int index)
Initializes a new instance of the [CalibrationPoint](#) class.

Properties

- int [Index](#) [get]
The index of the calibration point.
- bool [HasData](#) [get, set]
Flag to indicate whether data has been collected for this calibration point.
- Point [Position](#) [get, set]
The position of the calibration point.

- Point [GazePositionAverage](#) [get, set]
The average between the left and the right gaze point.
- Point [GazePositionLeft](#) [get, set]
The left gaze point.
- Point [GazePositionRight](#) [get, set]
The right gaze point.

Events

- PropertyChangedEventHandler? [PropertyChanged](#)
Event to trigger property changes.

8.12.1 Detailed Description

A calibration point class holding several metrics connected to a calibration point.

8.12.2 Constructor & Destructor Documentation

8.12.2.1 CalibrationPoint()

```
GazeUtilityLibrary.DataStructs.CalibrationPoint.CalibrationPoint (
    Point position,
    int index ) [inline]
```

Initializes a new instance of the [CalibrationPoint](#) class.

Parameters

<i>position</i>	The position of the calibration point.
<i>index</i>	The index of the calibration point.

8.12.3 Property Documentation

8.12.3.1 GazePositionAverage

```
Point GazeUtilityLibrary.DataStructs.CalibrationPoint.GazePositionAverage [get], [set]
```

The average between the left and the right gaze point.

8.12.3.2 GazePositionLeft

```
Point GazeUtilityLibrary.DataStructs.CalibrationPoint.GazePositionLeft [get], [set]
```

The left gaze point.

8.12.3.3 GazePositionRight

```
Point GazeUtilityLibrary.DataStructs.CalibrationPoint.GazePositionRight [get], [set]
```

The right gaze point.

8.12.3.4 HasData

```
bool GazeUtilityLibrary.DataStructs.CalibrationPoint.HasData [get], [set]
```

Flag to indicate whether data has been collected for this calibration point.

8.12.3.5 Index

```
int GazeUtilityLibrary.DataStructs.CalibrationPoint.Index [get]
```

The index of the calibration point.

8.12.3.6 Position

```
Point GazeUtilityLibrary.DataStructs.CalibrationPoint.Position [get], [set]
```

The position of the calibration point.

8.12.4 Event Documentation

8.12.4.1 PropertyChanged

```
PropertyChangedEventHandler? GazeUtilityLibrary.DataStructs.CalibrationPoint.PropertyChanged
```

Event to trigger property changes.


The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/DataStructs/CalibrationPoint.cs`


8.13 CustomCalibrationLibrary.Views.CalibrationPoint Class Reference

Interaction logic for CalibrationPoint.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.CalibrationPoint:

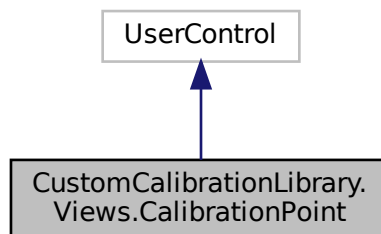

Collaboration diagram for CustomCalibrationLibrary.Views.CalibrationPoint:

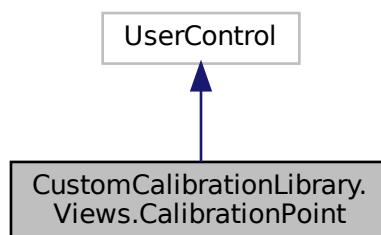

Public Member Functions

- [CalibrationPoint](#) ()
Initializes a new instance of the [CalibrationPoint](#) class.

8.13.1 Detailed Description

Interaction logic for CalibrationPoint.xaml

8.13.2 Constructor & Destructor Documentation

8.13.2.1 CalibrationPoint()

```
CustomCalibrationLibrary.Views.CalibrationPoint.CalibrationPoint ( ) [inline]
```

Initializes a new instance of the [CalibrationPoint](#) class.


The documentation for this class was generated from the following file:

- `source/CustomCalibrationLibrary/Views/CalibrationPoint.xaml.cs`


8.14 CustomCalibrationLibrary.ViewModels.CalibrationPointViewModel Class Reference

The view model for a calibration point.

Inheritance diagram for CustomCalibrationLibrary.ViewModels.CalibrationPointViewModel:

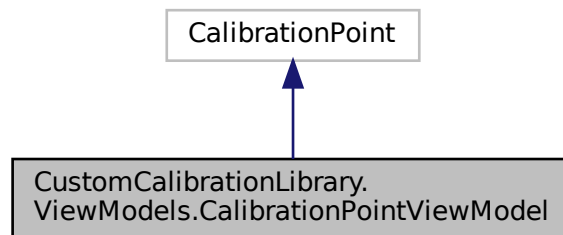

Collaboration diagram for CustomCalibrationLibrary.ViewModels.CalibrationPointViewModel:

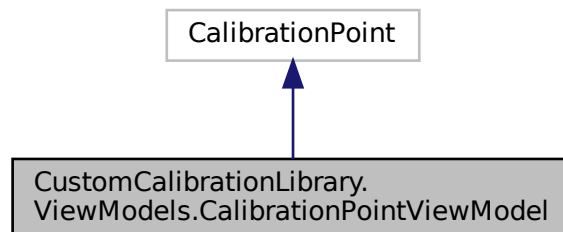

Public Member Functions

- [CalibrationPointViewModel](#) (Point point, int index)
Initializes a new instance of the [CalibrationPointViewModel](#) class.
- [CalibrationPointViewModel](#) (CalibrationPoint point)
Initializes a new instance of the [CalibrationPointViewModel](#) class.

8.14.1 Detailed Description

The view model for a calibration point.

8.14.2 Constructor & Destructor Documentation

8.14.2.1 CalibrationPointViewModel() [1/2]

```
CustomCalibrationLibrary.ViewModels.CalibrationPointViewModel.CalibrationPointViewModel (
    Point point,
    int index ) [inline]
```

Initializes a new instance of the [CalibrationPointViewModel](#) class.

Parameters

<i>point</i>	The position of the calibration point.
<i>index</i>	The index of the calibration point.

8.14.2.2 CalibrationPointViewModel() [2/2]

```
CustomCalibrationLibrary.ViewModels.CalibrationPointViewModel.CalibrationPointViewModel (
    CalibrationPoint point ) [inline]
```

Initializes a new instance of the [CalibrationPointViewModel](#) class.

Parameters

<i>point</i>	The calibration point object.
--------------	-------------------------------


The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/ViewModels/CalibrationPointViewModel.cs


8.15 CustomCalibrationLibrary.Views.CalibrationResult Class Reference

Interaction logic for CalibrationResult.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.CalibrationResult:

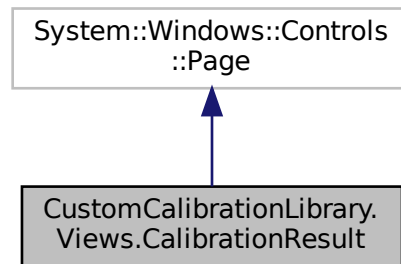

Collaboration diagram for CustomCalibrationLibrary.Views.CalibrationResult:

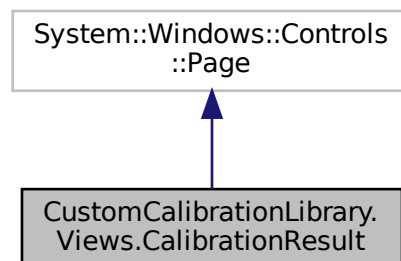

Public Member Functions

- [CalibrationResult](#) (CalibrationModel model)
Initializes a new instance of the [CalibrationResult](#) class.

8.15.1 Detailed Description

Interaction logic for CalibrationResult.xaml

8.15.2 Constructor & Destructor Documentation

8.15.2.1 CalibrationResult()

```
CustomCalibrationLibrary.Views.CalibrationResult.CalibrationResult (
    CalibrationModel model ) [inline]
```

Initializes a new instance of the [CalibrationResult](#) class.

Parameters

<i>model</i>	The calibration model.
--------------	------------------------


The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Views/CalibrationResult.xaml.cs


8.16 CustomCalibrationLibrary.Views.CalibrationResultPoint Class Reference

Interaction logic for CalibrationResultPoint.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.CalibrationResultPoint:

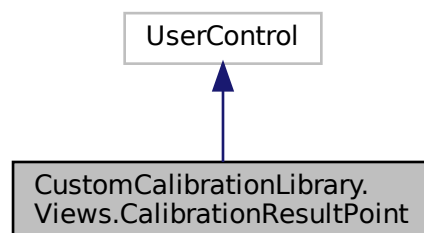

Collaboration diagram for CustomCalibrationLibrary.Views.CalibrationResultPoint:

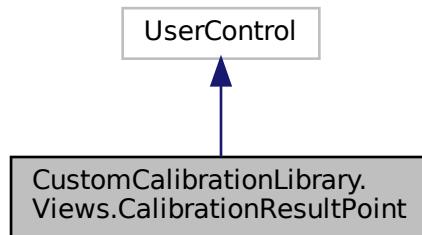

Public Member Functions

- [CalibrationResultPoint](#) ()
Initializes a new instance of the [CalibrationResultPoint](#) class.

8.16.1 Detailed Description

Interaction logic for CalibrationResultPoint.xaml

8.16.2 Constructor & Destructor Documentation

8.16.2.1 CalibrationResultPoint()

```
CustomCalibrationLibrary.Views.CalibrationResultPoint.CalibrationResultPoint ( ) [inline]
```

Initializes a new instance of the [CalibrationResultPoint](#) class.


The documentation for this class was generated from the following file:

- `source/CustomCalibrationLibrary/Views/CalibrationResultPoint.xaml.cs`


8.17 CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel Class Reference

View model class of the gaze calibration result.

Inheritance diagram for CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel:

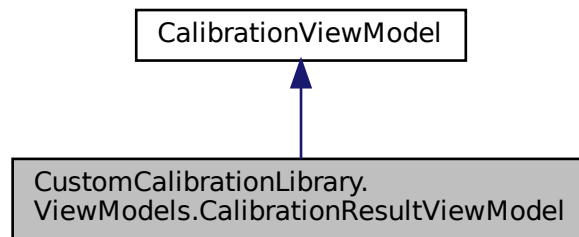

Collaboration diagram for CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel:

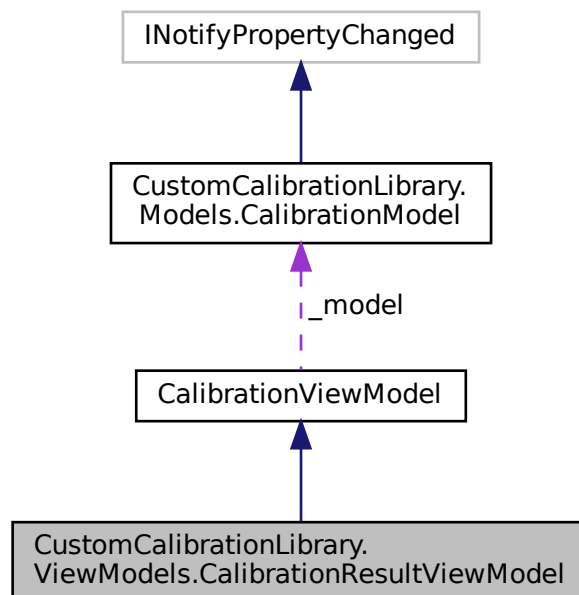

Public Member Functions

- [CalibrationResultViewModel](#) ([CalibrationModel](#) model)

Constructor

- void [OnGazeToggle](#) ()

Toggle the visibility of the live gaze point.

Properties

- ICommand [CalibrationRestartCommand](#) [get]
Command to restart the calibration
- ICommand [CalibrationAcceptCommand](#) [get]
Command to accept the calibration
- ICommand [GazeVisibilityCommand](#) [get]
Command to toggle the visibility of the live gaze point
- [LiveGazePoint](#) [GazePoint](#) [get]
The position of the live gaze point

Additional Inherited Members

8.17.1 Detailed Description

View model class of the gaze calibration result.

8.17.2 Constructor & Destructor Documentation

8.17.2.1 CalibrationResultViewModel()

```
CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel.CalibrationResultViewModel (
    CalibrationModel model ) [inline]
```

Constructor

Parameters

<i>model</i>	The claibration model
--------------	-----------------------

8.17.3 Member Function Documentation

8.17.3.1 OnGazeToggle()

```
void CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel.OnGazeToggle ( ) [inline]
```

Toggle the visibility of the live gaze point.

8.17.4 Property Documentation

8.17.4.1 CalibrationAcceptCommand

```
 ICommand CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel.CalibrationAccept↵  
 Command [get]
```

Command to accept the calibration

8.17.4.2 CalibrationRestartCommand

```
 ICommand CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel.CalibrationRestart↵  
 Command [get]
```

Command to restart the calibration

8.17.4.3 GazePoint

```
 LiveGazePoint CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel.GazePoint [get]
```

The position of the live gaze point

8.17.4.4 GazeVisibilityCommand

```
 ICommand CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel.GazeVisibilityCommand  
 [get]
```

Command to toggle the visibility of the live gaze point

The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/ViewModels/CalibrationResultViewModel.cs

8.18 Tobii.Research.Addons.CalibrationValidationPoint Class Reference

Represents a collected point that goes into the calibration validation. It contains calculated values for accuracy and precision as well as the original gaze samples collected for the point.

Public Member Functions

- override string [ToString](#) ()
Convert validation values to a string.

Properties

- NormalizedPoint2D [Coordinates](#) [get]
The 2D coordinates of this point (in Active Display Coordinate System).
- float [AccuracyLeftEye](#) [get]
The accuracy in degrees for the left eye.
- float [PrecisionLeftEye](#) [get]
The precision (standard deviation) in degrees for the left eye.
- float [PrecisionRMSLeftEye](#) [get]
The precision (root mean square of sample-to-sample error) in degrees for the left eye.
- float [AccuracyRightEye](#) [get]
The accuracy in degrees for the right eye.
- float [PrecisionRightEye](#) [get]
The precision (standard deviation) in degrees for the right eye.
- float [PrecisionRMSRightEye](#) [get]
The precision (root mean square of sample-to-sample error) in degrees for the right eye.
- bool [TimedOut](#) [get]
A boolean indicating if there was a timeout while collecting data for this point.
- GazeDataEventArgs[] [GazeData](#) [get]
The gaze data samples collected for this point. These samples are the base for the calculated accuracy and precision.

8.18.1 Detailed Description

Represents a collected point that goes into the calibration validation. It contains calculated values for accuracy and precision as well as the original gaze samples collected for the point.

8.18.2 Member Function Documentation

8.18.2.1 ToString()

```
override string Tobii.Research.Addons.CalibrationValidationPoint.ToString ( ) [inline]
```

Convert validation values to a string.

Returns

The validation string.

8.18.3 Property Documentation

8.18.3.1 AccuracyLeftEye

```
float Tobii.Research.Addons.CalibrationValidationPoint.AccuracyLeftEye [get]
```

The accuracy in degrees for the left eye.

8.18.3.2 AccuracyRightEye

```
float Tobii.Research.Addons.CalibrationValidationPoint.AccuracyRightEye [get]
```

The accuracy in degrees for the right eye.

8.18.3.3 Coordinates

```
NormalizedPoint2D Tobii.Research.Addons.CalibrationValidationPoint.Coordinates [get]
```

The 2D coordinates of this point (in Active Display Coordinate System).

8.18.3.4 GazeData

```
GazeDataEventArgs [] Tobii.Research.Addons.CalibrationValidationPoint.GazeData [get]
```

The gaze data samples collected for this point. These samples are the base for the calculated accuracy and precision.

8.18.3.5 PrecisionLeftEye

```
float Tobii.Research.Addons.CalibrationValidationPoint.PrecisionLeftEye [get]
```

The precision (standard deviation) in degrees for the left eye.

8.18.3.6 PrecisionRightEye

```
float Tobii.Research.Addons.CalibrationValidationPoint.PrecisionRightEye [get]
```

The precision (standard deviation) in degrees for the right eye.

8.18.3.7 PrecisionRMSLeftEye

`float Tobii.Research.Addons.CalibrationValidationPoint.PrecisionRMSLeftEye [get]`

The precision (root mean square of sample-to-sample error) in degrees for the left eye.

8.18.3.8 PrecisionRMSRightEye

`float Tobii.Research.Addons.CalibrationValidationPoint.PrecisionRMSRightEye [get]`

The precision (root mean square of sample-to-sample error) in degrees for the right eye.

8.18.3.9 TimedOut

`bool Tobii.Research.Addons.CalibrationValidationPoint.TimedOut [get]`

A boolean indicating if there was a timeout while collecting data for this point.

The documentation for this class was generated from the following file:

- `source/TobiiProSdkAddons/ScreenBasedCalibrationValidation.cs`

8.19 Tobii.Research.Addons.CalibrationValidationResult Class Reference

Contains the result of the calibration validation.

Public Member Functions

- override string [ToString](#) ()
Convert validation values to a string.

Properties

- List< [CalibrationValidationPoint](#) > [Points](#) [get]
The results of the calibration validation per point (same points as were collected).
- float [AverageAccuracyLeftEye](#) [get]
The accuracy in degrees averaged over all collected points for the left eye.
- float [AveragePrecisionLeftEye](#) [get]
The precision (standard deviation) in degrees averaged over all collected points for the left eye.
- float [AveragePrecisionRMSLeftEye](#) [get]
The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the left eye.
- float [AverageAccuracyRightEye](#) [get]
The accuracy in degrees averaged over all collected points for the right eye.
- float [AveragePrecisionRightEye](#) [get]
The precision (standard deviation) in degrees averaged over all collected points for the right eye.
- float [AveragePrecisionRMSRightEye](#) [get]
The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the right eye.

8.19.1 Detailed Description

Contains the result of the calibration validation.

8.19.2 Member Function Documentation

8.19.2.1 ToString()

```
override string Tobii.Research.Addons.CalibrationValidationResult.ToString ( ) [inline]
```

Convert validation values to a string.

Returns

The validation string.

8.19.3 Property Documentation

8.19.3.1 AverageAccuracyLeftEye

```
float Tobii.Research.Addons.CalibrationValidationResult.AverageAccuracyLeftEye [get]
```

The accuracy in degrees averaged over all collected points for the left eye.

8.19.3.2 AverageAccuracyRightEye

```
float Tobii.Research.Addons.CalibrationValidationResult.AverageAccuracyRightEye [get]
```

The accuracy in degrees averaged over all collected points for the right eye.

8.19.3.3 AveragePrecisionLeftEye

```
float Tobii.Research.Addons.CalibrationValidationResult.AveragePrecisionLeftEye [get]
```

The precision (standard deviation) in degrees averaged over all collected points for the left eye.

8.19.3.4 AveragePrecisionRightEye

```
float Tobii.Research.Addons.CalibrationValidationResult.AveragePrecisionRightEye [get]
```

The precision (standard deviation) in degrees averaged over all collected points for the right eye.

8.19.3.5 AveragePrecisionRMSLeftEye

```
float Tobii.Research.Addons.CalibrationValidationResult.AveragePrecisionRMSLeftEye [get]
```

The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the left eye.

8.19.3.6 AveragePrecisionRMSRightEye

```
float Tobii.Research.Addons.CalibrationValidationResult.AveragePrecisionRMSRightEye [get]
```

The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the right eye.

8.19.3.7 Points

```
List<CalibrationValidationPoint> Tobii.Research.Addons.CalibrationValidationResult.Points  
[get]
```

The results of the calibration validation per point (same points as were collected).


The documentation for this class was generated from the following file:

- source/TobiiProSdkAddons/ScreenBasedCalibrationValidation.cs


8.20 CustomCalibrationLibrary.ViewModels.CalibrationViewModel Class Reference

The view model class of the calibration view

Inheritance diagram for CustomCalibrationLibrary.ViewModels.CalibrationViewModel:

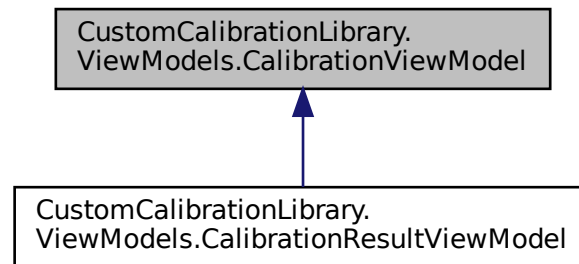

Collaboration diagram for CustomCalibrationLibrary.ViewModels.CalibrationViewModel:

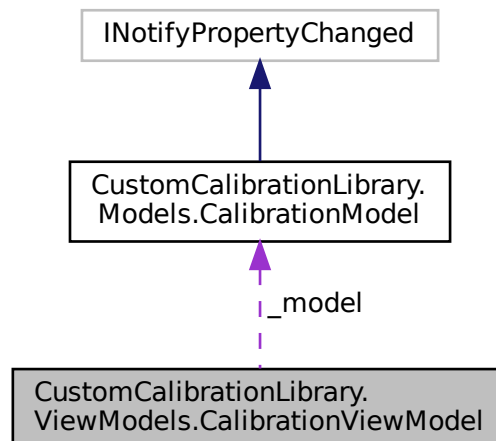

Public Member Functions

- [CalibrationViewModel](#) ([CalibrationModel](#) model)
Constructor

Protected Attributes

- [CalibrationModel _model](#)
The claibration model.

Properties

- ObservableCollection< [CalibrationPointViewModel](#) > [CalibrationPoints](#) [get]
The collection of calibration points to be shown on the view
- Cursor **CursorType** [get]

8.20.1 Detailed Description

The view model class of the calibration view

8.20.2 Constructor & Destructor Documentation

8.20.2.1 CalibrationViewModel()

```
CustomCalibrationLibrary.ViewModels.CalibrationViewModel.CalibrationViewModel (
    CalibrationModel model ) [inline]
```

Constructor

Parameters

<i>model</i>	The calibration model
--------------	-----------------------

8.20.3 Member Data Documentation

8.20.3.1 _model

```
CalibrationModel CustomCalibrationLibrary.ViewModels.CalibrationViewModel._model [protected]
```

The claibration model.

8.20.4 Property Documentation

8.20.4.1 CalibrationPoints

```
ObservableCollection<CalibrationPointViewModel> CustomCalibrationLibrary.ViewModels.CalibrationWindowViewModel.CalibrationPoints [get]
```

The collection of calibration points to be shown on the view


The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/ViewModels/CalibrationViewModel.cs


8.21 CustomCalibrationLibrary.Views.CalibrationWindow Class Reference

Interaction logic for MainWindow.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.CalibrationWindow:

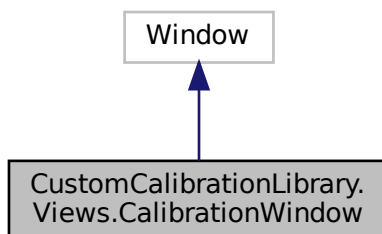

Collaboration diagram for CustomCalibrationLibrary.Views.CalibrationWindow:

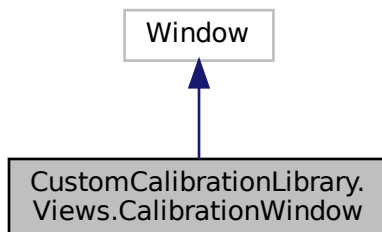

8.21.1 Detailed Description

Interaction logic for MainWindow.xaml


The documentation for this class was generated from the following file:

- `source/CustomCalibrationLibrary/Views/CalibrationWindow.xaml.cs`


8.22 CustomCalibrationLibrary.Views.Computing Class Reference

Interaction logic for Computing.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.Computing:

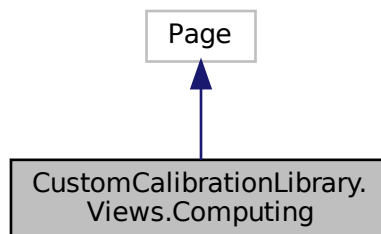

Collaboration diagram for CustomCalibrationLibrary.Views.Computing:

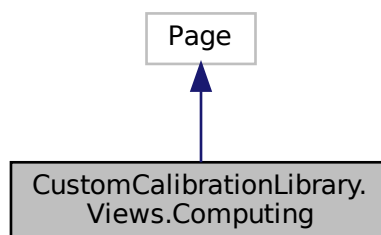

Public Member Functions

- [Computing](#) ()

Initializes a new instance of the [Computing](#) class.

8.22.1 Detailed Description

Interaction logic for Computing.xaml

8.22.2 Constructor & Destructor Documentation

8.22.2.1 Computing()

```
CustomCalibrationLibrary.Views.Computing.Computing ( ) [inline]
```

Initializes a new instance of the [Computing](#) class.

The documentation for this class was generated from the following file:

- `source/CustomCalibrationLibrary/Views/Computing.xaml.cs`

8.23 GazeUtilityLibrary.ConfigItem Class Reference

configuration file class

Public Member Functions

- [ConfigItem](#) ()
Initializes a new instance of the [ConfigItem](#) class.

Properties

- string [CalibrationLogColumnOrder](#) [get, set]
Allows to define the order and the delimiters between the different calibration data values.
- string[] [CalibrationLogColumnName](#) [get, set]
Defines the titles of the calibration data log value columns.
- bool [CalibrationLogWriteOutput](#) [get, set]
Defines whether gaze calibration data is written to a log file.
- double[][] [CalibrationPoints](#) [get, set]
Define the calibration points to be shown during the calibration process.
- double [DriftCompensationDispersionThreshold](#) [get, set]
In order to detect a fixation with the I-DT algorithm a dispersion threshold is required. Provide an angle in degrees.
- double [DriftCompensationDispersionThresholdMax](#) [get, set]
In order to prevent drift compensation from getting out of hand limit the maximal allowed dispersion. If the drift compensation angle is larger than the here defined degrees, no compensation is applied. Provide an angle in degrees.
- int [DriftCompensationDurationThreshold](#) [get, set]
Specifies the amount of time (in milliseconds) required to fixate the target during drift compensation.
- int [DriftCompensationTimer](#) [get, set]
Specifies the amount of time (in milliseconds) to wait for a fixation point during drift compensation.

- bool [DriftCompensationWindowShow](#) [get, set]
If set to true the drift compensation window is shown on the drift compensation command. Otherwise only the drift compensation process is done without showing the window.
- int [ValidationDurationThreshold](#) [get, set]
Specifies the amount of time (in milliseconds) required to fixate the target during validation.
- string [ValidationLogColumnOrder](#) [get, set]
Allows to define the order and the delimiters between the different validation data values.
- string[] [ValidationLogColumnNameTitle](#) [get, set]
Defines the titles of the validation data log value columns.
- bool [ValidationLogWriteOutput](#) [get, set]
Defines whether gaze validation data is written to a log file.
- double[][] [ValidationPoints](#) [get, set]
Define the validation points to be shown during the validation process.
- int [ValidationTimer](#) [get, set]
Specifies the amount of time (in milliseconds) to wait for a fixation point during validation.
- string [DataLogColumnOrder](#) [get, set]
Allows to define the order and the delimiters between the different gaze data values.
- string[] [DataLogColumnNameTitle](#) [get, set]
Defines the titles of the gaze data log value columns.
- int [DataLogCount](#) [get, set]
Number of maximal allowed output data files in the output path. Oldest files are deleted first.
- string [DataLogFormatDiameter](#) [get, set]
Allows to define the format of how the pupil diameter (in millimetres) will be logged.
- string [DataLogFormatOrigin](#) [get, set]
Allows to define the format of how the gaze origin values (in millimetres) will be logged.
- string [DataLogFormatNormalizedPoint](#) [get, set]
Allows to define the format of how normalized data points will be logged.
- string [DataLogFormatTimeStamp](#) [get, set]
Allows to define the format of the timestamp.
- string [DataLogFormatTimeStampRelative](#) [get, set]
Allows to define the format of the relative timestamp in milliseconds.
- string [DataLogFormatValidation](#) [get, set]
Allows to define the format of the validation values.
- string [DataLogPath](#) [get, set]
Defines the location of the output file. It must be the path to a folder (not a file).
- bool [DataLogWriteOutput](#) [get, set]
Defines whether gaze data is written to a log file.
- bool [DataLogDisabledOnStartup](#) [get, set]
Defines whether gaze data storing is disabled on Gaze application start.
- string? [ConfigName](#) [get, set]
The name of the experiment.
- string? [LicensePath](#) [get, set]
Defines the location of the license files. It must be the path to a folder (not a file).
- int [ReadyTimer](#) [get, set]
Specifies the amount of time (in milliseconds) to wait for the eye tracker to become ready while it is in any other state.
- int [TrackerDevice](#) [get, set]
Choose the tracker device (1: [Tobii Pro SDK](#), 2: [Mouse Tracker](#)).
- bool [MouseControl](#) [get, set]
Defines whether the mouse cursor shall be controlled by the gaze of the subject during the experiment.
- bool [MouseControlHide](#) [get, set]
Defines whether the mouse cursor shall be hidden during the experiment.

- bool [MouseCalibrationHide](#) [get, set]
Defines whether the mouse cursor shall be hidden on the calibration window.
- string [MouseStandardIconPath](#) [get, set]
Defines the Path to the standard mouse pointer icon.
- string [TobiiApplicationPath](#) [get, set]
*Defines the *Tobii* installation path. It must be the path to a folder (not a file).*
- string [TobiiCalibrate](#) [get, set]
*The *Tobii* application to run a calibration.*
- string [TobiiCalibrateArguments](#) [get, set]
The arguments to pass to the calibration application. Use S as a placeholder for the device serial number and A as a placeholder for the device address.
- [ConfigScreenArea](#) [ScreenArea](#) [get, set]
Hold the screen area once the config file is dumped during experimentation.

8.23.1 Detailed Description

configuration file class

8.23.2 Constructor & Destructor Documentation

8.23.2.1 ConfigItem()

```
GazeUtilityLibrary.ConfigItem.ConfigItem ( ) [inline]
```

Initializes a new instance of the [ConfigItem](#) class.

8.23.3 Property Documentation

8.23.3.1 CalibrationLogColumnOrder

```
string GazeUtilityLibrary.ConfigItem.CalibrationLogColumnOrder [get], [set]
```

Allows to define the order and the delimiters between the different calibration data values.

8.23.3.2 CalibrationLogColumnNameTitle

```
string [] GazeUtilityLibrary.ConfigItem.CalibrationLogColumnNameTitle [get], [set]
```

Defines the titles of the calibration data log value columns.

8.23.3.3 CalibrationLogWriteOutput

```
bool GazeUtilityLibrary.ConfigItem.CalibrationLogWriteOutput [get], [set]
```

Defines whether gaze calibration data is written to a log file.

8.23.3.4 CalibrationPoints

```
double [][] GazeUtilityLibrary.ConfigItem.CalibrationPoints [get], [set]
```

Define the calibration points to be shown during the calibration process.

8.23.3.5 ConfigName

```
string? GazeUtilityLibrary.ConfigItem.ConfigName [get], [set]
```

The name of the experiment.

8.23.3.6 DataLogColumnOrder

```
string GazeUtilityLibrary.ConfigItem.DataLogColumnOrder [get], [set]
```

Allows to define the order and the delimiters between the different gaze data values.

8.23.3.7 DataLogColumnNameTitle

```
string [] GazeUtilityLibrary.ConfigItem.DataLogColumnNameTitle [get], [set]
```

Defines the titles of the gaze data log value columns.

8.23.3.8 DataLogCount

```
int GazeUtilityLibrary.ConfigItem.DataLogCount [get], [set]
```

Number of maximal allowed output data files in the output path. Oldest files are deleted first.

8.23.3.9 DataLogDisabledOnStartup

```
bool GazeUtilityLibrary.ConfigItem.DataLogDisabledOnStartup [get], [set]
```

Defines whether gaze data storing is disabled on Gaze application start.

8.23.3.10 DataLogFormatDiameter

```
string GazeUtilityLibrary.ConfigItem.DataLogFormatDiameter [get], [set]
```

Allows to define the format of how the pupil diameter (in millimetres) will be logged.

8.23.3.11 DataLogFormatNormalizedPoint

```
string GazeUtilityLibrary.ConfigItem.DataLogFormatNormalizedPoint [get], [set]
```

Allows to define the format of how normalized data points will be logged.

8.23.3.12 DataLogFormatOrigin

```
string GazeUtilityLibrary.ConfigItem.DataLogFormatOrigin [get], [set]
```

Allows to define the format of how the gaze origin values (in millimetres) will be logged.

8.23.3.13 DataLogFormatTimeStamp

```
string GazeUtilityLibrary.ConfigItem.DataLogFormatTimeStamp [get], [set]
```

Allows to define the format of the timestamp.

8.23.3.14 DataLogFormatTimeStampRelative

```
string GazeUtilityLibrary.ConfigItem.DataLogFormatTimeStampRelative [get], [set]
```

Allows to define the format of the relative timestamp in milliseconds.

8.23.3.15 DataLogFormatValidation

```
string GazeUtilityLibrary.ConfigItem.DataLogFormatValidation [get], [set]
```

Allows to define the format of the validation values.

8.23.3.16 DataLogPath

```
string GazeUtilityLibrary.ConfigItem.DataLogPath [get], [set]
```

Defines the location of the output file. It must be the path to a folder (not a file).

8.23.3.17 DataLogWriteOutput

```
bool GazeUtilityLibrary.ConfigItem.DataLogWriteOutput [get], [set]
```

Defines whether gaze data is written to a log file.

8.23.3.18 DriftCompensationDispersionThreshold

```
double GazeUtilityLibrary.ConfigItem.DriftCompensationDispersionThreshold [get], [set]
```

In order to detect a fixation with the I-DT algorithm a dispersion threshold is required. Provide an angle in degrees.

8.23.3.19 DriftCompensationDispersionThresholdMax

```
double GazeUtilityLibrary.ConfigItem.DriftCompensationDispersionThresholdMax [get], [set]
```

In order to prevent drift compensation from getting out of hand limit the maximal allowed dispersion. If the drift compensation angle is larger than the here defined degrees, no compensation is applied. Provide an angle in degrees.

8.23.3.20 DriftCompensationDurationThreshold

```
int GazeUtilityLibrary.ConfigItem.DriftCompensationDurationThreshold [get], [set]
```

Specifies the amount of time (in milliseconds) required to fixate the target during drift compensation.

8.23.3.21 DriftCompensationTimer

```
int GazeUtilityLibrary.ConfigItem.DriftCompensationTimer [get], [set]
```

Specifies the amount of time (in milliseconds) to wait for a fixation point during drift compensation.

8.23.3.22 DriftCompensationWindowShow

```
bool GazeUtilityLibrary.ConfigItem.DriftCompensationWindowShow [get], [set]
```

If set to true the drift compensation window is shown on the drift compensation command. Otherwise only the drift compensation process is done without showing the window.

8.23.3.23 LicensePath

```
string? GazeUtilityLibrary.ConfigItem.LicensePath [get], [set]
```

Defines the location of the license files. It must be the path to a folder (not a file).

8.23.3.24 MouseCalibrationHide

```
bool GazeUtilityLibrary.ConfigItem.MouseCalibrationHide [get], [set]
```

Defines whether the mouse cursor shall be hidden on the calibration window.

8.23.3.25 MouseControl

```
bool GazeUtilityLibrary.ConfigItem.MouseControl [get], [set]
```

Defines whether the mouse cursor shall be controlled by the gaze of the subject during the experiment.

8.23.3.26 MouseControlHide

```
bool GazeUtilityLibrary.ConfigItem.MouseControlHide [get], [set]
```

Defines whether the mouse cursor shall be hidden during the experiment.

8.23.3.27 MouseStandardIconPath

```
string GazeUtilityLibrary.ConfigItem.MouseStandardIconPath [get], [set]
```

Defines the Path to the standard mouse pointer icon.

8.23.3.28 ReadyTimer

```
int GazeUtilityLibrary.ConfigItem.ReadyTimer [get], [set]
```

Specifies the amount of time (in milliseconds) to wait for the eye tracker to become ready while it is in any other state.

8.23.3.29 ScreenArea

```
ConfigScreenArea GazeUtilityLibrary.ConfigItem.ScreenArea [get], [set]
```

Hold the screen area once the config file is dumped during experimentation.

8.23.3.30 TobiiApplicationPath

```
string GazeUtilityLibrary.ConfigItem.TobiiApplicationPath [get], [set]
```

Defines the [Tobii](#) installation path. It must be the path to a folder (not a file).

8.23.3.31 TobiiCalibrate

```
string GazeUtilityLibrary.ConfigItem.TobiiCalibrate [get], [set]
```

The [Tobii](#) application to run a calibration.

8.23.3.32 TobiiCalibrateArguments

```
string GazeUtilityLibrary.ConfigItem.TobiiCalibrateArguments [get], [set]
```

The arguments to pass to the calibration application. Use S as a placeholder for the device serial number and A as a placeholder for the device address.

8.23.3.33 TrackerDevice

```
int GazeUtilityLibrary.ConfigItem.TrackerDevice [get], [set]
```

Choose the tracker device (1: [Tobii Pro SDK](#), 2: Mouse [Tracker](#)).

8.23.3.34 ValidationDurationThreshold

```
int GazeUtilityLibrary.ConfigItem.ValidationDurationThreshold [get], [set]
```

Specifies the amount of time (in milliseconds) required to fixate the target during validation.

8.23.3.35 ValidationLogColumnOrder

```
string GazeUtilityLibrary.ConfigItem.ValidationLogColumnOrder [get], [set]
```

Allows to define the order and the delimiters between the different validation data values.

8.23.3.36 ValidationLogColumnNameTitle

```
string [] GazeUtilityLibrary.ConfigItem.ValidationLogColumnNameTitle [get], [set]
```

Defines the titles of the validation data log value columns.

8.23.3.37 ValidationLogWriteOutput

```
bool GazeUtilityLibrary.ConfigItem.ValidationLogWriteOutput [get], [set]
```

Defines whether gaze validation data is written to a log file.

8.23.3.38 ValidationPoints

```
double [][] GazeUtilityLibrary.ConfigItem.ValidationPoints [get], [set]
```

Define the validation points to be shown during the validation process.

8.23.3.39 ValidationTimer

```
int GazeUtilityLibrary.ConfigItem.ValidationTimer [get], [set]
```

Specifies the amount of time (in milliseconds) to wait for a fixation point during validation.

The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/GazeConfiguration.cs

8.24 GazeUtilityLibrary.ConfigScreenArea Class Reference

The JSON structure of the screen area.

Public Member Functions

- [ConfigScreenArea](#) ()
Initializes a new instance of the [ConfigScreenArea](#) class.
- [ConfigScreenArea](#) ([ScreenArea](#) screenArea)
Initializes a new instance of the [ConfigScreenArea](#) class.

Properties

- double [Width](#) [get, set]
The width of the screen.
- double [Height](#) [get, set]
The height of the screen.
- double[] [Center](#) [get, set]
The coordinates of the center point of the screen.
- double[] [TopLeft](#) [get, set]
The coordinates of the top left point of the screen.
- double[] [TopRight](#) [get, set]
The coordinates of the to right point of the screen.
- double[] [BottomLeft](#) [get, set]
The coordinates of the bottom left point of the screen.
- double[] [BottomRight](#) [get, set]
The coordinates of the bottom right point of the screen.

8.24.1 Detailed Description

The JSON structure of the screen area.

8.24.2 Constructor & Destructor Documentation

8.24.2.1 ConfigScreenArea() [1/2]

```
GazeUtilityLibrary.ConfigScreenArea.ConfigScreenArea ( ) [inline]
```

Initializes a new instance of the [ConfigScreenArea](#) class.

8.24.2.2 ConfigScreenArea() [2/2]

```
GazeUtilityLibrary.ConfigScreenArea.ConfigScreenArea (
    ScreenArea screenArea ) [inline]
```

Initializes a new instance of the [ConfigScreenArea](#) class.

Parameters

<i>screenArea</i>	A screen area object.
-------------------	-----------------------

8.24.3 Property Documentation

8.24.3.1 BottomLeft

```
double [ ] GazeUtilityLibrary.ConfigScreenArea.BottomLeft [get], [set]
```

The coordinates of the bottom left point of the screen.

8.24.3.2 BottomRight

```
double [ ] GazeUtilityLibrary.ConfigScreenArea.BottomRight [get], [set]
```

The coordinates of the bottom right point of the screen.

8.24.3.3 Center

```
double [ ] GazeUtilityLibrary.ConfigScreenArea.Center [get], [set]
```

The coordinates of the center point of the screen.

8.24.3.4 Height

```
double GazeUtilityLibrary.ConfigScreenArea.Height [get], [set]
```

The height of the screen.

8.24.3.5 TopLeft

```
double [] GazeUtilityLibrary.ConfigScreenArea.TopLeft [get], [set]
```

The coordinates of the top left point of the screen.

8.24.3.6 TopRight

```
double [] GazeUtilityLibrary.ConfigScreenArea.TopRight [get], [set]
```

The coordinates of the to right point of the screen.

8.24.3.7 Width

```
double GazeUtilityLibrary.ConfigScreenArea.Width [get], [set]
```

The width of the screen.


The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/GazeConfiguration.cs`


8.25 CustomCalibrationLibrary.Views.Disconnect Class Reference

Interaction logic for Disconnect.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.Disconnect:

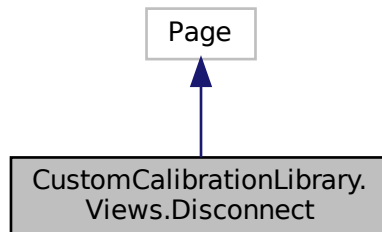

Collaboration diagram for CustomCalibrationLibrary.Views.Disconnect:

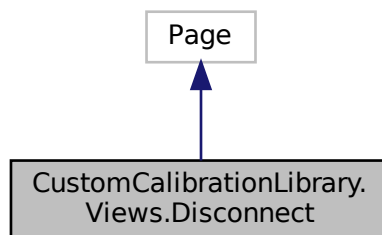

Public Member Functions

- [Disconnect](#) ([CalibrationModel](#) model)
Initializes a new instance of the [Disconnect](#) class.

Properties

- ICommand [CalibrationAbortCommand](#) [get]
Command to abort the calibration

8.25.1 Detailed Description

Interaction logic for Disconnect.xaml

8.25.2 Constructor & Destructor Documentation

8.25.2.1 Disconnect()

```
CustomCalibrationLibrary.Views.Disconnect.Disconnect (
    CalibrationModel model ) [inline]
```

Initializes a new instance of the [Disconnect](#) class.

Parameters

<i>model</i>	The calibration model
--------------	-----------------------

8.25.3 Property Documentation

8.25.3.1 CalibrationAbortCommand

```
ICommand CustomCalibrationLibrary.Views.Disconnect.CalibrationAbortCommand [get]
```

Command to abort the calibration

The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Views/Disconnect.xaml.cs

8.26 GazeUtilityLibrary.DriftCompensation Class Reference

The class to handle drift compensation.

Public Member Functions

- [DriftCompensation](#) (Vector3 fixationPoint, int fixationFrameCount, double dispersionThreshold, double dispersionThresholdMax)
Initializes a new instance of the [DriftCompensation](#) class.
- void [Reset](#) ()
Reset the drift compensation quaternion to the identity.
- void [Start](#) ()
Start the drift compensation.
- bool [Update](#) (GazeData gazeData)
Collect gaze data samples of a fixation and once enough samples are collected, compute the drift compensation quaternion.

Properties

- Quaternion [Q](#) [get]
The drift compensation quatrenion.
- double [DeviationAngle](#) [get]
The deviation angle of the drift compensation.
- double [Dispersion](#) [get]
The dispersion of the drift compensation fixation.

8.26.1 Detailed Description

The class to handle drift compensation.

8.26.2 Constructor & Destructor Documentation

8.26.2.1 DriftCompensation()

```
GazeUtilityLibrary.DriftCompensation.DriftCompensation (
    Vector3 fixationPoint,
    int fixationFrameCount,
    double dispersionThreshold,
    double dispersionThresholdMax ) [inline]
```

Initializes a new instance of the [DriftCompensation](#) class.

Parameters

<i>fixationPoint</i>	The target fixation point.
<i>fixationFrameCount</i>	The required number of frames during fixation.
<i>dispersionThreshold</i>	The dispersion threshold for the fixation.
<i>dispersionThresholdMax</i>	The maximal allowed deviation angle.

8.26.3 Member Function Documentation

8.26.3.1 Reset()

```
void GazeUtilityLibrary.DriftCompensation.Reset ( ) [inline]
```

Reset the drift compensation quaternion to the identity.

8.26.3.2 Start()

```
void GazeUtilityLibrary.DriftCompensation.Start ( ) [inline]
```

Start the drift compensation.

8.26.3.3 Update()

```
bool GazeUtilityLibrary.DriftCompensation.Update (
    GazeData gazeData ) [inline]
```

Collect gaze data samples of a fixation and once enough samples are collected, compute the drift compensation quaternion.

Parameters

<i>gazeData</i>	The gaze data sample to collect if it belongs to a fixation.
-----------------	--

Returns

True if new drift compensation is computed, false if the process is ongoing.

8.26.4 Property Documentation

8.26.4.1 DeviationAngle

```
double GazeUtilityLibrary.DriftCompensation.DeviationAngle [get]
```

The deviation angle of the drift compensation.

8.26.4.2 Dispersion

```
double GazeUtilityLibrary.DriftCompensation.Dispersion [get]
```

The dispersion of the drift compensation fixation.

8.26.4.3 Q

Quaternion GazeUtilityLibrary.DriftCompensation.Q [get]

The drift compensation quatrenion.

The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/DriftCompensation.cs

8.27 GazeUtilityLibrary.DataStructs.DriftCompensationData Class Reference

The drift compensation data structure

Public Member Functions

- [DriftCompensationData](#) ([ScreenArea](#) screen, Quaternion driftCompensation, [GazeData3d](#) gazeData)
Constructor

Properties

- Vector2 [GazePosition2d](#) [get]
The drift compensated 2d gaze position
- Vector3 [GazePosition3d](#) [get]
The drift compensated 3d gaze position
- Quaternion [Compensation](#) [get]
The drift compensation quaternion

8.27.1 Detailed Description

The drift compensation data structure

8.27.2 Constructor & Destructor Documentation

8.27.2.1 DriftCompensationData()

```
GazeUtilityLibrary.DataStructs.DriftCompensationData.DriftCompensationData (
    ScreenArea screen,
    Quaternion driftCompensation,
    GazeData3d gazeData ) [inline]
```

Constructor

Parameters

<i>screen</i>	The screen area
<i>driftCompensation</i>	The drift compensation quaternion
<i>gazeData</i>	The 3d gaze data structure

8.27.3 Property Documentation

8.27.3.1 Compensation

`Quaternion GazeUtilityLibrary.DataStructs.DriftCompensationData.Compensation` [get]

The drift compensation quaternion

8.27.3.2 GazePosition2d

`Vector2 GazeUtilityLibrary.DataStructs.DriftCompensationData.GazePosition2d` [get]

The drift compensated 2d gaze position

8.27.3.3 GazePosition3d

`Vector3 GazeUtilityLibrary.DataStructs.DriftCompensationData.GazePosition3d` [get]

The drift compensated 3d gaze position

The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/DataStructs/DriftCompensationData.cs`

8.28 CustomCalibrationLibrary.ViewModels.DriftCompensationView↔ Model Class Reference

The view model class of the drift compensation view.

Public Member Functions

- [DriftCompensationViewModel](#) ()

Constructor

Properties

- CalibrationPoint [FixationPoint](#) [get, set]

The point on the screen which the participant is supposed to fixate.

8.28.1 Detailed Description

The view model class of the drift compensation view.

8.28.2 Constructor & Destructor Documentation

8.28.2.1 DriftCompensationViewModel()

```
CustomCalibrationLibrary.ViewModels.DriftCompensationViewModel.DriftCompensationViewModel ( )  
[inline]
```

Constructor

8.28.3 Property Documentation

8.28.3.1 FixationPoint

```
CalibrationPoint CustomCalibrationLibrary.ViewModels.DriftCompensationViewModel.FixationPoint  
[get], [set]
```

The point on the screen which the participant is supposed to fixate.


The documentation for this class was generated from the following file:

- `source/CustomCalibrationLibrary/ViewModels/DriftCompensationViewModel.cs`


8.29 CustomCalibrationLibrary.Views.DriftCompensationWindow Class Reference

Interaction logic for DriftCompensation.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.DriftCompensationWindow:

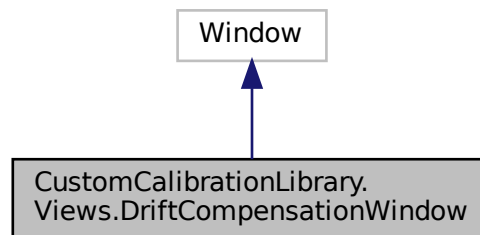

Collaboration diagram for CustomCalibrationLibrary.Views.DriftCompensationWindow:

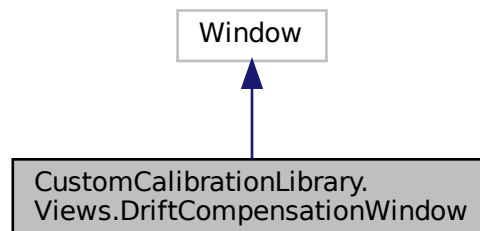

Public Member Functions

- [DriftCompensationWindow](#) ()
Initializes a new instance of the [DriftCompensationWindow](#) class.

8.29.1 Detailed Description

Interaction logic for DriftCompensation.xaml

8.29.2 Constructor & Destructor Documentation

8.29.2.1 DriftCompensationWindow()

```
CustomCalibrationLibrary.Views.DriftCompensationWindow.DriftCompensationWindow ( ) [inline]
```

Initializes a new instance of the [DriftCompensationWindow](#) class.

The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Views/DriftCompensationWindow.xaml.cs

8.30 GazeUtilityLibrary.DataStructs.EyeData Class Reference

The eye data set, including pupil information.

Public Member Functions

- [EyeData](#) (float pupilDiameter, bool isPupilDiameterValid)
Initializes a new instance of the [EyeData](#) class.

Properties

- float [PupilDiameter](#) [get]
The diameter of the pupil
- bool [IsPupilDiameterValid](#) [get]
The validity flag of the pupil diameter

8.30.1 Detailed Description

The eye data set, including pupil information.

8.30.2 Constructor & Destructor Documentation

8.30.2.1 EyeData()

```
GazeUtilityLibrary.DataStructs.EyeData.EyeData (
    float pupilDiameter,
    bool isPupilDiameterValid ) [inline]
```

Initializes a new instance of the [EyeData](#) class.

Parameters

<i>pupilDiameter</i>	The pupil diameter.
<i>isPupilDiameterValid</i>	The validity of the pupil diameter.

8.30.3 Property Documentation

8.30.3.1 IsPupilDiameterValid

```
bool GazeUtilityLibrary.DataStructs.EyeData.IsPupilDiameterValid [get]
```

The validity flag of the pupil diameter

8.30.3.2 PupilDiameter

```
float GazeUtilityLibrary.DataStructs.EyeData.PupilDiameter [get]
```

The diameter of the pupil


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/DataStructs/EyeData.cs


8.31 GazeUtilityLibrary.Tracker.EyeTrackerPro Class Reference

Interface to the [Tobii](#) SDK Pro engine

Inheritance diagram for GazeUtilityLibrary.Tracker.EyeTrackerPro:

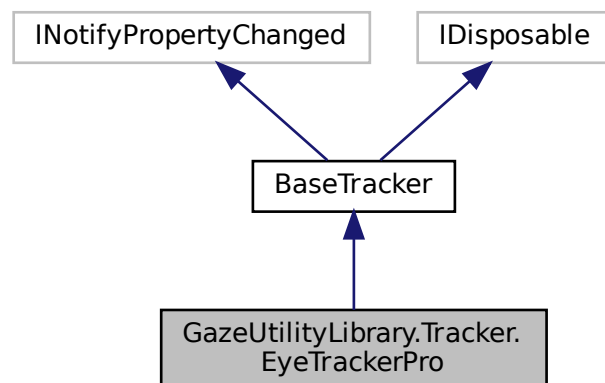

Collaboration diagram for GazeUtilityLibrary.Tracker.EyeTrackerPro:

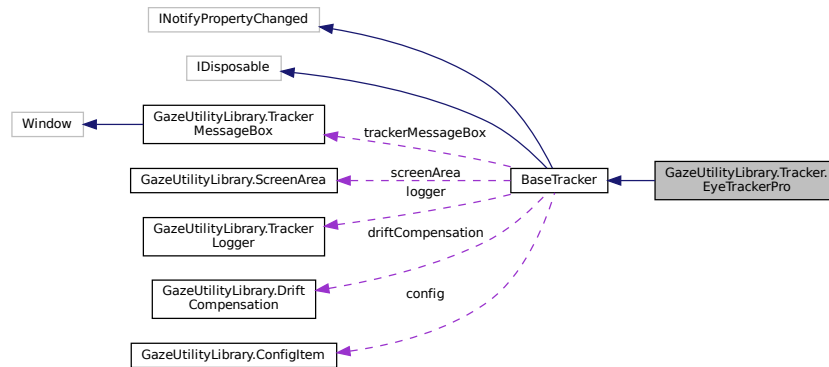

Public Member Functions

- [EyeTrackerPro](#) ([TrackerLogger](#) logger, [ConfigItem](#) config)
Initializes a new instance of the [EyeTrackerPro](#) class.
- override async Task [InitCalibrationAsync](#) ()
Initialise the screen based calibration.
- override void [InitCalibration](#) ()
Initialise the screen based calibration.
- override void [InitValidation](#) ()
Initialise the screen based calibration.
- override async Task< bool > [CollectCalibrationDataAsync](#) (Point point)
Collects gaze data of a calibration point.
- override async Task< bool > [CollectValidationDataAsync](#) (Point point)
Collects gaze data of a validation point.
- override async Task [FinishCalibrationAsync](#) ()
Finish the screen based async calibration process.
- override void [FinishCalibration](#) ()
Finish the screen based calibration process.
- override void [FinishValidation](#) ()
Finish the screen based validation process.
- override async Task< List< [GazeCalibrationData](#) > > [ApplyCalibration](#) ()
Compute and apply the calibration data. Transform the Tobii calibration result into the [GazeCalibrationData](#) structure.
- override? [GazeValidationData](#) [ComputeValidation](#) ()
Compute the validation data.
- bool [IsLicenseOk](#) ()
Determines whether the license is applied to the eyetracker device
- override bool [IsInitialised](#) ()
Checks if the tracker device exists.
- override string [PatternReplace](#) (string pattern)
Replaces a patten string with information from the eye tracker. Supported patterns are S for the serial number and A for the address.

Protected Member Functions

- override void [InitDriftCompensation](#) ()
Initialise the drift compensation.
- override int [GetFixationFrameCount](#) (int durationThreshold)
Get the number of required gaze samples to compute a fixation. This is based on the duration threshold and the sample rate of the device.
- override Vector3 [GetUnitDirection](#) ()
Get the unit vector pointing in the direction of the gaze vector.

Additional Inherited Members

8.31.1 Detailed Description

Interface to the [Tobii](#) SDK Pro engine

See also

[GazeHelper.TrackerHandler](#)

8.31.2 Constructor & Destructor Documentation

8.31.2.1 EyeTrackerPro()

```
GazeUtilityLibrary.Tracker.EyeTrackerPro.EyeTrackerPro (
    TrackerLogger logger,
    ConfigItem config ) [inline]
```

Initializes a new instance of the [EyeTrackerPro](#) class.

Parameters

<i>logger</i>	The logger.
<i>config</i>	The config item.

8.31.3 Member Function Documentation

8.31.3.1 ApplyCalibration()

```
override async Task<List<GazeCalibrationData> > GazeUtilityLibrary.Tracker.EyeTrackerPro.↔
ApplyCalibration ( ) [inline], [virtual]
```

Compute and apply the calibration data. Transform the Tobii calibration result into the [GazeCalibrationData](#) structure.

Returns

The calibration data result wrapped by an async handler.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.2 CollectCalibrationDataAsync()

```
override async Task<bool> GazeUtilityLibrary.Tracker.EyeTrackerPro.CollectCalibrationDataAsync  
(  
    Point point ) [inline], [virtual]
```

Collects gaze data of a calibration point.

Parameters

<i>point</i>	
--------------	--

Returns

True on success, false on failure, wrapped by an async handler.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.3 CollectValidationDataAsync()

```
override async Task<bool> GazeUtilityLibrary.Tracker.EyeTrackerPro.CollectValidationDataAsync  
(  
    Point point ) [inline], [virtual]
```

Collects gaze data of a validation point.

Parameters

<i>point</i>	
--------------	--

Returns

True on success, false on failure, wrapped by an async handler.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.4 ComputeValidation()

```
override? GazeValidationData GazeUtilityLibrary.Tracker.EyeTrackerPro.ComputeValidation ( )  
[inline], [virtual]
```

Compute the validation data.

Returns

The validation data result.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.5 FinishCalibration()

```
override void GazeUtilityLibrary.Tracker.EyeTrackerPro.FinishCalibration ( ) [inline], [virtual]
```

Finish the screen based calibration process.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.6 FinishCalibrationAsync()

```
override async Task GazeUtilityLibrary.Tracker.EyeTrackerPro.FinishCalibrationAsync ( ) [inline],  
[virtual]
```

Finish the screen based async calibration process.

Returns

An async handler

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.7 FinishValidation()

```
override void GazeUtilityLibrary.Tracker.EyeTrackerPro.FinishValidation ( ) [inline], [virtual]
```

Finish the screen based validation process.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.8 GetFixationFrameCount()

```
override int GazeUtilityLibrary.Tracker.EyeTrackerPro.GetFixationFrameCount (   
    int durationThreshold ) [inline], [protected], [virtual]
```

Get the number of required gaze samples to compute a fixation. This is based on the duration threshold and the sample rate of the device.

Parameters

<i>durationThreshold</i>	The required fixation duration in milliseconds.
--------------------------	---

Returns

The number of required samples.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.9 GetUnitDirection()

```
override Vector3 GazeUtilityLibrary.Tracker.EyeTrackerPro.GetUnitDirection ( ) [inline],  
[protected], [virtual]
```

Get the unit vector pointing in the direction of the gaze vector.

Returns

The unit vector pointing in the negative z direction.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.10 InitCalibration()

```
override void GazeUtilityLibrary.Tracker.EyeTrackerPro.InitCalibration ( ) [inline], [virtual]
```

Initialise the screen based calibration.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.11 InitCalibrationAsync()

```
override async Task GazeUtilityLibrary.Tracker.EyeTrackerPro.InitCalibrationAsync ( ) [inline],  
[virtual]
```

Initialise the screen based calibration.

Returns

An async handler

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.12 InitDriftCompensation()

```
override void GazeUtilityLibrary.Tracker.EyeTrackerPro.InitDriftCompensation ( ) [inline],  
[protected], [virtual]
```

Initialise the drift compensation.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.13 InitValidation()

```
override void GazeUtilityLibrary.Tracker.EyeTrackerPro.InitValidation ( ) [inline], [virtual]
```

Initialise the screen based calibration.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.14 IsInitialised()

```
override bool GazeUtilityLibrary.Tracker.EyeTrackerPro.IsInitialised ( ) [inline], [virtual]
```

Checks if the tracker device exists.

Returns

True if the tracker device exists, false otherwise.

Reimplemented from [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.31.3.15 IsLicenseOk()

```
bool GazeUtilityLibrary.Tracker.EyeTrackerPro.IsLicenseOk ( ) [inline]
```

Determines whether the license is applied to the eyetracker device

Returns

true if [is license ok]; otherwise, false.

8.31.3.16 PatternReplace()

```
override string GazeUtilityLibrary.Tracker.EyeTrackerPro.PatternReplace (
    string pattern ) [inline], [virtual]
```

Replaces a patten string with information from the eye tracker. Supported patterns are S for the serial number and A for the address.

Returns

The string where patterns were replaced.

Reimplemented from [GazeUtilityLibrary.Tracker.BaseTracker](#).


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/Tracker/EyeTrackerPro.cs


8.32 CustomCalibrationLibrary.Views.FixationPoint Class Reference

Interaction logic for FixationPoint.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.FixationPoint:

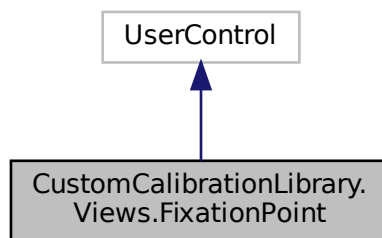

Collaboration diagram for CustomCalibrationLibrary.Views.FixationPoint:

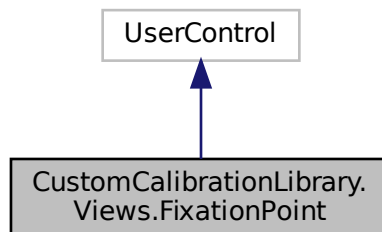

Public Member Functions

- [FixationPoint](#) ()

Initializes a new instance of the [FixationPoint](#) class.

8.32.1 Detailed Description

Interaction logic for FixationPoint.xaml

8.32.2 Constructor & Destructor Documentation

8.32.2.1 FixationPoint()

```
CustomCalibrationLibrary.Views.FixationPoint.FixationPoint ( ) [inline]
```

Initializes a new instance of the [FixationPoint](#) class.

The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Views/FixationPoint.xaml.cs

8.33 GazeUtilityLibrary.DataStructs.GazeCalibrationData Class Reference

The gaze calibration data structure

Public Member Functions

- [GazeCalibrationData](#) (double xCoord, double yCoord, double xCoordLeft, double yCoordLeft, bool validityLeft, double xCoordRight, double yCoordRight, bool validityRight)

Initializes a new instance of the [GazeDataArgs](#) class.

- string[] [Prepare](#) ([ConfigItem](#) config)

Prepare a list of formatted calibration data values

Properties

- double **XCoord** [get]
The x coordinate of the calibration point.
- double **YCoord** [get]
The y coordinate of the calibration point.
- double **XCoordLeft** [get]
The x coord of the gaze point of the left eye.
- double **YCoordLeft** [get]
The y coord of the gaze point of the left eye.
- bool **ValidityLeft** [get]
The validity of gaze point coordinate of the left eye.
- double **XCoordRight** [get]
The x coord of the gaze point of the right eye.
- double **YCoordRight** [get]
The y coord of the gaze point of the right eye.
- bool **ValidityRight** [get]
The validity of gaze point coordinate of the right eye.

8.33.1 Detailed Description

The gaze calibration data structure

8.33.2 Constructor & Destructor Documentation

8.33.2.1 GazeCalibrationData()

```
GazeUtilityLibrary.DataStructs.GazeCalibrationData.GazeCalibrationData (
    double xCoord,
    double yCoord,
    double xCoordLeft,
    double yCoordLeft,
    bool validityLeft,
    double xCoordRight,
    double yCoordRight,
    bool validityRight ) [inline]
```

Initializes a new instance of the GazeDataArgs class.

Parameters

<i>xCoord</i>	The x coord of the calibration point.
<i>yCoord</i>	The y coord of the calibration point.
<i>xCoordLeft</i>	The x coord of the gaze point of the left eye.
<i>yCoordLeft</i>	The y coord of the gaze point of the left eye.
<i>validityLeft</i>	The validity of gaze point coordinate of the left eye.
<i>xCoordRight</i>	The x coord of the gaze point of the right eye.
<i>yCoordRight</i>	The y coord of the gaze point of the right eye.
<i>validityRight</i>	The validity of gaze point coordinate of the right eye.

8.33.3 Member Function Documentation

8.33.3.1 Prepare()

```
string [ ] GazeUtilityLibrary.DataStructs.GazeCalibrationData.Prepare (
    ConfigItem config ) [inline]
```

Prepare a list of formatted calibration data values

Parameters

<i>config</i>	The gaze configuration structure
---------------	----------------------------------

Returns

A list of formatted values. Each index corresponds to a specific value. This allows to reorder the list according to a format string.

8.33.4 Property Documentation

8.33.4.1 ValidityLeft

```
bool GazeUtilityLibrary.DataStructs.GazeCalibrationData.ValidityLeft [get]
```

The validity of gaze point coordinate of the left eye.

8.33.4.2 ValidityRight

```
bool GazeUtilityLibrary.DataStructs.GazeCalibrationData.ValidityRight [get]
```

The validity of gaze point coordinate of the right eye.

8.33.4.3 XCoord

```
double GazeUtilityLibrary.DataStructs.GazeCalibrationData.XCoord [get]
```

The x coordinate of the calibration point.

8.33.4.4 XCoordLeft

```
double GazeUtilityLibrary.DataStructs.GazeCalibrationData.XCoordLeft [get]
```

The x coord of the gaze point of the left eye.

8.33.4.5 XCoordRight

```
double GazeUtilityLibrary.DataStructs.GazeCalibrationData.XCoordRight [get]
```

The x coord of the gaze point of the right eye.

8.33.4.6 YCoord

```
double GazeUtilityLibrary.DataStructs.GazeCalibrationData.YCoord [get]
```

The y coordinate of the calibration point.

8.33.4.7 YCoordLeft

```
double GazeUtilityLibrary.DataStructs.GazeCalibrationData.YCoordLeft [get]
```

The y coord of the gaze point of the left eye.

8.33.4.8 YCoordRight

```
double GazeUtilityLibrary.DataStructs.GazeCalibrationData.YCoordRight [get]
```

The y coord of the gaze point of the right eye.


The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/DataStructs/GazeCalibrationData.cs`


8.34 GazeUtilityLibrary.GazeConfigError Class Reference

The gaze config error class to convert error flags to binary strings.

Inheritance diagram for GazeUtilityLibrary.GazeConfigError:

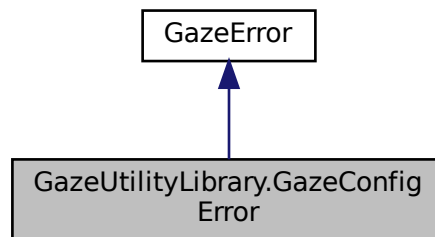

Collaboration diagram for GazeUtilityLibrary.GazeConfigError:

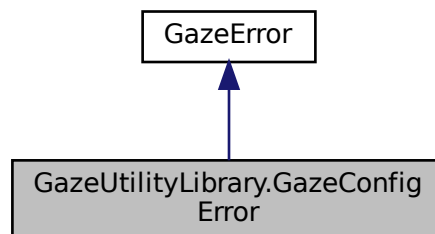

Public Member Functions

- `string GetGazeConfigErrorString ()`
Gets the gaze error string.

Properties

- `EGazeConfigError Error` [set]
The error flags.

Additional Inherited Members

8.34.1 Detailed Description

The gaze config error class to convert error flags to binary strings.

8.34.2 Member Function Documentation

8.34.2.1 GetGazeConfigErrorString()

```
string GazeUtilityLibrary.GazeConfigError.GetGazeConfigErrorString ( ) [inline]
```

Gets the gaze error string.

Returns

the error string with binary error values if errors occurred, the empty string otherwise

8.34.3 Property Documentation

8.34.3.1 Error

```
EGazeConfigError GazeUtilityLibrary.GazeConfigError.Error [set]
```

The error flags.

The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/GazeError.cs

8.35 GazeUtilityLibrary.GazeConfiguration Class Reference

The gaze configuration handler.

Public Member Functions

- [GazeConfiguration](#) ([TrackerLogger](#) logger)
Initializes a new instance of the [GazeConfiguration](#) class.
- bool [InitConfig](#) ()
Initialise the gaze configuration by parsing and checking the configuration file.
- bool [CleanupGazeOutputFile](#) (string error)
Close the gaze outputfile and rename it by appending error codes.
- bool [CleanupCalibrationOutputFile](#) (string error)
Close the calibration outputfile and rename it by appending error codes.
- bool [CleanupValidationOutputFile](#) (string error)
Close the validation outputfile and rename it by appending error codes.
- bool [DumpCurrentConfigurationFile](#) ()
Dump current configuration to the disk.
- bool [PrepareGazeOutputFile](#) (string? subjectCode, string? outputPath)
Prepare the gaze output file based on the configuration.
- bool [PrepareCalibrationOutputFile](#) (string? subjectCode)
Prepare the calibration output file based on the configuration.
- bool [PrepareValidationOutputFile](#) (string? subjectCode)
Prepare the validation output file based on the configuration.
- void [WriteToGazeOutput](#) (string[] formatted_values)
Write to the gaze output file
- void [WriteToCalibrationOutput](#) (string[] formatted_values)
Write to the calibration output file
- void [WriteToValidationOutput](#) (string[] formatted_values)
Write to the calibration output file

Properties

- [ConfigItem?? Config](#) [get]
The JSON structure holding the configuratyon options.

8.35.1 Detailed Description

The gaze configuration handler.

8.35.2 Constructor & Destructor Documentation

8.35.2.1 GazeConfiguration()

```
GazeUtilityLibrary.GazeConfiguration.GazeConfiguration (
    TrackerLogger logger ) [inline]
```

Initializes a new instance of the [GazeConfiguration](#) class.

Parameters

<i>logger</i>	The log handler.
---------------	------------------

8.35.3 Member Function Documentation

8.35.3.1 CleanupCalibrationOutputFile()

```
bool GazeUtilityLibrary.GazeConfiguration.CleanupCalibrationOutputFile (  
    string error ) [inline]
```

Close the calibration outputfile and rename it by appending error codes.

Parameters

<i>error</i>	
--------------	--

Returns

True on success, False on failure.

8.35.3.2 CleanupGazeOutputFile()

```
bool GazeUtilityLibrary.GazeConfiguration.CleanupGazeOutputFile (  
    string error ) [inline]
```

Close the gaze outputfile and rename it by appending error codes.

Parameters

<i>error</i>	
--------------	--

Returns

True on success, False on failure.

8.35.3.3 CleanupValidationOutputFile()

```
bool GazeUtilityLibrary.GazeConfiguration.CleanupValidationOutputFile (  
    string error ) [inline]
```

Close the validation outputfile and rename it by appending error codes.

Parameters

<i>error</i>	
--------------	--

Returns

True on success, False on failure.

8.35.3.4 DumpCurrentConfigurationFile()

```
bool GazeUtilityLibrary.GazeConfiguration.DumpCurrentConfigurationFile ( ) [inline]
```

Dump current configuration to the disk.

Returns

True on success, False on failure.

8.35.3.5 InitConfig()

```
bool GazeUtilityLibrary.GazeConfiguration.InitConfig ( ) [inline]
```

Initialise the gaze configuration by parsing and checking the configuration file.

Returns

True on success, False on failure.

8.35.3.6 PrepareCalibrationOutputFile()

```
bool GazeUtilityLibrary.GazeConfiguration.PrepareCalibrationOutputFile (
    string? subjectCode ) [inline]
```

Prepare the calibration output file based on the configuration.

Parameters

<i>subjectCode</i>	An optional subject code to be appended to the file name if set.
--------------------	--

Returns

True on success, False on failure.

8.35.3.7 PrepareGazeOutputFile()

```
bool GazeUtilityLibrary.GazeConfiguration.PrepareGazeOutputFile (
    string?  subjectCode,
    string?  outputPath ) [inline]
```

Prepare the gaze output file based on the configuration.

Parameters

<i>subjectCode</i>	An optional subject code to be appended to the file name if set.
<i>outputPath</i>	An optional output path where the file will be stored.

Returns

True on success, False on failure.

8.35.3.8 PrepareValidationOutputFile()

```
bool GazeUtilityLibrary.GazeConfiguration.PrepareValidationOutputFile (
    string?  subjectCode ) [inline]
```

Prepare the validation output file based on the configuration.

Parameters

<i>subjectCode</i>	An optional subject code to be appended to the file name if set.
--------------------	--

Returns

True on success, False on failure.

8.35.3.9 WriteToCalibrationOutput()

```
void GazeUtilityLibrary.GazeConfiguration.WriteToCalibrationOutput (
    string[] formatted_values ) [inline]
```

Write to the calibration output file

Parameters

<i>formatted_values</i>	The list of formatted values to be written to the file.
-------------------------	---

8.35.3.10 WriteToGazeOutput()

```
void GazeUtilityLibrary.GazeConfiguration.WriteToGazeOutput (
    string[] formatted_values ) [inline]
```

Write to the gaze output file

Parameters

<i>formatted_values</i>	The list of formatted values to be written to the file.
-------------------------	---

8.35.3.11 WriteToValidationOutput()

```
void GazeUtilityLibrary.GazeConfiguration.WriteToValidationOutput (
    string[] formatted_values ) [inline]
```

Write to the calibration output file

Parameters

<i>formatted_values</i>	The list of formatted values to be written to the file.
-------------------------	---

8.35.4 Property Documentation**8.35.4.1 Config**

```
ConfigItem?? GazeUtilityLibrary.GazeConfiguration.Config [get]
```

The JSON structure holding the configuratyion options.

The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/GazeConfiguration.cs

8.36 GazeUtilityLibrary.DataStructs.GazeData Class Reference

The class definition of a gaze data set

Public Member Functions

- [GazeData](#) (TimeSpan timestamp, TimeSpan timestampReceived, Vector2 gazePoint2d, bool isGazePoint2dValid)
Initializes a new instance of the GazeDataArgs class.
- [GazeData](#) (TimeSpan timestamp, TimeSpan timestampReceived, Vector2 gazePoint2dLeft, bool isGazePoint2dValidLeft, Vector2 gazePoint2dRight, bool isGazePoint2dValidRight)
Initializes a new instance of the GazeDataArgs class.
- [GazeData](#) (TimeSpan timestamp, TimeSpan timestampReceived, Vector2 gazePoint2dLeft, bool isGazePoint2dValidLeft, Vector2 gazePoint2dRight, bool isGazePoint2dValidRight, Vector3 gazePoint3dLeft, bool isGazePoint3dValidLeft, Vector3 gazePoint3dRight, bool isGazePoint3dValidRight, Vector3 gazeOrigin3dLeft, bool isGazeOrigin3dValidLeft, Vector3 gazeOrigin3dRight, bool isGazeOrigin3dValidRight, float pupilDiameterLeft, bool isPupilDiameterValidLeft, float pupilDiameterRight, bool isPupilDiameterValidRight)
Initializes a new instance of the GazeDataArgs class.
- string[] [Prepare](#) ([ConfigItem](#) config, int trialId, string tag, TimeSpan startTime)
Prepare a list of formatted gaze data values

Properties

- TimeSpan [Timestamp](#) [get]
The timestamp of the data sample.
- TimeSpan [TimestampReceived](#) [get]
The device timestamp of the data sample.
- [GazeDataCollection? Left](#) [get]
The gaze data set, including 2d and (optionally) 3d gaze data as well as optional eye data of the left eye.
- [GazeDataCollection? Right](#) [get]
The gaze data set, including 2d and (optionally) 3d gaze data as well as optional eye data of the right eye.
- [GazeDataCollection Combined](#) [get]
The gaze data set, including 2d and (optionally) 3d gaze data as well as optional eye data of the combined eyes.
- [DriftCompensationData? DriftCompensation](#) [get, set]
The drift compensation information.

8.36.1 Detailed Description

The class definition of a gaze data set

8.36.2 Constructor & Destructor Documentation

8.36.2.1 GazeData() [1/3]

```
GazeUtilityLibrary.DataStructs.GazeData.GazeData (
    TimeSpan timestamp,
    TimeSpan timestampReceived,
    Vector2 gazePoint2d,
    bool isGazePoint2dValid ) [inline]
```

Initializes a new instance of the GazeDataArgs class.

Parameters

<i>timestamp</i>	The timestamp when the data was captured by the device.
<i>timestampReceived</i>	The timestamp when the data was received by the system.
<i>gazePoint2d</i>	The 2d coordinates of the combined gaze point.
<i>isGazePoint2dValid</i>	The validity of the combined 2d gaze point.

8.36.2.2 GazeData() [2/3]

```
GazeUtilityLibrary.DataStructs.GazeData.GazeData (
    TimeSpan timestamp,
    TimeSpan timestampReceived,
    Vector2 gazePoint2dLeft,
    bool isGazePoint2dValidLeft,
    Vector2 gazePoint2dRight,
    bool isGazePoint2dValidRight ) [inline]
```

Initializes a new instance of the GazeDataArgs class.

Parameters

<i>timestamp</i>	The timestamp when the data was captured by the device.
<i>timestampReceived</i>	The timestamp when the data was received by the system.
<i>gazePoint2dLeft</i>	The 2d coordinates of the left gaze point.
<i>isGazePoint2dValidLeft</i>	The validity of the left 2d gaze point.
<i>gazePoint2dRight</i>	The 2d coordinates of the right gaze point.
<i>isGazePoint2dValidRight</i>	The validity of the right 2d gaze point.

8.36.2.3 GazeData() [3/3]

```
GazeUtilityLibrary.DataStructs.GazeData.GazeData (
    TimeSpan timestamp,
    TimeSpan timestampReceived,
    Vector2 gazePoint2dLeft,
    bool isGazePoint2dValidLeft,
    Vector2 gazePoint2dRight,
    bool isGazePoint2dValidRight,
    Vector3 gazePoint3dLeft,
    bool isGazePoint3dValidLeft,
    Vector3 gazePoint3dRight,
    bool isGazePoint3dValidRight,
    Vector3 gazeOrigin3dLeft,
    bool isGazeOrigin3dValidLeft,
    Vector3 gazeOrigin3dRight,
    bool isGazeOrigin3dValidRight,
    float pupilDiameterLeft,
```

```

bool isPupilDiameterValidLeft,
float pupilDiameterRight,
bool isPupilDiameterValidRight ) [inline]

```

Initializes a new instance of the GazeDataArgs class.

Parameters

<i>timestamp</i>	The timestamp when the data was captured by the device.
<i>timestampReceived</i>	The timestamp when the data was received by the system.
<i>gazePoint2dLeft</i>	The 2d coordinates of the left gaze point.
<i>isGazePoint2dValidLeft</i>	The validity of the left 2d gaze point.
<i>gazePoint2dRight</i>	The 2d coordinates of the right gaze point.
<i>isGazePoint2dValidRight</i>	The validity of the right 2d gaze point.
<i>gazePoint3dLeft</i>	The 3d coordinates of the left gaze point.
<i>isGazePoint3dValidLeft</i>	The validity of the left 3d gaze point.
<i>gazePoint3dRight</i>	The 3d coordinates of the right gaze point.
<i>isGazePoint3dValidRight</i>	The validity of the right 3d gaze point.
<i>gazeOrigin3dLeft</i>	The 3d coordinates of the left gaze origin.
<i>isGazeOrigin3dValidLeft</i>	The validity of the left 3d gaze origin.
<i>gazeOrigin3dRight</i>	The 3d coordinates of the right gaze origin.
<i>isGazeOrigin3dValidRight</i>	The validity of the right 3d gaze origin.
<i>pupilDiameterLeft</i>	The pupil diameter the left eye.
<i>isPupilDiameterValidLeft</i>	The validity of the left pupil diameter.
<i>pupilDiameterRight</i>	The pupil diameter the left eye.
<i>isPupilDiameterValidRight</i>	The validity of the left pupil diameter.

8.36.3 Member Function Documentation

8.36.3.1 Prepare()

```

string [] GazeUtilityLibrary.DataStructs.GazeData.Prepare (
    ConfigItem config,
    int trialId,
    string tag,
    TimeSpan startTime ) [inline]

```

Prepare a list of formatted gaze data values

Parameters

<i>config</i>	The gaze configuration structure
<i>trialId</i>	The ID of the current trial.
<i>tag</i>	An arbitrary tag to associate with the data sample.
<i>startTime</i>	The system time to use to compute the relative timestamp

Returns

A list of formatted values. Each index corresponds to a specific value. This allows to reorder the list according to a format string.

8.36.4 Property Documentation

8.36.4.1 Combined

`GazeDataCollection` `GazeUtilityLibrary.DataStructs.GazeData.Combined` [get]

The gaze data set, including 2d and (optionally) 3d gaze data as well as optional eye data of the combined eyes.

8.36.4.2 DriftCompensation

`DriftCompensationData?` `GazeUtilityLibrary.DataStructs.GazeData.DriftCompensation` [get], [set]

The drift compensation information.

8.36.4.3 Left

`GazeDataCollection?` `GazeUtilityLibrary.DataStructs.GazeData.Left` [get]

The gaze data set, including 2d and (optionally) 3d gaze data as well as optional eye data of the left eye.

8.36.4.4 Right

`GazeDataCollection?` `GazeUtilityLibrary.DataStructs.GazeData.Right` [get]

The gaze data set, including 2d and (optionally) 3d gaze data as well as optional eye data of the right eye.

8.36.4.5 Timestamp

`TimeSpan` `GazeUtilityLibrary.DataStructs.GazeData.Timestamp` [get]

The timestamp of the data sample.

8.36.4.6 TimestampReceived

TimeSpan GazeUtilityLibrary.DataStructs.GazeData.TimestampReceived [get]

The device timestamp of the data sample.

The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/DataStructs/GazeData.cs

8.37 GazeUtilityLibrary.DataStructs.GazeData2d Class Reference

The 2d gaze data set.

Public Member Functions

- [GazeData2d](#) (Vector2 gazePoint, bool isGazePointValid)

Initializes a new instance of the [GazeData2d](#) class.

Properties

- Vector2 [GazePoint](#) [get]

The 2d gaze point.

- bool [IsGazePointValid](#) [get]

The validity flag of the 2d gaze point.

8.37.1 Detailed Description

The 2d gaze data set.

8.37.2 Constructor & Destructor Documentation

8.37.2.1 GazeData2d()

```
GazeUtilityLibrary.DataStructs.GazeData2d.GazeData2d (  
    Vector2 gazePoint,  
    bool isGazePointValid ) [inline]
```

Initializes a new instance of the [GazeData2d](#) class.

Parameters

<i>gazePoint</i>	The 2d coordinates of the gaze point.
<i>isGazePointValid</i>	The validity of the 2d gaze point.

8.37.3 Property Documentation

8.37.3.1 GazePoint

`Vector2 GazeUtilityLibrary.DataStructs.GazeData2d.GazePoint [get]`

The 2d gaze point.

8.37.3.2 IsGazePointValid

`bool GazeUtilityLibrary.DataStructs.GazeData2d.IsGazePointValid [get]`

The validity flag of the 2d gaze point.

The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/DataStructs/GazeData2d.cs`

8.38 GazeUtilityLibrary.DataStructs.GazeData3d Class Reference

The 3d gaze data set.

Public Member Functions

- [GazeData3d](#) (Vector3 gazePoint, bool isGazePointValid, Vector3 gazeOrigin, bool isGazeOriginValid)
Initializes a new instance of the [GazeData3d](#) class.

Properties

- Vector3 [GazePoint](#) [get]
The 3d gaze point.
- bool [IsGazePointValid](#) [get]
The validity of the 3d gaze point.
- Vector3 [GazeOrigin](#) [get]
The 3d origin of the gaze.
- Vector3 [GazeDirection](#) [get]
The 3d gaze direction vector.
- float [GazeDistance](#) [get]
The gaze distance from the origin to the gaze point.
- bool [IsGazeOriginValid](#) [get]
The validity of the 3d origin.

8.38.1 Detailed Description

The 3d gaze data set.

8.38.2 Constructor & Destructor Documentation

8.38.2.1 GazeData3d()

```
GazeUtilityLibrary.DataStructs.GazeData3d.GazeData3d (  
    Vector3 gazePoint,  
    bool isGazePointValid,  
    Vector3 gazeOrigin,  
    bool isGazeOriginValid ) [inline]
```

Initializes a new instance of the [GazeData3d](#) class.

Parameters

<i>gazePoint</i>	The 3d coordinates of the gaze point.
<i>isGazePointValid</i>	The validity of the 3d gaze point.
<i>gazeOrigin</i>	The 3d coordinates of the gaze origin.
<i>isGazeOriginValid</i>	The validity of the 3d gaze origin.

8.38.3 Property Documentation

8.38.3.1 GazeDirection

```
Vector3 GazeUtilityLibrary.DataStructs.GazeData3d.GazeDirection [get]
```

The 3d gaze direction vector.

8.38.3.2 GazeDistance

```
float GazeUtilityLibrary.DataStructs.GazeData3d.GazeDistance [get]
```

The gaze distance from the origin to the gaze point.

8.38.3.3 GazeOrigin

```
Vector3 GazeUtilityLibrary.DataStructs.GazeData3d.GazeOrigin [get]
```

The 3d origin of the gaze.

8.38.3.4 GazePoint

```
Vector3 GazeUtilityLibrary.DataStructs.GazeData3d.GazePoint [get]
```

The 3d gaze point.

8.38.3.5 IsGazeOriginValid

```
bool GazeUtilityLibrary.DataStructs.GazeData3d.IsGazeOriginValid [get]
```

The validity of the 3d origin.

8.38.3.6 IsGazePointValid

```
bool GazeUtilityLibrary.DataStructs.GazeData3d.IsGazePointValid [get]
```

The validity of the 3d gaze point.

The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/DataStructs/GazeData3d.cs`

8.39 GazeUtilityLibrary.DataStructs.GazeDataCollection Class Reference

The gaze data set, including 2d and (optionally) 3d gaze data as well as optional eye data.

Public Member Functions

- [GazeDataCollection](#) (Vector2 gazePoint2d, bool isGazePoint2dValid)
Initializes a new instance of the GazeDataItem class.
- [GazeDataCollection](#) (Vector2 gazePoint2d, bool isGazePoint2dValid, Vector3 gazePoint3d, bool isGazePoint3dValid, Vector3 gazeOrigin3d, bool isGazeOrigin3dValid, float pupilDiameter, bool isPupilDiameterValid)
Initializes a new instance of the GazeDataItem class.

Properties

- [GazeData2d](#) [GazeData2d](#) [get]
The 2d gaze data.
- [GazeData3d?](#) [GazeData3d](#) [get]
The 3d gaze data.
- [EyeData?](#) [EyeData](#) [get]
Pupil data of the eye.

8.39.1 Detailed Description

The gaze data set, including 2d and (optionally) 3d gaze data as well as optional eye data.

8.39.2 Constructor & Destructor Documentation

8.39.2.1 GazeDataCollection() [1/2]

```
GazeUtilityLibrary.DataStructs.GazeDataCollection.GazeDataCollection (
    Vector2 gazePoint2d,
    bool isGazePoint2dValid ) [inline]
```

Initializes a new instance of the GazeDataItem class.

Parameters

<i>gazePoint2d</i>	The 2d coordinates of the gaze point.
<i>isGazePoint2dValid</i>	The validity of the 2d gaze point.

8.39.2.2 GazeDataCollection() [2/2]

```
GazeUtilityLibrary.DataStructs.GazeDataCollection.GazeDataCollection (
    Vector2 gazePoint2d,
    bool isGazePoint2dValid,
    Vector3 gazePoint3d,
    bool isGazePoint3dValid,
    Vector3 gazeOrigin3d,
    bool isGazeOrigin3dValid,
    float pupilDiameter,
    bool isPupilDiameterValid ) [inline]
```

Initializes a new instance of the GazeDataItem class.

Parameters

<i>gazePoint2d</i>	The 2d coordinates of the gaze point.
<i>isGazePoint2dValid</i>	The validity of the 2d gaze point.
<i>gazePoint3d</i>	The 3d coordinates of the gaze point.
<i>isGazePoint3dValid</i>	The validity of the 3d gaze point.
<i>gazeOrigin3d</i>	The 3d coordinates of the gaze origin.
<i>isGazeOrigin3dValid</i>	The validity of the 3d gaze origin.
<i>pupilDiameter</i>	The pupil diameter.
<i>isPupilDiameterValid</i>	The validity of the pupil diameter.

8.39.3 Property Documentation

8.39.3.1 EyeData

`EyeData?` `GazeUtilityLibrary.DataStructs.GazeDataCollection.EyeData` [get]

Pupil data of the eye.

8.39.3.2 GazeData2d

`GazeData2d` `GazeUtilityLibrary.DataStructs.GazeDataCollection.GazeData2d` [get]

The 2d gaze data.

8.39.3.3 GazeData3d

`GazeData3d?` `GazeUtilityLibrary.DataStructs.GazeDataCollection.GazeData3d` [get]

The 3d gaze data.


The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/DataStructs/GazeDataCollection.cs`


8.40 GazeUtilityLibrary.GazeDataError Class Reference

The gaze data error class to convert error flags to binary strings.

Inheritance diagram for GazeUtilityLibrary.GazeDataError:

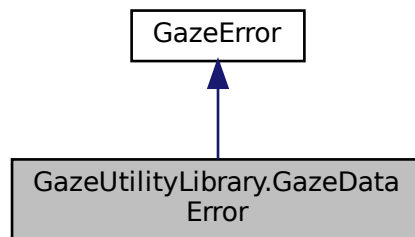

Collaboration diagram for GazeUtilityLibrary.GazeDataError:

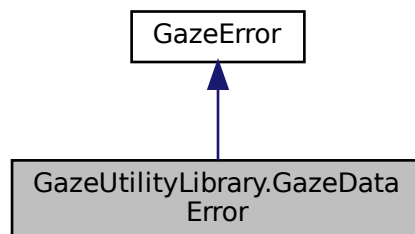

Public Member Functions

- string [GetGazeDataErrorString](#) ()
Gets the gaze error string.

Properties

- [EGazeDataError Error](#) [set]
The error flags.

Additional Inherited Members

8.40.1 Detailed Description

The gaze data error class to convert error flags to binary strings.

8.40.2 Member Function Documentation

8.40.2.1 GetGazeDataErrorString()

```
string GazeUtilityLibrary.GazeDataError.GetGazeDataErrorString ( ) [inline]
```

Gets the gaze error string.

Returns

the error string with binary error values if errors occurred, the empty string otherwise

8.40.3 Property Documentation

8.40.3.1 Error

```
EGazeDataError GazeUtilityLibrary.GazeDataError.Error [set]
```

The error flags.


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/GazeError.cs

8.41 GazeUtilityLibrary.GazeError Class Reference

The base error class to convert error flags to binary strings.

Inheritance diagram for GazeUtilityLibrary.GazeError:

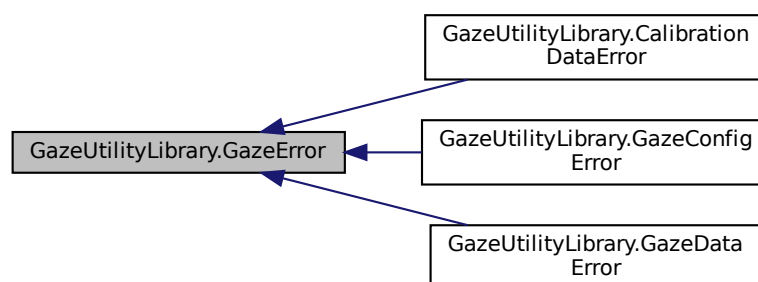

Protected Member Functions

- string [ConvertToBinString](#) (int val, int len)
Converts a integer value to a binary string.

8.41.1 Detailed Description

The base error class to convert error flags to binary strings.

8.41.2 Member Function Documentation

8.41.2.1 ConvertToBinString()

```
string GazeUtilityLibrary.GazeError.ConvertToBinString (
    int val,
    int len ) [inline], [protected]
```

Converts a integer value to a binary string.

Parameters

<i>val</i>	The value.
<i>len</i>	The length of the binary string.

Returns

a binary string of specified length, left-padded with '0'

The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/GazeError.cs

8.42 GazeUtilityLibrary.DataStructs.GazeValidationData Class Reference

The gaze validation data structure

Public Member Functions

- [GazeValidationData](#) ()
Initializes a new instance of the [GazeValidationData](#) class.
- [GazeValidationData](#) (float accuracyLeft, float accuracyRight, float precisionLeft, float precisionRight, float precisionRmsLeft, float precisionRmsRight)
Initializes a new instance of the [GazeValidationData](#) class.
- void [AddPoint](#) (Vector2 point, float accuracyLeft, float accuracyRight, float precisionLeft, float precisionRight, float precisionRmsLeft, float precisionRmsRight)
Add a new validation point to the list.

Properties

- float [AccuracyLeft](#) [get]
The accuracy in degrees averaged over all collected points for the left eye.
- float [AccuracyRight](#) [get]
The accuracy in degrees averaged over all collected points for the right eye.
- float [PrecisionLeft](#) [get]
The precision (standard deviation) in degrees averaged over all collected points for the left eye.
- float [PrecisionRight](#) [get]
The precision (standard deviation) in degrees averaged over all collected points for the right eye.
- float [PrecisionRmsLeft](#) [get]
The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the left eye.
- float [PrecisionRmsRight](#) [get]
The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the right eye.
- List< [GazeValidationPoint](#) > [Points](#) [get]
The list of all

8.42.1 Detailed Description

The gaze validation data structure

8.42.2 Constructor & Destructor Documentation

8.42.2.1 GazeValidationData() [1/2]

```
GazeUtilityLibrary.DataStructs.GazeValidationData.GazeValidationData ( ) [inline]
```

Initializes a new instance of the [GazeValidationData](#) class.

8.42.2.2 GazeValidationData() [2/2]

```
GazeUtilityLibrary.DataStructs.GazeValidationData.GazeValidationData (
    float accuracyLeft,
    float accuracyRight,
    float precisionLeft,
    float precisionRight,
    float precisionRmsLeft,
    float precisionRmsRight ) [inline]
```

Initializes a new instance of the [GazeValidationData](#) class.

Parameters

<i>accuracyLeft</i>	The accuracy in degrees averaged over all collected points for the left eye.
<i>accuracyRight</i>	The accuracy in degrees averaged over all collected points for the right eye.
<i>precisionLeft</i>	The precision (standard deviation) in degrees averaged over all collected points for the left eye.
<i>precisionRight</i>	The precision (standard deviation) in degrees averaged over all collected points for the right eye.
<i>precisionRmsLeft</i>	The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the left eye.
<i>precisionRmsRight</i>	The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the right eye.

8.42.3 Member Function Documentation

8.42.3.1 AddPoint()

```
void GazeUtilityLibrary.DataStructs.GazeValidationData.AddPoint (
    Vector2 point,
    float accuracyLeft,
    float accuracyRight,
    float precisionLeft,
    float precisionRight,
    float precisionRmsLeft,
    float precisionRmsRight ) [inline]
```

Add a new validation point to the list.

Parameters

<i>point</i>	The validation point coordinates.
<i>accuracyLeft</i>	The accuracy in degrees averaged over all collected points for the left eye.
<i>accuracyRight</i>	The accuracy in degrees averaged over all collected points for the right eye.
<i>precisionLeft</i>	The precision (standard deviation) in degrees averaged over all collected points for the left eye.
<i>precisionRight</i>	The precision (standard deviation) in degrees averaged over all collected points for the right eye.
<i>precisionRmsLeft</i>	The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the left eye.
<i>precisionRmsRight</i>	The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the right eye.

8.42.4 Property Documentation

8.42.4.1 AccuracyLeft

```
float GazeUtilityLibrary.DataStructs.GazeValidationData.AccuracyLeft [get]
```

The accuracy in degrees averaged over all collected points for the left eye.

8.42.4.2 AccuracyRight

```
float GazeUtilityLibrary.DataStructs.GazeValidationData.AccuracyRight [get]
```

The accuracy in degrees averaged over all collected points for the right eye.

8.42.4.3 Points

```
List<GazeValidationPoint> GazeUtilityLibrary.DataStructs.GazeValidationData.Points [get]
```

The list of all

8.42.4.4 PrecisionLeft

```
float GazeUtilityLibrary.DataStructs.GazeValidationData.PrecisionLeft [get]
```

The precision (standard deviation) in degrees averaged over all collected points for the left eye.

8.42.4.5 PrecisionRight

```
float GazeUtilityLibrary.DataStructs.GazeValidationData.PrecisionRight [get]
```

The precision (standard deviation) in degrees averaged over all collected points for the right eye.

8.42.4.6 PrecisionRmsLeft

```
float GazeUtilityLibrary.DataStructs.GazeValidationData.PrecisionRmsLeft [get]
```

The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the left eye.

8.42.4.7 PrecisionRmsRight

```
float GazeUtilityLibrary.DataStructs.GazeValidationData.PrecisionRmsRight [get]
```

The precision (root mean square of sample-to-sample error) in degrees averaged over all collected points for the right eye.

The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/DataStructs/GazeValidationData.cs`

8.43 GazeUtilityLibrary.DataStructs.GazeValidationPoint Class Reference

A validation point.

Public Member Functions

- [GazeValidationPoint](#) (Vector2 point, [GazeValidationData](#) result)
Initializes a new instance of the [GazeValidationPoint](#) class.
- string[] [Prepare](#) ([ConfigItem](#) config)
Prepare a list of formatted calibration data values

Properties

- Vector2 [Point](#) [get]
The validation point.
- [GazeValidationData Result](#) [get]
The validation result of this point.

8.43.1 Detailed Description

A validation point.

8.43.2 Constructor & Destructor Documentation

8.43.2.1 GazeValidationPoint()

```
GazeUtilityLibrary.DataStructs.GazeValidationPoint.GazeValidationPoint (
    Vector2 point,
    GazeValidationData result ) [inline]
```

Initializes a new instance of the [GazeValidationPoint](#) class.

Parameters

<i>point</i>	The validation point.
<i>result</i>	The validation result of this point.

8.43.3 Member Function Documentation

8.43.3.1 Prepare()

```
string [ ] GazeUtilityLibrary.DataStructs.GazeValidationPoint.Prepare (
    ConfigItem config ) [inline]
```

Prepare a list of formatted calibration data values

Parameters

<i>config</i>	The gaze configuration structure
---------------	----------------------------------

Returns

A list of formatted values. Each index corresponds to a specific value. This allows to reorder the list according to a format string.

8.43.4 Property Documentation

8.43.4.1 Point

```
Vector2 GazeUtilityLibrary.DataStructs.GazeValidationPoint.Point [get]
```

The validation point.

8.43.4.2 Result

```
GazeValidationData GazeUtilityLibrary.DataStructs.GazeValidationPoint.Result [get]
```

The validation result of this point.


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/DataStructs/GazeValidationData.cs


8.44 CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter Class Reference

Converts True to Hidden and False to Visible

Inheritance diagram for CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter:

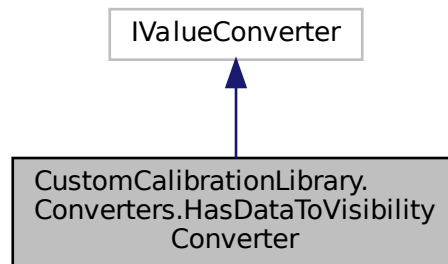

Collaboration diagram for CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter:

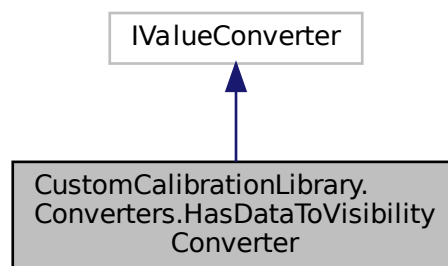

Public Member Functions

- object [Convert](#) (object value, Type targetType, object parameter, System.Globalization.CultureInfo culture)
Value converter.
- object [ConvertBack](#) (object value, Type targetType, object parameter, System.Globalization.CultureInfo culture)
Reverted value converter.

8.44.1 Detailed Description

Converts True to Hidden and False to Visible

8.44.2 Member Function Documentation

8.44.2.1 Convert()

```
object CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter.Convert (
    object value,
    Type targetType,
    object parameter,
    System.Globalization.CultureInfo culture ) [inline]
```

Value converter.

Parameters

<i>value</i>	The value to convert.
<i>targetType</i>	The type of the target value.
<i>parameter</i>	The conversion parameter.
<i>culture</i>	The language localisation.

Returns

The converted value object

Exceptions

<i>InvalidOperationException</i>	
----------------------------------	--

8.44.2.2 ConvertBack()

```
object CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter.ConvertBack (
    object value,
    Type targetType,
    object parameter,
    System.Globalization.CultureInfo culture ) [inline]
```

Reverted value converter.

Parameters

<i>value</i>	The value to convert.
<i>targetType</i>	The type of the target value.
<i>parameter</i>	The conversion parameter.
<i>culture</i>	The language localisation.

Returns

The converted value object

Exceptions

<i>NotSupportedException</i>	
------------------------------	--

The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Converters/HasDataToVisibilityConverter.cs

8.45 GazeUtilityLibrary.JsonConfigParser Class Reference

The config file "config.json" is parsed and its values are attributed to the [ConfigItem](#) class.

Public Member Functions

- [JsonConfigParser](#) ([TrackerLogger](#) logger)
Initializes a new instance of the [JsonConfigParser](#) class.
- [ConfigItem?](#) [ParseJsonConfig](#) (ref [GazeConfigError](#) error)
Parses the json configuration.
- void [SerializeJsonConfig](#) ([ConfigItem](#) item, string path)
Serializes the json configuration object to a string and writes it to a file.
- [ConfigItem](#) [GetDefaultConfig](#) ()
Gets the default configuration values.

8.45.1 Detailed Description

The config file "config.json" is parsed and its values are attributed to the [ConfigItem](#) class.

8.45.2 Constructor & Destructor Documentation

8.45.2.1 JsonConfigParser()

```
GazeUtilityLibrary.JsonConfigParser.JsonConfigParser (
    TrackerLogger logger ) [inline]
```

Initializes a new instance of the [JsonConfigParser](#) class.

Parameters

<i>logger</i>	The logger.
---------------	-------------

8.45.3 Member Function Documentation

8.45.3.1 GetDefaultConfig()

```
ConfigItem GazeUtilityLibrary.JsonConfigParser.GetDefaultConfig ( ) [inline]
```

Gets the default configuration values.

Returns

the default configuration values.

8.45.3.2 ParseJsonConfig()

```
ConfigItem? GazeUtilityLibrary.JsonConfigParser.ParseJsonConfig (
    ref GazeConfigError error ) [inline]
```

Parses the json configuration.

Returns

the updated [ConfigItem](#) class.

8.45.3.3 SerializeJsonConfig()

```
void GazeUtilityLibrary.JsonConfigParser.SerializeJsonConfig (
    ConfigItem item,
    string path ) [inline]
```

Serializes the json configuration object to a string and writes it to a file.

Parameters

<i>item</i>	The json configuration item.
<i>path</i>	The path where the file will be written.


The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/GazeConfiguration.cs`


8.46 GazeUtilityLibrary.DataStructs.LiveGazePoint Class Reference

The live gaze point used for verification during the calibration process.

Inheritance diagram for GazeUtilityLibrary.DataStructs.LiveGazePoint:

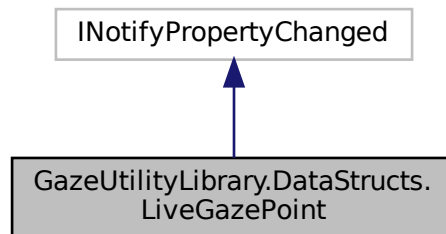

Collaboration diagram for GazeUtilityLibrary.DataStructs.LiveGazePoint:

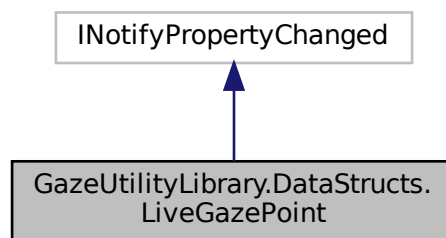

Properties

- double `X` [get, set]
The normalized x coordinate on the screen
- double `Y` [get, set]
The normalized y coordinate on the screen
- bool `Visibility` [get, set]
The visibility flag.

Events

- PropertyChangedEventHandler? `PropertyChanged`
Event to trigger property changes.

8.46.1 Detailed Description

The live gaze point used for verification during the calibration process.

8.46.2 Property Documentation

8.46.2.1 Visibility

```
bool GazeUtilityLibrary.DataStructs.LiveGazePoint.Visibility [get], [set]
```

The visibility flag.

8.46.2.2 X

```
double GazeUtilityLibrary.DataStructs.LiveGazePoint.X [get], [set]
```

The normalized x coordinate on the screen

8.46.2.3 Y

```
double GazeUtilityLibrary.DataStructs.LiveGazePoint.Y [get], [set]
```

The normalized y coordinate on the screen

8.46.3 Event Documentation

8.46.3.1 PropertyChanged

```
PropertyChangedEventHandler? GazeUtilityLibrary.DataStructs.LiveGazePoint.PropertyChanged
```

Event to trigger property changes.

The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/DataStructs/LiveGazePoint.cs`

8.47 CustomCalibrationLibrary.ViewModels.Monitor Class Reference

A representation of the screen.

Public Member Functions

- [Monitor](#) (int index, string name)
Initializes a new instance of the [Monitor](#) class.

Properties

- string [Name](#) [get]
The name of the screen.
- int [Index](#) [get]
The screen index.

8.47.1 Detailed Description

A representation of the screen.

8.47.2 Constructor & Destructor Documentation

8.47.2.1 Monitor()

```
CustomCalibrationLibrary.ViewModels.Monitor.Monitor (  
    int index,  
    string name ) [inline]
```

Initializes a new instance of the [Monitor](#) class.

Parameters

<i>index</i>	The screen index.
<i>name</i>	The name of the screen.

8.47.3 Property Documentation

8.47.3.1 Index

```
int CustomCalibrationLibrary.ViewModels.Monitor.Index [get]
```

The screen index.

8.47.3.2 Name

```
string CustomCalibrationLibrary.ViewModels.Monitor.Name [get]
```

The name of the screen.

The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/ViewModels/ScreenSelectionViewModel.cs

8.48 GazeUtilityLibrary.MouseHider Class Reference

hide standard mouse pointer and restore it

Public Member Functions

- [MouseHider](#) ([TrackerLogger](#) logger)
Initializes a new instance of the [MouseHider](#) class.
- void [HideCursor](#) ()
Hides the cursor.
- void [ShowCursor](#) (string? pathToCur)
Shows the cursor.

8.48.1 Detailed Description

hide standard mouse pointer and restore it

8.48.2 Constructor & Destructor Documentation

8.48.2.1 MouseHider()

```
GazeUtilityLibrary.MouseHider.MouseHider (  
    TrackerLogger logger ) [inline]
```

Initializes a new instance of the [MouseHider](#) class.

Parameters

<i>logger</i>	The logger.
---------------	-------------

8.48.3 Member Function Documentation

8.48.3.1 HideCursor()

```
void GazeUtilityLibrary.MouseHider.HideCursor ( ) [inline]
```

Hides the cursor.

Hides the standard mouse pointer by replacing the current icon with a transparent icon.

8.48.3.2 ShowCursor()

```
void GazeUtilityLibrary.MouseHider.ShowCursor (
    string? pathToCur ) [inline]
```

Shows the cursor.

the standard mouse pointer by replacing the current icon with the standard mouse pointer icon

Parameters

<i>pathToCur</i>	The path to the standard mouse pointer icon.
------------------	--


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/MouseHider.cs


8.49 GazeUtilityLibrary.Tracker.MouseTracker Class Reference

This class is used to hook into the system mouse events and track the position

Inheritance diagram for GazeUtilityLibrary.Tracker.MouseTracker:

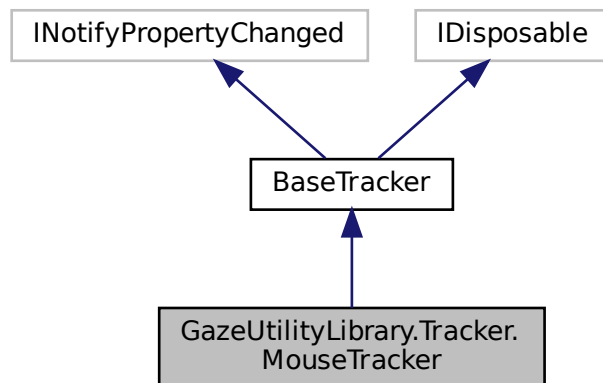

Collaboration diagram for GazeUtilityLibrary.Tracker.MouseTracker:

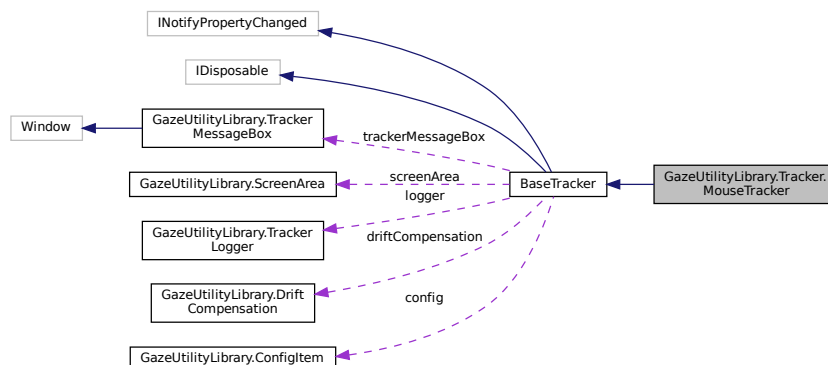

Public Member Functions

- `MouseTracker (TrackerLogger logger, ConfigItem config)`
Initializes a new instance of the `MouseTracker` class.
- `override Task< List< GazeCalibrationData > > ApplyCalibration ()`
Apply the calibration data. This is device specific and must be overwritten by the device class.
- `void Start ()`
Hooks the callback function `HookCallback(int, IntPtr, IntPtr)` to mouse events.
- `void Stop ()`
Removes to mouse event hook.
- `override Task InitCalibrationAsync ()`
Initialise the async calibration process. This is device specific and must be overwritten by the device class.
- `override void InitValidation ()`

- Initialise the validation process. This is device specific and must be overwritten by the device class.*
- override Task [FinishCalibrationAsync](#) ()
 - Finish the async calibration process. This is device specific and must be overwritten by the device class.*
- override void [FinishValidation](#) ()
 - Finish the validation process. This is device specific and must be overwritten by the device class.*
- override Task< bool > [CollectCalibrationDataAsync](#) (Point point)
 - Collect calibration data on a calibration point. This is device specific and must be overwritten by the device class.*
- override Task< bool > [CollectValidationDataAsync](#) (Point point)
 - Collect validation data on a validation point. This is device specific and must be overwritten by the device class.*
- override void [InitCalibration](#) ()
 - Initialise the calibration process. This is device specific and must be overwritten by the device class.*
- override void [FinishCalibration](#) ()
 - Finish the calibration process. This is device specific and must be overwritten by the device class.*
- override? [GazeValidationData ComputeValidation](#) ()
 - Apply the validation data. This is device specific and must be overwritten by the device class.*

Protected Member Functions

- override void [Dispose](#) (bool disposing)
 - Releases unmanaged and - optionally - managed resources.*
- override int [GetFixationFrameCount](#) (int durationThreshold)
 - Get the number of required gaze samples to compute a fixation. This is device specific and must be overwritten by the device because the duration of fixation point detection depends on the frame rate of the device.*
- override Vector3 [GetUnitDirection](#) ()
 - Get the unit vector pointing in the direction of the gaze vector. This is device specific as the gaze data are represented in a coordinate system as defined by the device.*
- override void [InitDriftCompensation](#) ()
 - Initialise the drift compensation. This is device specific and must be overwritten by the device class.*

Additional Inherited Members

8.49.1 Detailed Description

This class is used to hook into the system mouse events and track the position

See also

GazeHelper.TrackerHandler

8.49.2 Constructor & Destructor Documentation

8.49.2.1 MouseTracker()

```
GazeUtilityLibrary.Tracker.MouseTracker.MouseTracker (
    TrackerLogger logger,
    ConfigItem config ) [inline]
```

Initializes a new instance of the [MouseTracker](#) class.

Parameters

<i>logger</i>	The logger.
<i>config</i>	The config item.

8.49.3 Member Function Documentation

8.49.3.1 ApplyCalibration()

```
override Task<List<GazeCalibrationData> > GazeUtilityLibrary.Tracker.MouseTracker.ApplyCalibration ( ) [inline], [virtual]
```

Apply the calibration data. This is device specific and must be overwritten by the device class.

Returns

The calibration data result wrapped by an async handler.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.2 CollectCalibrationDataAsync()

```
override Task<bool> GazeUtilityLibrary.Tracker.MouseTracker.CollectCalibrationDataAsync ( Point point ) [inline], [virtual]
```

Collect calibration data on a calibration point. This is device specific and must be overwritten by the device class.

Parameters

<i>point</i>	The calibration point for which to collect data
--------------	---

Returns

True on success, false on failure, wrapped by an async handler.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.3 CollectValidationDataAsync()

```
override Task<bool> GazeUtilityLibrary.Tracker.MouseTracker.CollectValidationDataAsync ( Point point ) [inline], [virtual]
```

Collect validation data on a validation point. This is device specific and must be overwritten by the device class.

Parameters

<i>point</i>	The calibration point for which to collect data
--------------	---

Returns

True on success, false on failure, wrapped by an async handler.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.4 ComputeValidation()

```
override? GazeValidationData GazeUtilityLibrary.Tracker.MouseTracker.ComputeValidation ( )  
[inline], [virtual]
```

Apply the validation data. This is device specific and must be overwritten by the device class.

Returns

The validation data result.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.5 Dispose()

```
override void GazeUtilityLibrary.Tracker.MouseTracker.Dispose (  
    bool disposing ) [inline], [protected], [virtual]
```

Releases unmanaged and - optionally - managed resources.

Parameters

<i>disposing</i>	true to release both managed and unmanaged resources; false to release only unmanaged resources.
------------------	--

Reimplemented from [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.6 FinishCalibration()

```
override void GazeUtilityLibrary.Tracker.MouseTracker.FinishCalibration ( ) [inline], [virtual]
```

Finish the calibration process. This is device specific and must be overwritten by the device class.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.7 FinishCalibrationAsync()

```
override Task GazeUtilityLibrary.Tracker.MouseTracker.FinishCalibrationAsync ( ) [inline],  
[virtual]
```

Finish the async calibration process. This is device specific and must be overwritten by the device class.

Returns

An async handler

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.8 FinishValidation()

```
override void GazeUtilityLibrary.Tracker.MouseTracker.FinishValidation ( ) [inline], [virtual]
```

Finish the validation process. This is device specific and must be overwritten by the device class.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.9 GetFixationFrameCount()

```
override int GazeUtilityLibrary.Tracker.MouseTracker.GetFixationFrameCount (   
    int durationThreshold ) [inline], [protected], [virtual]
```

Get the number of required gaze samples to compute a fixation. This is device specific and must be overwritten by the device because the duration of fixation point detection depends on the frame rate of the device.

Parameters

<i>durationThreshold</i>	The required fixation duration in milliseconds.
--------------------------	---

Returns

The number of gaze samples to require for fixation detection.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.10 GetUnitDirection()

```
override Vector3 GazeUtilityLibrary.Tracker.MouseTracker.GetUnitDirection ( ) [inline], [protected],  
[virtual]
```

Get the unit vector pointing in the direction of the gaze vector. This is device specific as the gaze data are represented in a coordinate system as defined by the device.

Returns

The unit vector

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.11 InitCalibration()

```
override void GazeUtilityLibrary.Tracker.MouseTracker.InitCalibration ( ) [inline], [virtual]
```

Initialise the calibration process. This is device specific and must be overwritten by the device class.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.12 InitCalibrationAsync()

```
override Task GazeUtilityLibrary.Tracker.MouseTracker.InitCalibrationAsync ( ) [inline],  
[virtual]
```

Initialise the async calibration process. This is device specific and must be overwritten by the device class.

Returns

An async handler

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.13 InitDriftCompensation()

```
override void GazeUtilityLibrary.Tracker.MouseTracker.InitDriftCompensation ( ) [inline],  
[protected], [virtual]
```

Initialise the drift compensation. This is device specific and must be overwritten by the device class.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.14 InitValidation()

```
override void GazeUtilityLibrary.Tracker.MouseTracker.InitValidation ( ) [inline], [virtual]
```

Initialise the validation process. This is device specific and must be overwritten by the device class.

Implements [GazeUtilityLibrary.Tracker.BaseTracker](#).

8.49.3.15 Start()

```
void GazeUtilityLibrary.Tracker.MouseTracker.Start ( ) [inline]
```

Hooks the callback function HookCallback(int, IntPtr, IntPtr) to mouse events.

8.49.3.16 Stop()

```
void GazeUtilityLibrary.Tracker.MouseTracker.Stop ( ) [inline]
```

Removes to mouse event hook.

The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/Tracker/MouseTracker.cs

8.50 GazeControlLibrary.PipeCommand Class Reference

The JSON structure of a pipe command.

Public Member Functions

- [PipeCommand](#) (string command, bool reset, int? trialId, string label)
Initializes a new instance of the [PipeCommand](#) class.

Properties

- string [Command](#) [get, set]
The optional pipe command to be sent.
- string [Label](#) [get, set]
An optional label to annotate gaze data.
- int? [TrialId](#) [get, set]
An optional trial ID to annotate gaze data.
- bool? [ResetStartTime](#) [get, set]
An optional flag to indicate whether the relative timestamp should be reset.

8.50.1 Detailed Description

The JSON structure of a pipe command.

8.50.2 Constructor & Destructor Documentation

8.50.2.1 PipeCommand()

```
GazeControlLibrary.PipeCommand.PipeCommand (
    string command,
    bool reset,
    int? trialId,
    string label ) [inline]
```

Initializes a new instance of the [PipeCommand](#) class.

Parameters

<i>command</i>	The pipe command to be sent.
<i>reset</i>	A flag to indicate whether the relative timestamp should be reset.
<i>trialId</i>	An optional trial ID to annotate gaze data.
<i>label</i>	An optional label to annotate gaze data.

8.50.3 Property Documentation

8.50.3.1 Command

```
string GazeControlLibrary.PipeCommand.Command [get], [set]
```

The optional pipe command to be sent.

8.50.3.2 Label

```
string GazeControlLibrary.PipeCommand.Label [get], [set]
```

An optional label to annotate gaze data.

8.50.3.3 ResetStartTime

```
bool? GazeControlLibrary.PipeCommand.ResetStartTime [get], [set]
```

An optional flag to indicate whether the relative timestamp should be reset.

8.50.3.4 TrialId

```
int? GazeControlLibrary.PipeCommand.TrialId [get], [set]
```

An optional trial ID to annotate gaze data.


The documentation for this class was generated from the following file:

- `source/GazeControlLibrary/PipeCommand.cs`


8.51 CustomCalibrationLibrary.Converters.PositionConverter Class Reference

Converter class to convert a normalized coordinate to a pixel coordinate.

Inheritance diagram for CustomCalibrationLibrary.Converters.PositionConverter:

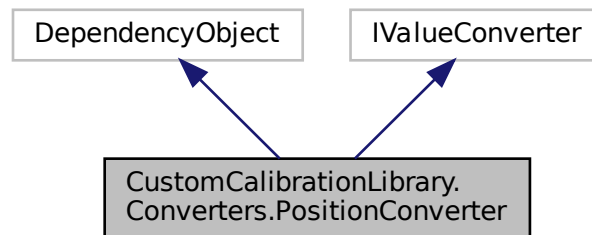

Collaboration diagram for CustomCalibrationLibrary.Converters.PositionConverter:

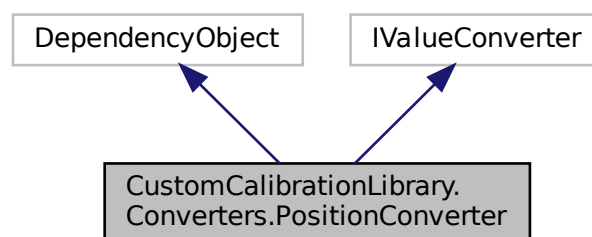

Public Member Functions

- object [Convert](#) (object value, Type targetType, object parameter, CultureInfo culture)
Value converter.
- object [ConvertBack](#) (object value, Type targetType, object parameter, CultureInfo culture)
Reverted value converter.

Static Public Attributes

- static readonly DependencyProperty [OffsetProperty](#)
The custom offset property of the value converter.

Properties

- string?? [Offset](#) [get, set]
The position offset.

8.51.1 Detailed Description

Converter class to convert a normalized coordinate to a pixel coordinate.

8.51.2 Member Function Documentation

8.51.2.1 Convert()

```
object CustomCalibrationLibrary.Converters.PositionConverter.Convert (
    object value,
    Type targetType,
    object parameter,
    CultureInfo culture ) [inline]
```

Value converter.

Parameters

<i>value</i>	The value to convert.
<i>targetType</i>	The type of the target value.
<i>parameter</i>	The conversion parameter.
<i>culture</i>	The language localisation.

Returns

The converted value object

8.51.2.2 ConvertBack()

```
object CustomCalibrationLibrary.Converters.PositionConverter.ConvertBack (
    object value,
    Type targetType,
    object parameter,
    CultureInfo culture ) [inline]
```

Reverted value converter.

Parameters

<i>value</i>	The value to convert.
<i>targetType</i>	The type of the target value.
<i>parameter</i>	The conversion parameter.
<i>culture</i>	The language localisation.

Returns

The converted value object

Exceptions

<i>NotSupportedException</i>	
------------------------------	--

8.51.3 Member Data Documentation**8.51.3.1 OffsetProperty**

```
readonly DependencyProperty CustomCalibrationLibrary.Converters.PositionConverter.Offset↔
Property [static]
```

Initial value:

```
=
    DependencyProperty.Register("Offset", typeof(string), typeof(PositionConverter), new
        PropertyMetadata(null))
```

The custom offset property of the value converter.

8.51.4 Property Documentation**8.51.4.1 Offset**

```
string?? CustomCalibrationLibrary.Converters.PositionConverter.Offset [get], [set]
```

The position offset.


The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Converters/PositionConverter.cs


8.52 CustomCalibrationLibrary.Converters.ProximityColorConverter Class Reference

Converter class to convert the proximito of a normailezed coordinate to the center point (0.5) into colors.

Inheritance diagram for CustomCalibrationLibrary.Converters.ProximityColorConverter:

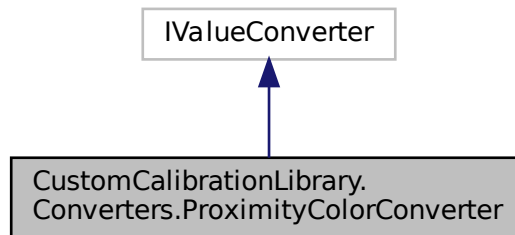

Collaboration diagram for CustomCalibrationLibrary.Converters.ProximityColorConverter:

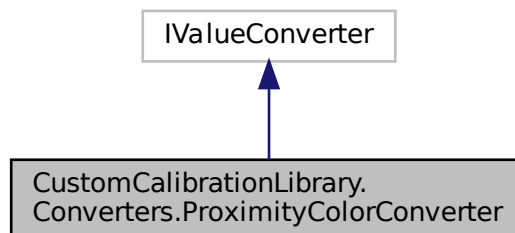

Public Member Functions

- object [Convert](#) (object value, Type targetType, object parameter, CultureInfo culture)
Value converter.
- object [ConvertBack](#) (object value, Type targetType, object parameter, CultureInfo culture)
Reverted value converter.

8.52.1 Detailed Description

Converter class to convert the proximito of a normailezed coordinate to the center point (0.5) into colors.

8.52.2 Member Function Documentation

8.52.2.1 Convert()

```
object CustomCalibrationLibrary.Converters.ProximityColorConverter.Convert (
    object value,
    Type targetType,
    object parameter,
    CultureInfo culture ) [inline]
```

Value converter.

Parameters

<i>value</i>	The value to convert.
<i>targetType</i>	The type of the target value.
<i>parameter</i>	The conversion parameter.
<i>culture</i>	The language localisation.

Returns

The converted value object

8.52.2.2 ConvertBack()

```
object CustomCalibrationLibrary.Converters.ProximityColorConverter.ConvertBack (
    object value,
    Type targetType,
    object parameter,
    CultureInfo culture ) [inline]
```

Reverted value converter.

Parameters

<i>value</i>	The value to convert.
<i>targetType</i>	The type of the target value.
<i>parameter</i>	The conversion parameter.
<i>culture</i>	The language localisation.

Returns

The converted value object

Exceptions

<code>NotSupportedException</code>	
------------------------------------	--

The documentation for this class was generated from the following file:

- `source/CustomCalibrationLibrary/Converters/ProximityColorConverter.cs`

8.53 GazeUtilityLibrary.ScreenArea Class Reference

The class describing the Screen area in 3d and 2d space.

Public Member Functions

- [ScreenArea](#) (Vector3 bottomLeft, Vector3 bottomRight, Vector3 topLeft, Vector3 topRight, float width, float height)
Constructor. Assigns parameters and computes the transformation matrix to transform a 3d point into a 2d point.
- Vector3? [GetIntersectionPoint](#) (Vector3 gazeOrigin, Vector3 gazeDirection)
Compute the intersection point with the screen plane given a gaze origin and a gaze direction. Note that this does not compute the intersection with the screen area but with the infinite plane which is co-aligned with the screen. Pass the here computed intersection point to the method [GetPoint2dNormalized](#) to get the normalized intersection point on the screen area.
- Vector2 [GetPoint2d](#) (Vector3 point)
Get the 2d point on the screen given a 3d point on the screen plane.
- Vector2 [GetPoint2dNormalized](#) (Vector3 point3d)
Get the normalized 2d point on the screen given a 3d point on the screen plane. Note that values outside of the interval [0, 1] indicate an intersection point outside of the screen area.
- bool [Dump](#) (string path, string prefix)
Dump the four screen corner points to a csv file

Properties

- float [Width](#) [get]
The width of the screen.
- float [Height](#) [get]
The height of the screen.
- Vector3 [BottomLeft](#) [get]
The coordinates of the bottom left point of the screen.
- Vector3 [BottomRight](#) [get]
The coordinates of the bottom right point of the screen.
- Vector3 [TopLeft](#) [get]
The coordinates of the top left point of the screen.
- Vector3 [TopRight](#) [get]
The coordinates of the top right point of the screen.
- Vector3 [Center](#) [get]
The coordinates of the center point of the screen.

8.53.1 Detailed Description

The class describing the Screen area in 3d and 2d space.

8.53.2 Constructor & Destructor Documentation

8.53.2.1 ScreenArea()

```
GazeUtilityLibrary.ScreenArea.ScreenArea (
    Vector3 bottomLeft,
    Vector3 bottomRight,
    Vector3 topLeft,
    Vector3 topRight,
    float width,
    float height ) [inline]
```

Constructor. Assigns parameters and computes the transformation matrix to transform a 3d point into a 2d point.

Parameters

<i>bottomLeft</i>	The bottom left 3d coordinate of the screen.
<i>bottomRight</i>	The bottom right 3d coordinate of the screen.
<i>topLeft</i>	The top left 3d coordinate of the screen.
<i>topRight</i>	The top right 3d coordinate of the screen
<i>width</i>	The width of the screen
<i>height</i>	The height of the screen

8.53.3 Member Function Documentation

8.53.3.1 Dump()

```
bool GazeUtilityLibrary.ScreenArea.Dump (
    string path,
    string prefix ) [inline]
```

Dump the four screen corner points to a csv file

Parameters

<i>path</i>	The folder to store the file.
<i>prefix</i>	The file prefix.

Returns

8.53.3.2 GetIntersectionPoint()

```
Vector3? GazeUtilityLibrary.ScreenArea.GetIntersectionPoint (
    Vector3 gazeOrigin,
    Vector3 gazeDirection ) [inline]
```

Compute the intersection point with the screen plane given a gaze origin and a gaze direction. Note that this does not compute the intersection with the screen area but with the infinite plane which is co-aligned with the screen. Pass the here computed intersection point to the method GetPoint2dNormalized to get the normalized intersection point on the screen area.

Parameters

<i>gazeOrigin</i>	The origin of the gaze.
<i>gazeDirection</i>	The direction of the gaze.

Returns

The intersection point with the screen or null if no intersection point exists.

8.53.3.3 GetPoint2d()

```
Vector2 GazeUtilityLibrary.ScreenArea.GetPoint2d (
    Vector3 point ) [inline]
```

Get the 2d point on the screen given given a 3d point on the screen plane.

Parameters

<i>point</i>	The 3d point on the screen plane to convert.
--------------	--

Returns

The 2d point on the screen plane

8.53.3.4 GetPoint2dNormalized()

```
Vector2 GazeUtilityLibrary.ScreenArea.GetPoint2dNormalized (
    Vector3 point3d ) [inline]
```

Get the normalized 2d point on the screen given a 3d point on the screen plane. Note that values outside of the interval $[0, 1]$ indicate an intersection point outside of the screen area.

Parameters

<i>point3d</i>	The 3d point on the screen plane to convert.
----------------	--

Returns

The normalized 2d point on the screen plane

8.53.4 Property Documentation

8.53.4.1 BottomLeft

`Vector3 GazeUtilityLibrary.ScreenArea.BottomLeft` [get]

The coordinates of the bottom left point of the screen.

8.53.4.2 BottomRight

`Vector3 GazeUtilityLibrary.ScreenArea.BottomRight` [get]

The coordinates of the bottom right point of the screen.

8.53.4.3 Center

`Vector3 GazeUtilityLibrary.ScreenArea.Center` [get]

The coordinates of the center point of the screen.

8.53.4.4 Height

`float GazeUtilityLibrary.ScreenArea.Height` [get]

The height of the screen.

8.53.4.5 TopLeft

```
Vector3 GazeUtilityLibrary.ScreenArea.TopLeft [get]
```

The coordinates of the top left point of the screen.

8.53.4.6 TopRight

```
Vector3 GazeUtilityLibrary.ScreenArea.TopRight [get]
```

The coordinates of the to right point of the screen.

8.53.4.7 Width

```
float GazeUtilityLibrary.ScreenArea.Width [get]
```

The width of the screen.


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/ScreenArea.cs


8.54 Tobii.Research.Addons.ScreenBasedCalibrationValidation Class Reference

Provides methods and properties for managing calibration validation for screen based eye trackers.

Inheritance diagram for Tobii.Research.Addons.ScreenBasedCalibrationValidation:

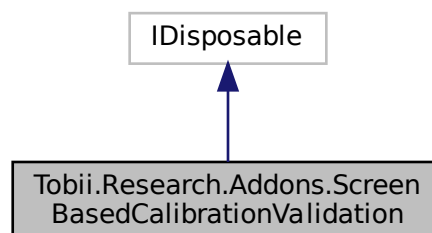

Collaboration diagram for Tobii.Research.Addons.ScreenBasedCalibrationValidation:

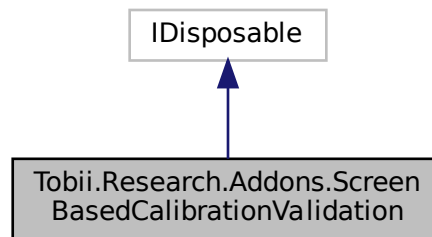

Public Types

- enum [ValidationState](#) { **NotInValidationMode**, **NotCollectingData**, **CollectingData** }

ValidationState.NotInValidationMode - [EnterValidationMode](#) must be called starting to collect data. *ValidationState.*↔
NotCollectingData - Ready to start collecting data or computing result. *ValidationState.CollectingData* - Currently collecting data. Will finish after the sample count is reached or a timeout.

Public Member Functions

- [ScreenBasedCalibrationValidation](#) (IEyeTracker eyeTracker, int sampleCount=30, int timeoutMS=1000)
 Create a calibration validation object for screen based eye trackers.
- void [StartCollectingData](#) (NormalizedPoint2D calibrationPointCoordinates)
 Starts collecting data for a calibration validation point. The argument used is the point the user is assumed to be looking at and is given in the active display area coordinate system. Please check State property to know when data collection is completed (or timed out).
- void [DiscardData](#) (NormalizedPoint2D calibrationPointCoordinates)
 Removes the collected data for a specific calibration validation point.
- void [EnterValidationMode](#) ()
 Enter the calibration validation mode and starts subscribing to gaze data from the eye tracker.
- void [LeaveValidationMode](#) ()
 Leaves the calibration validation mode, clears all collected data, and unsubscribes from the eye tracker.
- [CalibrationValidationResult Compute](#) ()
 Uses the collected data and tries to compute accuracy and precision values for all points. If the calculation is successful, the result is returned, and stored in the Result property of the CalibrationValidation object. If there is insufficient data to compute the results for a certain point that [CalibrationValidationPoint](#) will contain invalid data (NaN) for the results. Gaze data will still be untouched. If there is no valid data for any point, the average results of [CalibrationValidationResult](#) will be invalid (NaN) as well.
- void [Dispose](#) ()
 Dispose will unsubscribe to gaze data and exit validation mode, if the object is not already in *ValidationState.NotInValidationMode*↔
- override string [ToString](#) ()
 Convert validation values to a string.

Properties

- [ValidationState State](#) [get]
Get the current state of the validation object.
- [CalibrationValidationResult Result](#) [get]
Get the current [CalibrationValidationResult](#) with the computed accuracy and precision. [Compute](#) must have been called for this to contain valid data.

8.54.1 Detailed Description

Provides methods and properties for managing calibration validation for screen based eye trackers.

8.54.2 Member Enumeration Documentation

8.54.2.1 ValidationState

```
enum Tobii.Research.Addons.ScreenBasedCalibrationValidation.ValidationState [strong]
```

[ValidationState.NotInValidationMode](#) - [EnterValidationMode](#) must be called starting to collect data. [ValidationState.NotCollectingData](#) - Ready to start collecting data or computing result. [ValidationState.CollectingData](#) - Currently collecting data. Will finish after the sample count is reached or a timeout.

8.54.3 Constructor & Destructor Documentation

8.54.3.1 ScreenBasedCalibrationValidation()

```
Tobii.Research.Addons.ScreenBasedCalibrationValidation.ScreenBasedCalibrationValidation (  
    IEyeTracker eyeTracker,  
    int sampleCount = 30,  
    int timeoutMS = 1000 ) [inline]
```

Create a calibration validation object for screen based eye trackers.

Parameters

<i>eyeTracker</i>	An IEyeTracker instance.
<i>sampleCount</i>	The number of samples to collect. Default 30, minimum 10, maximum 3000.
<i>timeoutMS</i>	Timeout in milliseconds. Default 1000, minimum 100, maximum 3000.

8.54.4 Member Function Documentation

8.54.4.1 Compute()

```
CalibrationValidationResult Tobii.Research.Addons.ScreenBasedCalibrationValidation.Compute ( )
[inline]
```

Uses the collected data and tries to compute accuracy and precision values for all points. If the calculation is successful, the result is returned, and stored in the Result property of the CalibrationValidation object. If there is insufficient data to compute the results for a certain point that [CalibrationValidationPoint](#) will contain invalid data (NaN) for the results. Gaze data will still be untouched. If there is no valid data for any point, the average results of [CalibrationValidationResult](#) will be invalid (NaN) as well.

Returns

The [CalibrationValidationResult](#)

8.54.4.2 DiscardData()

```
void Tobii.Research.Addons.ScreenBasedCalibrationValidation.DiscardData (
    NormalizedPoint2D calibrationPointCoordinates ) [inline]
```

Removes the collected data for a specific calibration validation point.

Parameters

<i>calibrationPointCoordinates</i>	The calibration point to remove.
------------------------------------	----------------------------------

8.54.4.3 Dispose()

```
void Tobii.Research.Addons.ScreenBasedCalibrationValidation.Dispose ( ) [inline]
```

Dispose will unsubscribe to gaze data and exit validation mode, if the object is not already in ValidationState.NotInValidationMode

8.54.4.4 EnterValidationMode()

```
void Tobii.Research.Addons.ScreenBasedCalibrationValidation.EnterValidationMode ( ) [inline]
```

Enter the calibration validation mode and starts subscribing to gaze data from the eye tracker.

8.54.4.5 LeaveValidationMode()

```
void Tobii.Research.Addons.ScreenBasedCalibrationValidation.LeaveValidationMode ( ) [inline]
```

Leaves the calibration validation mode, clears all collected data, and unsubscribes from the eye tracker.

8.54.4.6 StartCollectingData()

```
void Tobii.Research.Addons.ScreenBasedCalibrationValidation.StartCollectingData (
    NormalizedPoint2D calibrationPointCoordinates ) [inline]
```

Starts collecting data for a calibration validation point. The argument used is the point the user is assumed to be looking at and is given in the active display area coordinate system. Please check State property to know when data collection is completed (or timed out).

Parameters

<i>calibrationPointCoordinates</i>	The normalized 2D point on the display area
------------------------------------	---

8.54.4.7 ToString()

```
override string Tobii.Research.Addons.ScreenBasedCalibrationValidation.ToString ( ) [inline]
```

Convert validation values to a string.

Returns

The validation string.

8.54.5 Property Documentation

8.54.5.1 Result

```
CalibrationValidationResult Tobii.Research.Addons.ScreenBasedCalibrationValidation.Result
[get]
```

Get the current [CalibrationValidationResult](#) with the computed accuracy and precision. [Compute](#) must have been called for this to contain valid data.

8.54.5.2 State

`ValidationState` Tobii.Research.Addons.ScreenBasedCalibrationValidation.State [get]

Get the current state of the validation object.


The documentation for this class was generated from the following file:

- source/TobiiProSdkAddons/ScreenBasedCalibrationValidation.cs


8.55 CustomCalibrationLibrary.Views.ScreenSelection Class Reference

Interaction logic for ScreenSelection.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.ScreenSelection:

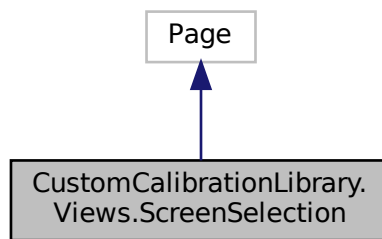

Collaboration diagram for CustomCalibrationLibrary.Views.ScreenSelection:

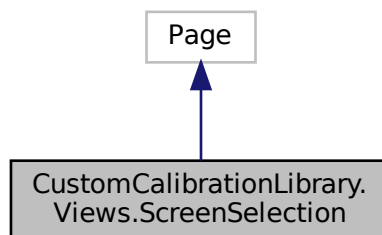

Public Member Functions

- `ScreenSelection` (`CalibrationModel` model, Window window)
Initializes a new instance of the `ScreenSelection` class.

8.55.1 Detailed Description

Interaction logic for ScreenSelection.xaml

8.55.2 Constructor & Destructor Documentation

8.55.2.1 ScreenSelection()

```
CustomCalibrationLibrary.Views.ScreenSelection.ScreenSelection (
    CalibrationModel model,
    Window window ) [inline]
```

Initializes a new instance of the [ScreenSelection](#) class.

Parameters

<i>model</i>	The calibration model.
<i>window</i>	The target window.

The documentation for this class was generated from the following file:

- `source/CustomCalibrationLibrary/Views/ScreenSelection.xaml.cs`

8.56 CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel Class Reference

The view model class for the screen selection view.

Public Member Functions

- [ScreenSelectionViewModel](#) ([CalibrationModel](#) model, [Window](#) window)
Initializes a new instance of the [ScreenSelectionViewModel](#) class.

Properties

- [ObservableCollection](#)< [Monitor](#) > [Monitors](#) [get]
The observable list of monitors to select from.
- [ICommand](#) [CalibrationStartCommand](#) [get]
Command to start the calibration
- [ICommand](#) [CalibrationAbortCommand](#) [get]
Command to abort the calibration
- [ICommand](#) [ScreenSwitchCommand](#) [get]
Command to switch the screen

8.56.1 Detailed Description

The view model class for the screen selection view.

8.56.2 Constructor & Destructor Documentation

8.56.2.1 ScreenSelectionViewModel()

```
CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel.ScreenSelectionViewModel (
    CalibrationModel model,
    Window window ) [inline]
```

Initializes a new instance of the [ScreenSelectionViewModel](#) class.

Parameters

<i>model</i>	The calibration model
<i>window</i>	The target window of the screen selection

8.56.3 Property Documentation

8.56.3.1 CalibrationAbortCommand

```
ICommand CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel.CalibrationAbortCommand
[get]
```

Command to abort the calibration

8.56.3.2 CalibrationStartCommand

```
ICommand CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel.CalibrationStartCommand
[get]
```

Command to start the calibration

8.56.3.3 Monitors

```
ObservableCollection<Monitor> CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel.↔
Monitors [get]
```

The observable list of monitors to select from.

8.56.3.4 ScreenSwitchCommand

```
 ICommand CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel.ScreenSwitchCommand
[get]
```

Command to switch the screen

The documentation for this class was generated from the following file:

- `source/CustomCalibrationLibrary/ViewModels/ScreenSelectionViewModel.cs`

8.57 GazeUtilityLibrary.ScreenTriangle Class Reference

A class to describe a triangle. This was supposed to be used to construct the [ScreenArea](#) but it turned out that it is simpler to work with the screen plane and use the normalised intersection points to check whether the gaze point is outside the screen area.

Public Member Functions

- [ScreenTriangle](#) (Vector3 v1, Vector3 v2, Vector3 v3)
Initializes a new instance of the [ScreenTriangle](#) class.
- Vector3? [GetIntersectionPoint](#) (Vector3 origin, Vector3 direction)
Compute the intersection point with the triangle with the Moller-Trumbore algorithm.

Properties

- Vector3 [V1](#) [get]
A corner point of the triangle.
- Vector3 [V2](#) [get]
A corner point of the triangle.
- Vector3 [V3](#) [get]
A corner point of the triangle.
- Vector3 [E1](#) [get]
The edge vector from v1 to v2.
- Vector3 [E2](#) [get]
The edge vector from v1 to v3.

8.57.1 Detailed Description

A class to describe a triangle. This was supposed to be used to construct the [ScreenArea](#) but it turned out that it is simpler to work with the screen plane and use the normalised intersection points to check wheter the gaze point is outside the screen area.

8.57.2 Constructor & Destructor Documentation

8.57.2.1 ScreenTriangle()

```
GazeUtilityLibrary.ScreenTriangle.ScreenTriangle (
    Vector3 v1,
    Vector3 v2,
    Vector3 v3 ) [inline]
```

Initializes a new instance of the [ScreenTriangle](#) class.

Parameters

<i>v1</i>	A corner point of the triangle.
<i>v2</i>	A corner point of the triangle.
<i>v3</i>	A corner point of the triangle.

8.57.3 Member Function Documentation

8.57.3.1 GetIntersectionPoint()

```
Vector3? GazeUtilityLibrary.ScreenTriangle.GetIntersectionPoint (
    Vector3 origin,
    Vector3 direction ) [inline]
```

Compute the intersection point with the triangle with the Moller-Trumbore algorithm.

Parameters

<i>origin</i>	The origin of the gaze point
<i>direction</i>	The direction of the gaze point

Returns

The intersection point or null if no intersection point could be computed.

8.57.4 Property Documentation

8.57.4.1 E1

`Vector3 GazeUtilityLibrary.ScreenTriangle.E1 [get]`

The edge vector from v1 to v2.

8.57.4.2 E2

`Vector3 GazeUtilityLibrary.ScreenTriangle.E2 [get]`

The edge vector from v1 to v3.

8.57.4.3 V1

`Vector3 GazeUtilityLibrary.ScreenTriangle.V1 [get]`

A corner point of the triangle.

8.57.4.4 V2

`Vector3 GazeUtilityLibrary.ScreenTriangle.V2 [get]`

A corner point of the triangle.

8.57.4.5 V3

`Vector3 GazeUtilityLibrary.ScreenTriangle.V3 [get]`

A corner point of the triangle.

The documentation for this class was generated from the following file:

- `source/GazeUtilityLibrary/ScreenTriangle.cs`

8.58 GazeUtilityLibrary.TrackerLogger Class Reference

Simple logger class.

Public Member Functions

- [TrackerLogger](#) (string? logPath, [EOutputType](#) type=EOutputType.gaze)
Initializes a new instance of the [TrackerLogger](#) class.
- void [DumpFatal](#) (Exception e)
Dumps exception to a new file if it is not possible to write to the main log file.
- void [Debug](#) (string message)
wrapper function for debug level logging.
- void [Info](#) (string message)
wrapper function for info level logging
- void [Warning](#) (string message)
wrapper function for warning level logging
- void [Error](#) (string message)
wrapper function for error level logging

8.58.1 Detailed Description

Simple logger class.

8.58.2 Constructor & Destructor Documentation

8.58.2.1 TrackerLogger()

```
GazeUtilityLibrary.TrackerLogger.TrackerLogger (  
    string? logPath,  
    EOutputType type = EOutputType.gaze ) [inline]
```

Initializes a new instance of the [TrackerLogger](#) class.

8.58.3 Member Function Documentation

8.58.3.1 Debug()

```
void GazeUtilityLibrary.TrackerLogger.Debug (  
    string message ) [inline]
```

wrapper function for debug level logging.

Parameters

<i>message</i>	The message.
----------------	--------------

8.58.3.2 DumpFatal()

```
void GazeUtilityLibrary.TrackerLogger.DumpFatal (
    Exception e ) [inline]
```

Dumps exception to a new file if it is not possible to write to the main log file.

Parameters

<i>e</i>	The exception.
----------	----------------

8.58.3.3 Error()

```
void GazeUtilityLibrary.TrackerLogger.Error (
    string message ) [inline]
```

wrapper function for error level logging

Parameters

<i>message</i>	The message.
----------------	--------------

8.58.3.4 Info()

```
void GazeUtilityLibrary.TrackerLogger.Info (
    string message ) [inline]
```

wrapper function for info level logging

Parameters

<i>message</i>	The message.
----------------	--------------

8.58.3.5 Warning()

```
void GazeUtilityLibrary.TrackerLogger.Warning (  
    string message ) [inline]
```

wrapper function for warning level logging

Parameters

<i>message</i>	The message.
----------------	--------------


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/Logger.cs


8.59 GazeUtilityLibrary.TrackerMessageBox Class Reference

Interaction logic for TrackerMessageBox.xaml

Inheritance diagram for GazeUtilityLibrary.TrackerMessageBox:

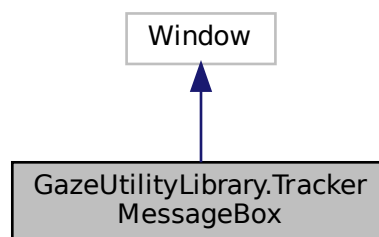

Collaboration diagram for GazeUtilityLibrary.TrackerMessageBox:

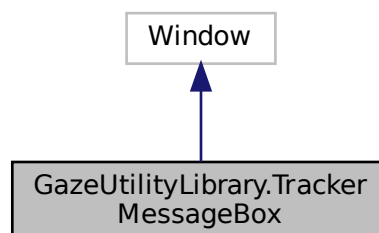

8.59.1 Detailed Description

Interaction logic for TrackerMessageBox.xaml


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/TrackerMessageBox.xaml.cs


8.60 GazeUtilityLibrary.DataStructs.UserPositionData Class Reference

The user position to be rendered on the screen.

Inheritance diagram for GazeUtilityLibrary.DataStructs.UserPositionData:

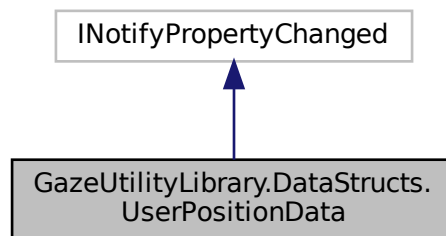

Collaboration diagram for GazeUtilityLibrary.DataStructs.UserPositionData:

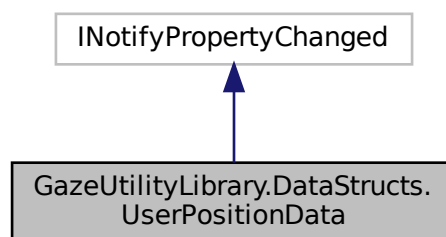

Public Member Functions

- [UserPositionData](#) ()
Initializes a new instance of the [UserPositionData](#) class.
- [UserPositionData](#) (double xCoordLeft, double yCoordLeft, double zCoordLeft, double xCoordRight, double yCoordRight, double zCoordRight)
Initializes a new instance of the [UserPositionData](#) class.

Properties

- double [XCoordLeft](#) [get, set]
The normalized x coordinate of the left eye.
- double [YCoordLeft](#) [get, set]
The normalized y coordinate of the left eye.
- double [ZCoordLeft](#) [get, set]
The normalized z coordinate of the left eye.
- double [XCoordRight](#) [get, set]
The normalized x coordinate of the right eye.
- double [YCoordRight](#) [get, set]
The normalized y coordinate of the right eye.
- double [ZCoordRight](#) [get, set]
The normalized z coordinate of the right eye.

Events

- PropertyChangedEventHandler? [PropertyChanged](#)
The property change event handler.

8.60.1 Detailed Description

The user position to be rendered on the screen.

8.60.2 Constructor & Destructor Documentation

8.60.2.1 [UserPositionData\(\)](#) [1/2]

```
GazeUtilityLibrary.DataStructs.UserPositionData.UserPositionData ( ) [inline]
```

Initializes a new instance of the [UserPositionData](#) class.

8.60.2.2 [UserPositionData\(\)](#) [2/2]

```
GazeUtilityLibrary.DataStructs.UserPositionData.UserPositionData (
    double xCoordLeft,
    double yCoordLeft,
    double zCoordLeft,
    double xCoordRight,
    double yCoordRight,
    double zCoordRight ) [inline]
```

Initializes a new instance of the [UserPositionData](#) class.

Parameters

<i>xCoordLeft</i>	The normalized x coordinate of the left eye.
<i>yCoordLeft</i>	The normalized y coordinate of the left eye.
<i>zCoordLeft</i>	The normalized z coordinate of the left eye.
<i>xCoordRight</i>	The normalized x coordinate of the right eye.
<i>yCoordRight</i>	The normalized y coordinate of the right eye.
<i>zCoordRight</i>	The normalized z coordinate of the right eye.

8.60.3 Property Documentation

8.60.3.1 XCoordLeft

```
double GazeUtilityLibrary.DataStructs.UserPositionData.XCoordLeft [get], [set]
```

The normalized x coordinate of the left eye.

8.60.3.2 XCoordRight

```
double GazeUtilityLibrary.DataStructs.UserPositionData.XCoordRight [get], [set]
```

The normalized x coordinate of the right eye.

8.60.3.3 YCoordLeft

```
double GazeUtilityLibrary.DataStructs.UserPositionData.YCoordLeft [get], [set]
```

The normalized y coordinate of the left eye.

8.60.3.4 YCoordRight

```
double GazeUtilityLibrary.DataStructs.UserPositionData.YCoordRight [get], [set]
```

The normalized y coordinate of the right eye.

8.60.3.5 ZCoordLeft

```
double GazeUtilityLibrary.DataStructs.UserPositionData.ZCoordLeft [get], [set]
```

The normalized z coordinate of the left eye.

8.60.3.6 ZCoordRight

```
double GazeUtilityLibrary.DataStructs.UserPositionData.ZCoordRight [get], [set]
```

The normalized z coordinate of the right eye.

8.60.4 Event Documentation

8.60.4.1 PropertyChanged

```
PropertyChangedEventHandler? GazeUtilityLibrary.DataStructs.UserPositionData.PropertyChanged
```

The property change event handler.


The documentation for this class was generated from the following file:

- source/GazeUtilityLibrary/DataStructs/UserPositionData.cs


8.61 CustomCalibrationLibrary.Views.UserPositionGuide Class Reference

Interaction logic for UserPositionGuide.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.UserPositionGuide:

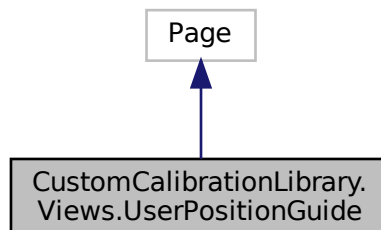

Collaboration diagram for CustomCalibrationLibrary.Views.UserPositionGuide:

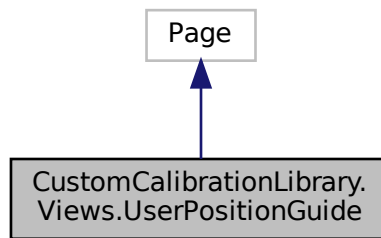

Public Member Functions

- [UserPositionGuide](#) (CalibrationModel model)
Initializes a new instance of the [UserPositionGuide](#) class.

8.61.1 Detailed Description

Interaction logic for UserPositionGuide.xaml

8.61.2 Constructor & Destructor Documentation

8.61.2.1 UserPositionGuide()

```
CustomCalibrationLibrary.Views.UserPositionGuide.UserPositionGuide (
    CalibrationModel model ) [inline]
```

Initializes a new instance of the [UserPositionGuide](#) class.

Parameters

<i>model</i>	The calibration model.
--------------	------------------------

The documentation for this class was generated from the following file:

- `source/CustomCalibrationLibrary/Views/UserPositionGuide.xaml.cs`

8.62 CustomCalibrationLibrary.ViewModels.UserPositionGuideView↵ Model Class Reference

The view model class for the user position guide view.

Public Member Functions

- [UserPositionGuideViewModel](#) (CalibrationModel model)
Constructor

Properties

- [UserPositionData](#) [UserPosition](#) [get]
The user position to be represented on the view
- ICommand [CalibrationStartCommand](#) [get]
Command to start the calibration
- ICommand [CalibrationAbortCommand](#) [get]
Command to abort the calibration

8.62.1 Detailed Description

The view model class for the user position guide view.

8.62.2 Constructor & Destructor Documentation

8.62.2.1 UserPositionGuideViewModel()

```
CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel.UserPositionGuideViewModel (
    CalibrationModel model ) [inline]
```

Constructor

Parameters

<i>model</i>	The calibartion model
--------------	-----------------------

8.62.3 Property Documentation

8.62.3.1 CalibrationAbortCommand

```
 ICommand CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel.CalibrationAbort↔  
 Command [get]
```

Command to abort the calibration

8.62.3.2 CalibrationStartCommand

```
 ICommand CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel.CalibrationStart↔  
 Command [get]
```

Command to start the calibration

8.62.3.3 UserPosition

```
 UserPositionData CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel.UserPosition  
 [get]
```

The user position to be represented on the view


The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/ViewModels/UserPositionGuideViewModel.cs


8.63 CustomCalibrationLibrary.Views.ValidationResult Class Reference

Interaction logic for ValidationResult.xaml

Inheritance diagram for CustomCalibrationLibrary.Views.ValidationResult:

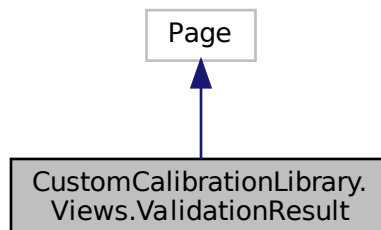

Collaboration diagram for CustomCalibrationLibrary.Views.ValidationResult:

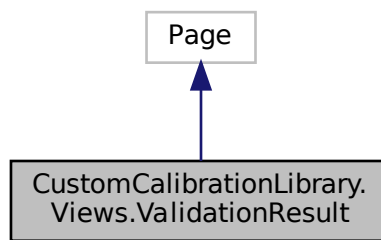

Public Member Functions

- [ValidationResult](#) (CalibrationModel model)
Initializes a new instance of the [ValidationResult](#) class.

8.63.1 Detailed Description

Interaction logic for ValidationResult.xaml

8.63.2 Constructor & Destructor Documentation

8.63.2.1 ValidationResult()

```
CustomCalibrationLibrary.Views.ValidationResult.ValidationResult (
    CalibrationModel model ) [inline]
```

Initializes a new instance of the [ValidationResult](#) class.

Parameters

<i>model</i>	The calibration model.
--------------	------------------------

The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/Views/ValidationResult.xaml.cs

8.64 CustomCalibrationLibrary.ViewModels.ValidationResultViewModel Class Reference

View model class of the gaze validation result.

Public Member Functions

- [ValidationResultViewModel](#) ([CalibrationModel](#) model)

Constructor

Properties

- ICommand [ValidationRestartCommand](#) [get]
Command to restart the validation
- ICommand [ValidationCloseCommand](#) [get]
Command to close the validation window
- [GazeValidationData](#) [ValidationData](#) [get]
The validation result

8.64.1 Detailed Description

View model class of the gaze validation result.

8.64.2 Constructor & Destructor Documentation

8.64.2.1 ValidationResultViewModel()

```
CustomCalibrationLibrary.ViewModels.ValidationResultViewModel.ValidationResultViewModel (
    CalibrationModel model ) [inline]
```

Constructor

Parameters

<i>model</i>	The claibration model
--------------	-----------------------

8.64.3 Property Documentation

8.64.3.1 ValidationCloseCommand

```
 ICommand CustomCalibrationLibrary.ViewModels.ValidationResultViewModel.ValidationCloseCommand  
 [get]
```

Command to close the validation window

8.64.3.2 ValidationData

```
 GazeValidationData CustomCalibrationLibrary.ViewModels.ValidationResultViewModel.Validation↵  
 Data [get]
```

The validation result

8.64.3.3 ValidationRestartCommand

```
 ICommand CustomCalibrationLibrary.ViewModels.ValidationResultViewModel.ValidationRestart↵  
 Command [get]
```

Command to restart the validation

The documentation for this class was generated from the following file:

- source/CustomCalibrationLibrary/ViewModels/ValidationResultViewModel.cs

Index

- [_model](#)
 - [CustomCalibrationLibrary.ViewModels.CalibrationViewModel, 90](#)
- [AccuracyLeft](#)
 - [GazeUtilityLibrary.DataStructs.GazeValidationData, 147](#)
- [AccuracyLeftEye](#)
 - [Tobii.Research.Addons.CalibrationValidationPoint, 84](#)
- [AccuracyRight](#)
 - [GazeUtilityLibrary.DataStructs.GazeValidationData, 148](#)
- [AccuracyRightEye](#)
 - [Tobii.Research.Addons.CalibrationValidationPoint, 85](#)
- [AddPoint](#)
 - [GazeUtilityLibrary.DataStructs.GazeValidationData, 147](#)
- [App](#)
 - [GazeToMouse.App, 35](#)
- [ApplyCalibration](#)
 - [GazeUtilityLibrary.Tracker.BaseTracker, 44](#)
 - [GazeUtilityLibrary.Tracker.EyeTrackerPro, 116](#)
 - [GazeUtilityLibrary.Tracker.MouseTracker, 162](#)
- [AverageAccuracyLeftEye](#)
 - [Tobii.Research.Addons.CalibrationValidationResult, 87](#)
- [AverageAccuracyRightEye](#)
 - [Tobii.Research.Addons.CalibrationValidationResult, 87](#)
- [AveragePrecisionLeftEye](#)
 - [Tobii.Research.Addons.CalibrationValidationResult, 87](#)
- [AveragePrecisionRightEye](#)
 - [Tobii.Research.Addons.CalibrationValidationResult, 87](#)
- [AveragePrecisionRMSLeftEye](#)
 - [Tobii.Research.Addons.CalibrationValidationResult, 88](#)
- [AveragePrecisionRMSRightEye](#)
 - [Tobii.Research.Addons.CalibrationValidationResult, 88](#)
- [BaseTracker](#)
 - [GazeUtilityLibrary.Tracker.BaseTracker, 43](#)
- [BottomLeft](#)
 - [GazeUtilityLibrary.ConfigScreenArea, 103](#)
 - [GazeUtilityLibrary.ScreenArea, 177](#)
- [BottomRight](#)
 - [GazeUtilityLibrary.ConfigScreenArea, 103](#)
 - [GazeUtilityLibrary.ScreenArea, 177](#)
- [Calibration](#)
 - [CustomCalibrationLibrary.Views.Calibration, 56](#)
- [CalibrationAbortCommand](#)
 - [CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel, 185](#)
- [CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel, 197](#)
- [CustomCalibrationLibrary.Views.CalibrationFailed, 62](#)
- [CustomCalibrationLibrary.Views.Disconnect, 106](#)
- [CalibrationAcceptCommand](#)
 - [CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel, 83](#)
- [CalibrationCommand](#)
 - [CustomCalibrationLibrary.Commands.CalibrationCommand, 58](#)
- [CalibrationEvent](#)
 - [CustomCalibrationLibrary.Models.CalibrationModel, 71](#)
- [CalibrationEventType](#)
 - [CustomCalibrationLibrary.Models, 24](#)
- [CalibrationFailed](#)
 - [CustomCalibrationLibrary.Views.CalibrationFailed, 62](#)
- [CalibrationFrame](#)
 - [CustomCalibrationLibrary.Views.CalibrationFrame, 64](#)
- [CalibrationLogColumnOrder](#)
 - [GazeUtilityLibrary.ConfigItem, 95](#)
- [CalibrationLogColumnName](#)
 - [GazeUtilityLibrary.ConfigItem, 95](#)
- [CalibrationLogWriteOutput](#)
 - [GazeUtilityLibrary.ConfigItem, 95](#)
- [CalibrationModel](#)
 - [CustomCalibrationLibrary.Models.CalibrationModel, 67](#)
- [CalibrationOutputValue](#)
 - [GazeUtilityLibrary.DataStructs, 29](#)
- [CalibrationPoint](#)
 - [CustomCalibrationLibrary.Views.CalibrationPoint, 75](#)
 - [GazeUtilityLibrary.DataStructs.CalibrationPoint, 73](#)
- [CalibrationPoints](#)
 - [CustomCalibrationLibrary.Models.CalibrationModel, 69](#)
 - [CustomCalibrationLibrary.ViewModels.CalibrationViewModel, 90](#)

- GazeUtilityLibrary.ConfigItem, 96
- CalibrationPointViewModel
 - CustomCalibrationLibrary.ViewModels.CalibrationPointViewModel, 77
- CalibrationRestartCommand
 - CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel, 83
 - CustomCalibrationLibrary.Views.CalibrationFailed, 63
- CalibrationResult
 - CustomCalibrationLibrary.Views.CalibrationResult, 79
- CalibrationResultPoint
 - CustomCalibrationLibrary.Views.CalibrationResultPoint, 80
- CalibrationResultViewModel
 - CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel, 82
- CalibrationStartCommand
 - CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel, 185
 - CustomCalibrationLibrary.ViewModels.UserPositionGuidanceViewModel, 198
- CalibrationStatus
 - CustomCalibrationLibrary.Models, 24
- CalibrationValidate
 - GazeToMouse.App, 36
- CalibrationViewModel
 - CustomCalibrationLibrary.ViewModels.CalibrationViewModel, 90
- CanExecute
 - CustomCalibrationLibrary.Commands.CalibrationCommand, 58
- CanExecuteChanged
 - CustomCalibrationLibrary.Commands.CalibrationCommand, 59
- Center
 - GazeUtilityLibrary.ConfigScreenArea, 103
 - GazeUtilityLibrary.ScreenArea, 177
- CleanupCalibrationOutputFile
 - GazeUtilityLibrary.GazeConfiguration, 129
- CleanupGazeOutputFile
 - GazeUtilityLibrary.GazeConfiguration, 129
- CleanupValidationOutputFile
 - GazeUtilityLibrary.GazeConfiguration, 129
- CollectCalibrationDataAsync
 - GazeUtilityLibrary.Tracker.BaseTracker, 44
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, 117
 - GazeUtilityLibrary.Tracker.MouseTracker, 162
- CollectValidationDataAsync
 - GazeUtilityLibrary.Tracker.BaseTracker, 44
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, 117
 - GazeUtilityLibrary.Tracker.MouseTracker, 162
- Combined
 - GazeUtilityLibrary.DataStructs.GazeData, 136
- Command
 - GazeControlLibrary.PipeCommand, 167
- CompensateDrift
 - GazeToMouse.App, 36
- Compensation
 - GazeUtilityLibrary.DataStructs.DriftCompensationData, 110
- Compute
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 181
- ComputeValidation
 - GazeUtilityLibrary.Tracker.BaseTracker, 46
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, 117
 - GazeUtilityLibrary.Tracker.MouseTracker, 163
- Computing
 - CustomCalibrationLibrary.Views.Computing, 93
- Config
 - GazeUtilityLibrary.GazeConfiguration, 132
- config
 - GazeUtilityLibrary.Tracker.BaseTracker, 52
- ConfigItem
 - GazeUtilityLibrary.ConfigItem, 95
 - GazeUtilityLibrary.ConfigItem, 96
 - GazeUtilityLibrary.ConfigScreenArea, 102, 103
- Convert
 - CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter, 152
 - CustomCalibrationLibrary.Converters.PositionConverter, 169
 - CustomCalibrationLibrary.Converters.ProximityColorConverter, 172
- ConvertBack
 - CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter, 152
 - CustomCalibrationLibrary.Converters.PositionConverter, 169
 - CustomCalibrationLibrary.Converters.ProximityColorConverter, 172
- ConvertToBinString
 - GazeUtilityLibrary.GazeError, 145
- Coordinates
 - Tobii.Research.Addons.CalibrationValidationPoint, 85
- CursorType
 - CustomCalibrationLibrary.Models.CalibrationModel, 69
- CustomCalibrate
 - GazeToMouse.App, 36
- CustomCalibrationLibrary, 23
- CustomCalibrationLibrary.Commands, 23
- CustomCalibrationLibrary.Commands.CalibrationCommand, 57
 - CalibrationCommand, 58
 - CanExecute, 58
 - CanExecuteChanged, 59
 - Execute, 58
- CustomCalibrationLibrary.Converters, 23
- CustomCalibrationLibrary.Converters.HasDataToVisibilityConverter, 151

- Convert, 152
- ConvertBack, 152
- CustomCalibrationLibrary.Converters.PositionConverter, 168
 - Convert, 169
 - ConvertBack, 169
 - Offset, 170
 - OffsetProperty, 170
- CustomCalibrationLibrary.Converters.ProximityColorConverter, 171
 - Convert, 172
 - ConvertBack, 172
- CustomCalibrationLibrary.Models, 23
 - CalibrationEventType, 24
 - CalibrationStatus, 24
- CustomCalibrationLibrary.Models.CalibrationModel, 65
 - CalibrationEvent, 71
 - CalibrationModel, 67
 - CalibrationPoints, 69
 - CursorType, 69
 - Error, 69
 - GazeDataCollected, 67
 - GazePoint, 69
 - GazePointChanged, 71
 - Index, 70
 - InitCalibration, 68
 - LastStatus, 70
 - NextCalibrationPoint, 68
 - OnCalibrationEvent, 68
 - Points, 70
 - PropertyChanged, 71
 - RedoCalibrationPoint, 68
 - SetCalibrationResult, 68
 - Status, 70
 - UpdateGazePoint, 69
 - UserPositionGuide, 70
 - UserPositionGuideChanged, 71
 - ValidationData, 70
- CustomCalibrationLibrary.ViewModels, 24
- CustomCalibrationLibrary.ViewModels.CalibrationPointViewModel, 76
 - CalibrationPointViewModel, 77
- CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel, 81
 - CalibrationAcceptCommand, 83
 - CalibrationRestartCommand, 83
 - CalibrationResultViewModel, 82
 - GazePoint, 83
 - GazeVisibilityCommand, 83
 - OnGazeToggle, 82
- CustomCalibrationLibrary.ViewModels.CalibrationViewModel, 89
 - _model, 90
 - CalibrationPoints, 90
 - CalibrationViewModel, 90
- CustomCalibrationLibrary.ViewModels.DriftCompensationViewModel, 110
 - DriftCompensationViewModel, 111
- FixationPoint, 111
- CustomCalibrationLibrary.ViewModels.Monitor, 157
 - Index, 157
 - Monitor, 157
 - Name, 158
- CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel, 184
 - CalibrationAbortCommand, 185
 - CalibrationStartCommand, 185
 - Monitors, 185
 - ScreenSelectionViewModel, 185
 - ScreenSwitchCommand, 186
- CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel, 197
 - CalibrationAbortCommand, 197
 - CalibrationStartCommand, 198
 - UserPosition, 198
 - UserPositionGuideViewModel, 197
- CustomCalibrationLibrary.ViewModels.ValidationResultViewModel, 200
 - ValidationCloseCommand, 200
 - ValidationData, 201
 - ValidationRestartCommand, 201
 - ValidationResultViewModel, 200
- CustomCalibrationLibrary.Views, 25
- CustomCalibrationLibrary.Views.Calibration, 55
 - Calibration, 56
- CustomCalibrationLibrary.Views.CalibrationFailed, 61
 - CalibrationAbortCommand, 62
 - CalibrationFailed, 62
 - CalibrationRestartCommand, 63
 - Error, 63
 - PropertyChanged, 63
- CustomCalibrationLibrary.Views.CalibrationFrame, 64
 - CalibrationFrame, 64
- CustomCalibrationLibrary.Views.CalibrationPoint, 75
 - CalibrationPoint, 75
- CustomCalibrationLibrary.Views.CalibrationResult, 78
 - CalibrationResult, 79
- CustomCalibrationLibrary.Views.CalibrationResultPoint, 79
 - CalibrationResultPoint, 80
- CustomCalibrationLibrary.Views.CalibrationWindow, 91
 - CalibrationWindow, 91
- CustomCalibrationLibrary.Views.Computing, 92
 - Computing, 93
- CustomCalibrationLibrary.Views.Disconnect, 105
 - CalibrationAbortCommand, 106
 - Disconnect, 106
- CustomCalibrationLibrary.Views.DriftCompensationWindow, 112
 - DriftCompensationWindow, 112
- CustomCalibrationLibrary.Views.FixationPoint, 121
 - FixationPoint, 122
- CustomCalibrationLibrary.Views.ScreenSelection, 183
 - ScreenSelection, 184
- CustomCalibrationLibrary.Views.UserPositionGuide, 195
 - UserPositionGuide, 196

- CustomCalibrationLibrary.Views.ValidationResult, 198
 - ValidationResult, 199
- DataLogColumnOrder
 - GazeUtilityLibrary.ConfigItem, 96
- DataLogColumnName
 - GazeUtilityLibrary.ConfigItem, 96
- DataLogCount
 - GazeUtilityLibrary.ConfigItem, 96
- DataLogDisabledOnStartup
 - GazeUtilityLibrary.ConfigItem, 96
- DataLogFormatDiameter
 - GazeUtilityLibrary.ConfigItem, 97
- DataLogFormatNormalizedPoint
 - GazeUtilityLibrary.ConfigItem, 97
- DataLogFormatOrigin
 - GazeUtilityLibrary.ConfigItem, 97
- DataLogFormatTimeStamp
 - GazeUtilityLibrary.ConfigItem, 97
- DataLogFormatTimeStampRelative
 - GazeUtilityLibrary.ConfigItem, 97
- DataLogFormatValidation
 - GazeUtilityLibrary.ConfigItem, 97
- DataLogPath
 - GazeUtilityLibrary.ConfigItem, 98
- DataLogWriteOutput
 - GazeUtilityLibrary.ConfigItem, 98
- Debug
 - GazeUtilityLibrary.TrackerLogger, 189
- DeviationAngle
 - GazeUtilityLibrary.DriftCompensation, 108
- DeviceName
 - GazeUtilityLibrary.Tracker.BaseTracker, 53
- DeviceStatus
 - GazeUtilityLibrary.Tracker.BaseTracker, 43
- dialogBoxTimer
 - GazeUtilityLibrary.Tracker.BaseTracker, 53
- DiscardData
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 181
- Disconnect
 - CustomCalibrationLibrary.Views.Disconnect, 106
- Dispersion
 - GazeUtilityLibrary.DriftCompensation, 108
- Dispose
 - GazeUtilityLibrary.Tracker.BaseTracker, 46
 - GazeUtilityLibrary.Tracker.MouseTracker, 163
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 181
- DriftCompensation
 - GazeUtilityLibrary.DataStructs.GazeData, 136
 - GazeUtilityLibrary.DriftCompensation, 107
- driftCompensation
 - GazeUtilityLibrary.Tracker.BaseTracker, 53
- DriftCompensationComputed
 - GazeUtilityLibrary.Tracker.BaseTracker, 54
- DriftCompensationData
 - GazeUtilityLibrary.DataStructs.DriftCompensationData, 109
- DriftCompensationDispersionThreshold
 - GazeUtilityLibrary.ConfigItem, 98
- DriftCompensationDispersionThresholdMax
 - GazeUtilityLibrary.ConfigItem, 98
- DriftCompensationDurationThreshold
 - GazeUtilityLibrary.ConfigItem, 98
- DriftCompensationEventHandler
 - GazeUtilityLibrary.Tracker.BaseTracker, 47
- DriftCompensationTimer
 - GazeUtilityLibrary.ConfigItem, 98
- DriftCompensationViewModel
 - CustomCalibrationLibrary.ViewModels.DriftCompensationViewModel, 111
- DriftCompensationWindow
 - CustomCalibrationLibrary.Views.DriftCompensationWindow, 112
- DriftCompensationWindowShow
 - GazeUtilityLibrary.ConfigItem, 99
- Dump
 - GazeUtilityLibrary.ScreenArea, 174
- DumpCurrentConfigurationFile
 - GazeUtilityLibrary.GazeConfiguration, 130
- DumpFatal
 - GazeUtilityLibrary.TrackerLogger, 190
- E1
 - GazeUtilityLibrary.ScreenTriangle, 188
- E2
 - GazeUtilityLibrary.ScreenTriangle, 188
- ECalibrationDataError
 - GazeUtilityLibrary, 27
- EGazeConfigError
 - GazeUtilityLibrary, 27
- EGazeDataError
 - GazeUtilityLibrary, 28
- EnterValidationMode
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 181
- EOutputType
 - GazeUtilityLibrary, 28
- Error
 - CustomCalibrationLibrary.Models.CalibrationModel, 69
 - CustomCalibrationLibrary.Views.CalibrationFailed, 63
 - GazeUtilityLibrary.CalibrationDataError, 61
 - GazeUtilityLibrary.GazeConfigError, 127
 - GazeUtilityLibrary.GazeDataError, 144
 - GazeUtilityLibrary.TrackerLogger, 190
- Execute
 - CustomCalibrationLibrary.Commands.CalibrationCommand, 58
- EyeData
 - GazeUtilityLibrary.DataStructs.EyeData, 113
 - GazeUtilityLibrary.DataStructs.GazeDataCollection, 142
- EyeTrackerPro
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, 116

- FinishCalibration
 - GazeUtilityLibrary.Tracker.BaseTracker, [47](#)
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, [118](#)
 - GazeUtilityLibrary.Tracker.MouseTracker, [163](#)
- FinishCalibrationAsync
 - GazeUtilityLibrary.Tracker.BaseTracker, [47](#)
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, [118](#)
 - GazeUtilityLibrary.Tracker.MouseTracker, [163](#)
- FinishValidation
 - GazeUtilityLibrary.Tracker.BaseTracker, [47](#)
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, [118](#)
 - GazeUtilityLibrary.Tracker.MouseTracker, [164](#)
- FixationPoint
 - CustomCalibrationLibrary.ViewModels.DriftCompensationViewModel, [111](#)
 - CustomCalibrationLibrary.Views.FixationPoint, [122](#)
- GazeCalibrationData
 - GazeUtilityLibrary.DataStructs.GazeCalibrationData, [123](#)
- GazeConfiguration
 - GazeUtilityLibrary.GazeConfiguration, [128](#)
- GazeControl, [25](#)
- GazeControl.App, [33](#)
- GazeControlLibrary, [25](#)
 - LogLevel, [26](#)
- GazeControlLibrary.PipeCommand, [166](#)
 - Command, [167](#)
 - Label, [167](#)
 - PipeCommand, [166](#)
 - ResetStartTime, [167](#)
 - TrialId, [167](#)
- GazeData
 - GazeUtilityLibrary.DataStructs.GazeData, [133](#), [134](#)
 - Tobii.Research.Addons.CalibrationValidationPoint, [85](#)
- GazeData2d
 - GazeUtilityLibrary.DataStructs.GazeData2d, [137](#)
 - GazeUtilityLibrary.DataStructs.GazeDataCollection, [142](#)
- GazeData3d
 - GazeUtilityLibrary.DataStructs.GazeData3d, [139](#)
 - GazeUtilityLibrary.DataStructs.GazeDataCollection, [142](#)
- GazeDataCollected
 - CustomCalibrationLibrary.Models.CalibrationModel, [67](#)
- GazeDataCollection
 - GazeUtilityLibrary.DataStructs.GazeDataCollection, [141](#)
- GazeDataHandler
 - GazeUtilityLibrary.Tracker.BaseTracker, [47](#)
- GazeDataReceived
 - GazeUtilityLibrary.Tracker.BaseTracker, [54](#)
- GazeDirection
 - GazeUtilityLibrary.DataStructs.GazeData3d, [139](#)
- GazeDistance
 - GazeUtilityLibrary.DataStructs.GazeData3d, [139](#)
- GazeOrigin
 - GazeUtilityLibrary.DataStructs.GazeData3d, [139](#)
- GazeOutputValue
 - GazeUtilityLibrary.DataStructs, [29](#)
- GazePoint
 - CustomCalibrationLibrary.Models.CalibrationModel, [69](#)
 - CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel, [83](#)
 - GazeUtilityLibrary.DataStructs.GazeData2d, [138](#)
 - GazeUtilityLibrary.DataStructs.GazeData3d, [140](#)
- GazePointChanged
 - CustomCalibrationLibrary.Models.CalibrationModel, [71](#)
- GazePosition
 - GazeUtilityLibrary.DataStructs.DriftCompensationData, [110](#)
- GazePosition3d
 - GazeUtilityLibrary.DataStructs.DriftCompensationData, [110](#)
- GazePositionAverage
 - GazeUtilityLibrary.DataStructs.CalibrationPoint, [73](#)
- GazePositionLeft
 - GazeUtilityLibrary.DataStructs.CalibrationPoint, [73](#)
- GazePositionRight
 - GazeUtilityLibrary.DataStructs.CalibrationPoint, [74](#)
- GazeRecordingDisable
 - GazeToMouse.App, [36](#)
- GazeRecordingEnable
 - GazeToMouse.App, [36](#)
- GazeToMouse, [26](#)
- GazeToMouse.App, [34](#)
 - App, [35](#)
 - CalibrationValidate, [36](#)
 - CompensateDrift, [36](#)
 - CustomCalibrate, [36](#)
 - GazeRecordingDisable, [36](#)
 - GazeRecordingEnable, [36](#)
 - LastTag, [37](#)
 - MouseTrackingDisable, [37](#)
 - MouseTrackingEnable, [37](#)
 - ResetDriftCompensation, [37](#)
 - StartTime, [37](#)
 - Tag, [37](#)
 - TrialId, [38](#)
- GazeUtilityLibrary, [26](#)
 - ECalibrationDataError, [27](#)
 - EGazeConfigError, [27](#)
 - EGazeDataError, [28](#)
 - EOutputType, [28](#)
- GazeUtilityLibrary.CalibrationDataError, [59](#)
 - Error, [61](#)
 - GetCalibrationDataErrorString, [60](#)
- GazeUtilityLibrary.ConfigItem, [93](#)
 - CalibrationLogColumnOrder, [95](#)
 - CalibrationLogColumnName, [95](#)
 - CalibrationLogWriteOutput, [95](#)
 - CalibrationPoints, [96](#)
 - ConfigItem, [95](#)

- ConfigName, 96
- DataLogColumnOrder, 96
- DataLogColumnTitle, 96
- DataLogCount, 96
- DataLogDisabledOnStartup, 96
- DataLogFormatDiameter, 97
- DataLogFormatNormalizedPoint, 97
- DataLogFormatOrigin, 97
- DataLogFormatTimeStamp, 97
- DataLogFormatTimeStampRelative, 97
- DataLogFormatValidation, 97
- DataLogPath, 98
- DataLogWriteOutput, 98
- DriftCompensationDispersionThreshold, 98
- DriftCompensationDispersionThresholdMax, 98
- DriftCompensationDurationThreshold, 98
- DriftCompensationTimer, 98
- DriftCompensationWindowShow, 99
- LicensePath, 99
- MouseCalibrationHide, 99
- MouseControl, 99
- MouseControlHide, 99
- MouseStandardIconPath, 99
- ReadyTimer, 100
- ScreenArea, 100
- TobiiApplicationPath, 100
- TobiiCalibrate, 100
- TobiiCalibrateArguments, 100
- TrackerDevice, 100
- ValidationDurationThreshold, 101
- ValidationLogColumnOrder, 101
- ValidationLogColumnTitle, 101
- ValidationLogWriteOutput, 101
- ValidationPoints, 101
- ValidationTimer, 101
- GazeUtilityLibrary.ConfigScreenArea, 102
 - BottomLeft, 103
 - BottomRight, 103
 - Center, 103
 - ConfigScreenArea, 102, 103
 - Height, 103
 - TopLeft, 104
 - TopRight, 104
 - Width, 104
- GazeUtilityLibrary.DataStructs, 28
 - CalibrationOutputValue, 29
 - GazeOutputValue, 29
 - ValidationOutputValue, 29
- GazeUtilityLibrary.DataStructs.CalibrationPoint, 72
 - CalibrationPoint, 73
 - GazePositionAverage, 73
 - GazePositionLeft, 73
 - GazePositionRight, 74
 - HasData, 74
 - Index, 74
 - Position, 74
 - PropertyChanged, 74
- GazeUtilityLibrary.DataStructs.DriftCompensationData, 109
 - Compensation, 110
 - DriftCompensationData, 109
 - GazePosition2d, 110
 - GazePosition3d, 110
- GazeUtilityLibrary.DataStructs.EyeData, 113
 - EyeData, 113
 - IsPupilDiameterValid, 114
 - PupilDiameter, 114
- GazeUtilityLibrary.DataStructs.GazeCalibrationData, 122
 - GazeCalibrationData, 123
 - Prepare, 124
 - ValidityLeft, 124
 - ValidityRight, 124
 - XCoord, 124
 - XCoordLeft, 124
 - XCoordRight, 125
 - YCoord, 125
 - YCoordLeft, 125
 - YCoordRight, 125
- GazeUtilityLibrary.DataStructs.GazeData, 133
 - Combined, 136
 - DriftCompensation, 136
 - GazeData, 133, 134
 - Left, 136
 - Prepare, 135
 - Right, 136
 - Timestamp, 136
 - TimestampReceived, 136
- GazeUtilityLibrary.DataStructs.GazeData2d, 137
 - GazeData2d, 137
 - GazePoint, 138
 - IsGazePointValid, 138
- GazeUtilityLibrary.DataStructs.GazeData3d, 138
 - GazeData3d, 139
 - GazeDirection, 139
 - GazeDistance, 139
 - GazeOrigin, 139
 - GazePoint, 140
 - IsGazeOriginValid, 140
 - IsGazePointValid, 140
- GazeUtilityLibrary.DataStructs.GazeDataCollection, 140
 - EyeData, 142
 - GazeData2d, 142
 - GazeData3d, 142
 - GazeDataCollection, 141
- GazeUtilityLibrary.DataStructs.GazeValidationData, 145
 - AccuracyLeft, 147
 - AccuracyRight, 148
 - AddPoint, 147
 - GazeValidationData, 146
 - Points, 148
 - PrecisionLeft, 148
 - PrecisionRight, 148
 - PrecisionRmsLeft, 148
 - PrecisionRmsRight, 148

- GazeUtilityLibrary.DataStructs.GazeValidationPoint, 149
 - GazeValidationPoint, 149
 - Point, 150
 - Prepare, 150
 - Result, 150
- GazeUtilityLibrary.DataStructs.LiveGazePoint, 155
 - PropertyChanged, 156
 - Visibility, 156
 - X, 156
 - Y, 156
- GazeUtilityLibrary.DataStructs.UserPositionData, 192
 - PropertyChanged, 195
 - UserPositionData, 193
 - XCoordLeft, 194
 - XCoordRight, 194
 - YCoordLeft, 194
 - YCoordRight, 194
 - ZCoordLeft, 194
 - ZCoordRight, 195
- GazeUtilityLibrary.DriftCompensation, 106
 - DeviationAngle, 108
 - Dispersion, 108
 - DriftCompensation, 107
 - Q, 108
 - Reset, 107
 - Start, 107
 - Update, 108
- GazeUtilityLibrary.GazeConfigError, 126
 - Error, 127
 - GetGazeConfigErrorString, 127
- GazeUtilityLibrary.GazeConfiguration, 127
 - CleanupCalibrationOutputFile, 129
 - CleanupGazeOutputFile, 129
 - CleanupValidationOutputFile, 129
 - Config, 132
 - DumpCurrentConfigurationFile, 130
 - GazeConfiguration, 128
 - InitConfig, 130
 - PrepareCalibrationOutputFile, 130
 - PrepareGazeOutputFile, 131
 - PrepareValidationOutputFile, 131
 - WriteToCalibrationOutput, 131
 - WriteToGazeOutput, 132
 - WriteToValidationOutput, 132
- GazeUtilityLibrary.GazeDataError, 143
 - Error, 144
 - GetGazeDataErrorString, 144
- GazeUtilityLibrary.GazeError, 144
 - ConvertToBinString, 145
- GazeUtilityLibrary.JsonConfigParser, 153
 - GetDefaultConfig, 154
 - JsonConfigParser, 153
 - ParseJsonConfig, 154
 - SerializeJsonConfig, 154
- GazeUtilityLibrary.MouseHider, 158
 - HideCursor, 159
 - MouseHider, 158
 - ShowCursor, 159
- GazeUtilityLibrary.ScreenArea, 173
 - BottomLeft, 177
 - BottomRight, 177
 - Center, 177
 - Dump, 174
 - GetIntersectionPoint, 175
 - GetPoint2d, 175
 - GetPoint2dNormalized, 175
 - Height, 177
 - ScreenArea, 174
 - TopLeft, 177
 - TopRight, 178
 - Width, 178
- GazeUtilityLibrary.ScreenTriangle, 186
 - E1, 188
 - E2, 188
 - GetIntersectionPoint, 187
 - ScreenTriangle, 187
 - V1, 188
 - V2, 188
 - V3, 188
- GazeUtilityLibrary.Tracker, 30
- GazeUtilityLibrary.Tracker.BaseTracker, 40
 - ApplyCalibration, 44
 - BaseTracker, 43
 - CollectCalibrationDataAsync, 44
 - CollectValidationDataAsync, 44
 - ComputeValidation, 46
 - config, 52
 - DeviceName, 53
 - DeviceStatus, 43
 - dialogBoxTimer, 53
 - Dispose, 46
 - driftCompensation, 53
 - DriftCompensationComputed, 54
 - DriftCompensationEventHandler, 47
 - FinishCalibration, 47
 - FinishCalibrationAsync, 47
 - FinishValidation, 47
 - GazeDataHandler, 47
 - GazeDataReceived, 54
 - GetFixationFrameCount, 48
 - GetUnitDirection, 48
 - InitCalibration, 48
 - InitCalibrationAsync, 48
 - InitDriftCompensation, 49
 - InitValidation, 49
 - IsInitialised, 49
 - IsReady, 49
 - logger, 53
 - OnGazeDataReceived, 50
 - OnPropertyChanged, 50
 - OnTrackerDisabled, 50
 - OnTrackerDisabledTimeout, 51
 - OnTrackerEnabled, 51
 - OnUserPositionDataReceived, 51
 - PatternReplace, 51
 - PropertyChanged, 54

- ResetDriftCompensation, [52](#)
- ScreenArea, [54](#)
- screenArea, [53](#)
- StartDriftCompensation, [52](#)
- State, [54](#)
- TrackerDisabled, [54](#)
- TrackerEnabled, [55](#)
- trackerMessageBox, [53](#)
- UserPositionDataHandler, [52](#)
- UserPositionDataReceived, [55](#)
- GazeUtilityLibrary.Tracker.EyeTrackerPro, [114](#)
 - ApplyCalibration, [116](#)
 - CollectCalibrationDataAsync, [117](#)
 - CollectValidationDataAsync, [117](#)
 - ComputeValidation, [117](#)
 - EyeTrackerPro, [116](#)
 - FinishCalibration, [118](#)
 - FinishCalibrationAsync, [118](#)
 - FinishValidation, [118](#)
 - GetFixationFrameCount, [118](#)
 - GetUnitDirection, [119](#)
 - InitCalibration, [119](#)
 - InitCalibrationAsync, [119](#)
 - InitDriftCompensation, [119](#)
 - InitValidation, [120](#)
 - IsInitialised, [120](#)
 - IsLicenseOk, [120](#)
 - PatternReplace, [120](#)
- GazeUtilityLibrary.Tracker.MouseTracker, [159](#)
 - ApplyCalibration, [162](#)
 - CollectCalibrationDataAsync, [162](#)
 - CollectValidationDataAsync, [162](#)
 - ComputeValidation, [163](#)
 - Dispose, [163](#)
 - FinishCalibration, [163](#)
 - FinishCalibrationAsync, [163](#)
 - FinishValidation, [164](#)
 - GetFixationFrameCount, [164](#)
 - GetUnitDirection, [164](#)
 - InitCalibration, [165](#)
 - InitCalibrationAsync, [165](#)
 - InitDriftCompensation, [165](#)
 - InitValidation, [165](#)
 - MouseTracker, [161](#)
 - Start, [165](#)
 - Stop, [166](#)
- GazeUtilityLibrary.TrackerLogger, [189](#)
 - Debug, [189](#)
 - DumpFatal, [190](#)
 - Error, [190](#)
 - Info, [190](#)
 - TrackerLogger, [189](#)
 - Warning, [190](#)
- GazeUtilityLibrary.TrackerMessageBox, [191](#)
- GazeValidationData
 - GazeUtilityLibrary.DataStructs.GazeValidationData, [146](#)
- GazeValidationPoint
 - GazeUtilityLibrary.DataStructs.GazeValidationPoint, [149](#)
- GazeVisibilityCommand
 - CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel, [83](#)
- GetCalibrationDataErrorString
 - GazeUtilityLibrary.CalibrationDataError, [60](#)
- GetDefaultConfig
 - GazeUtilityLibrary.JsonConfigParser, [154](#)
- GetFixationFrameCount
 - GazeUtilityLibrary.Tracker.BaseTracker, [48](#)
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, [118](#)
 - GazeUtilityLibrary.Tracker.MouseTracker, [164](#)
- GetGazeConfigErrorString
 - GazeUtilityLibrary.GazeConfigError, [127](#)
- GetGazeDataErrorString
 - GazeUtilityLibrary.GazeDataError, [144](#)
- GetIntersectionPoint
 - GazeUtilityLibrary.ScreenArea, [175](#)
 - GazeUtilityLibrary.ScreenTriangle, [187](#)
- GetPoint2d
 - GazeUtilityLibrary.ScreenArea, [175](#)
- GetPoint2dNormalized
 - GazeUtilityLibrary.ScreenArea, [175](#)
- GetUnitDirection
 - GazeUtilityLibrary.Tracker.BaseTracker, [48](#)
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, [119](#)
 - GazeUtilityLibrary.Tracker.MouseTracker, [164](#)
- HasData
 - GazeUtilityLibrary.DataStructs.CalibrationPoint, [74](#)
- Height
 - GazeUtilityLibrary.ConfigScreenArea, [103](#)
 - GazeUtilityLibrary.ScreenArea, [177](#)
- HideCursor
 - GazeUtilityLibrary.MouseHider, [159](#)
- Index
 - CustomCalibrationLibrary.Models.CalibrationModel, [70](#)
 - CustomCalibrationLibrary.ViewModels.Monitor, [157](#)
 - GazeUtilityLibrary.DataStructs.CalibrationPoint, [74](#)
- Info
 - GazeUtilityLibrary.TrackerLogger, [190](#)
- InitCalibration
 - CustomCalibrationLibrary.Models.CalibrationModel, [68](#)
 - GazeUtilityLibrary.Tracker.BaseTracker, [48](#)
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, [119](#)
 - GazeUtilityLibrary.Tracker.MouseTracker, [165](#)
- InitCalibrationAsync
 - GazeUtilityLibrary.Tracker.BaseTracker, [48](#)
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, [119](#)
 - GazeUtilityLibrary.Tracker.MouseTracker, [165](#)
- InitConfig
 - GazeUtilityLibrary.GazeConfiguration, [130](#)
- InitDriftCompensation
 - GazeUtilityLibrary.Tracker.BaseTracker, [49](#)

- GazeUtilityLibrary.Tracker.EyeTrackerPro, 119
- GazeUtilityLibrary.Tracker.MouseTracker, 165
- InitValidation
 - GazeUtilityLibrary.Tracker.BaseTracker, 49
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, 120
 - GazeUtilityLibrary.Tracker.MouseTracker, 165
- IsGazeOriginValid
 - GazeUtilityLibrary.DataStructs.GazeData3d, 140
- IsGazePointValid
 - GazeUtilityLibrary.DataStructs.GazeData2d, 138
 - GazeUtilityLibrary.DataStructs.GazeData3d, 140
- IsInitialised
 - GazeUtilityLibrary.Tracker.BaseTracker, 49
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, 120
- IsLicenseOk
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, 120
- IsPupilDiameterValid
 - GazeUtilityLibrary.DataStructs.EyeData, 114
- IsReady
 - GazeUtilityLibrary.Tracker.BaseTracker, 49
- JsonConfigParser
 - GazeUtilityLibrary.JsonConfigParser, 153
- Label
 - GazeControlLibrary.PipeCommand, 167
- LastStatus
 - CustomCalibrationLibrary.Models.CalibrationModel, 70
- LastTag
 - GazeToMouse.App, 37
- LeaveValidationMode
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 181
- Left
 - GazeUtilityLibrary.DataStructs.GazeData, 136
- LicensePath
 - GazeUtilityLibrary.ConfigItem, 99
- logger
 - GazeUtilityLibrary.Tracker.BaseTracker, 53
- LogLevel
 - GazeControlLibrary, 26
- Monitor
 - CustomCalibrationLibrary.ViewModels.Monitor, 157
- Monitors
 - CustomCalibrationLibrary.ViewModels.ScreenSelectionView, 185
- MouseCalibrationHide
 - GazeUtilityLibrary.ConfigItem, 99
- MouseControl
 - GazeUtilityLibrary.ConfigItem, 99
- MouseControlHide
 - GazeUtilityLibrary.ConfigItem, 99
- MouseHider
 - GazeUtilityLibrary.MouseHider, 158
- MouseStandardIconPath
 - GazeUtilityLibrary.ConfigItem, 99
- MouseTracker
 - GazeUtilityLibrary.Tracker.MouseTracker, 161
- MouseTrackingDisable
 - GazeToMouse.App, 37
- MouseTrackingEnable
 - GazeToMouse.App, 37
- Name
 - CustomCalibrationLibrary.ViewModels.Monitor, 158
- NextCalibrationPoint
 - CustomCalibrationLibrary.Models.CalibrationModel, 68
- Offset
 - CustomCalibrationLibrary.Converters.PositionConverter, 170
- OffsetProperty
 - CustomCalibrationLibrary.Converters.PositionConverter, 170
- OnCalibrationEvent
 - CustomCalibrationLibrary.Models.CalibrationModel, 68
- OnGazeDataReceived
 - GazeUtilityLibrary.Tracker.BaseTracker, 50
- OnGazeToggle
 - CustomCalibrationLibrary.ViewModels.CalibrationResultViewModel, 82
- OnPropertyChanged
 - GazeUtilityLibrary.Tracker.BaseTracker, 50
- OnTrackerDisabled
 - GazeUtilityLibrary.Tracker.BaseTracker, 50
- OnTrackerDisabledTimeout
 - GazeUtilityLibrary.Tracker.BaseTracker, 51
- OnTrackerEnabled
 - GazeUtilityLibrary.Tracker.BaseTracker, 51
- OnUserPositionDataReceived
 - GazeUtilityLibrary.Tracker.BaseTracker, 51
- ParseJsonConfig
 - GazeUtilityLibrary.JsonConfigParser, 154
- PatternReplace
 - GazeUtilityLibrary.Tracker.BaseTracker, 51
 - GazeUtilityLibrary.Tracker.EyeTrackerPro, 120
- PipeCommand
 - GazeControlLibrary.PipeCommand, 166
- Point
 - GazeUtilityLibrary.DataStructs.GazeValidationPoint, 150
- Points
 - CustomCalibrationLibrary.Models.CalibrationModel, 70
 - GazeUtilityLibrary.DataStructs.GazeValidationData, 148
 - Tobii.Research.Addons.CalibrationValidationResult, 88
- Position
 - GazeUtilityLibrary.DataStructs.CalibrationPoint, 74
- PrecisionLeft

- GazeUtilityLibrary.DataStructs.GazeValidationData, 148
- PrecisionLeftEye
 - Tobii.Research.Addons.CalibrationValidationPoint, 85
- PrecisionRight
 - GazeUtilityLibrary.DataStructs.GazeValidationData, 148
- PrecisionRightEye
 - Tobii.Research.Addons.CalibrationValidationPoint, 85
- PrecisionRmsLeft
 - GazeUtilityLibrary.DataStructs.GazeValidationData, 148
- PrecisionRMSLeftEye
 - Tobii.Research.Addons.CalibrationValidationPoint, 85
- PrecisionRmsRight
 - GazeUtilityLibrary.DataStructs.GazeValidationData, 148
- PrecisionRMSRightEye
 - Tobii.Research.Addons.CalibrationValidationPoint, 86
- Prepare
 - GazeUtilityLibrary.DataStructs.GazeCalibrationData, 124
 - GazeUtilityLibrary.DataStructs.GazeData, 135
 - GazeUtilityLibrary.DataStructs.GazeValidationPoint, 150
- PrepareCalibrationOutputFile
 - GazeUtilityLibrary.GazeConfiguration, 130
- PrepareGazeOutputFile
 - GazeUtilityLibrary.GazeConfiguration, 131
- PrepareValidationOutputFile
 - GazeUtilityLibrary.GazeConfiguration, 131
- PropertyChanged
 - CustomCalibrationLibrary.Models.CalibrationModel, 71
 - CustomCalibrationLibrary.Views.CalibrationFailed, 63
 - GazeUtilityLibrary.DataStructs.CalibrationPoint, 74
 - GazeUtilityLibrary.DataStructs.LiveGazePoint, 156
 - GazeUtilityLibrary.DataStructs.UserPositionData, 195
 - GazeUtilityLibrary.Tracker.BaseTracker, 54
- PupilDiameter
 - GazeUtilityLibrary.DataStructs.EyeData, 114
- Q
 - GazeUtilityLibrary.DriftCompensation, 108
- ReadyTimer
 - GazeUtilityLibrary.ConfigItem, 100
- RedoCalibrationPoint
 - CustomCalibrationLibrary.Models.CalibrationModel, 68
- Reset
 - GazeUtilityLibrary.DriftCompensation, 107
- ResetDriftCompensation
 - GazeToMouse.App, 37
 - GazeUtilityLibrary.Tracker.BaseTracker, 52
- ResetStartTime
 - GazeControlLibrary.PipeCommand, 167
- Result
 - GazeUtilityLibrary.DataStructs.GazeValidationPoint, 150
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 182
- Right
 - GazeUtilityLibrary.DataStructs.GazeData, 136
- ScreenArea
 - GazeUtilityLibrary.ConfigItem, 100
 - GazeUtilityLibrary.ScreenArea, 174
 - GazeUtilityLibrary.Tracker.BaseTracker, 54
- screenArea
 - GazeUtilityLibrary.Tracker.BaseTracker, 53
- ScreenBasedCalibrationValidation
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 180
- ScreenSelection
 - CustomCalibrationLibrary.Views.ScreenSelection, 184
- ScreenSelectionViewModel
 - CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel, 185
- ScreenSwitchCommand
 - CustomCalibrationLibrary.ViewModels.ScreenSelectionViewModel, 186
- ScreenTriangle
 - GazeUtilityLibrary.ScreenTriangle, 187
- SerializeJsonConfig
 - GazeUtilityLibrary.JsonConfigParser, 154
- SetCalibrationResult
 - CustomCalibrationLibrary.Models.CalibrationModel, 68
- ShowCursor
 - GazeUtilityLibrary.MouseHider, 159
- ShowMouse, 30
- ShowMouse.App, 38
- Start
 - GazeUtilityLibrary.DriftCompensation, 107
 - GazeUtilityLibrary.Tracker.MouseTracker, 165
- StartCollectingData
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 182
- StartDriftCompensation
 - GazeUtilityLibrary.Tracker.BaseTracker, 52
- StartTime
 - GazeToMouse.App, 37
- State
 - GazeUtilityLibrary.Tracker.BaseTracker, 54
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 182
- Status
 - CustomCalibrationLibrary.Models.CalibrationModel, 70
- Stop

- GazeUtilityLibrary.Tracker.MouseTracker, 166
- Tag
 - GazeToMouse.App, 37
- TimedOut
 - Tobii.Research.Addons.CalibrationValidationPoint, 86
- Timestamp
 - GazeUtilityLibrary.DataStructs.GazeData, 136
- TimestampReceived
 - GazeUtilityLibrary.DataStructs.GazeData, 136
- Tobii, 30
- Tobii.Research, 30
- Tobii.Research.Addons, 30
- Tobii.Research.Addons.CalibrationValidationPoint, 83
 - AccuracyLeftEye, 84
 - AccuracyRightEye, 85
 - Coordinates, 85
 - GazeData, 85
 - PrecisionLeftEye, 85
 - PrecisionRightEye, 85
 - PrecisionRMSLeftEye, 85
 - PrecisionRMSRightEye, 86
 - TimedOut, 86
 - ToString, 84
- Tobii.Research.Addons.CalibrationValidationResult, 86
 - AverageAccuracyLeftEye, 87
 - AverageAccuracyRightEye, 87
 - AveragePrecisionLeftEye, 87
 - AveragePrecisionRightEye, 87
 - AveragePrecisionRMSLeftEye, 88
 - AveragePrecisionRMSRightEye, 88
 - Points, 88
 - ToString, 87
- Tobii.Research.Addons.ScreenBasedCalibrationValidation, 178
 - Compute, 181
 - DiscardData, 181
 - Dispose, 181
 - EnterValidationMode, 181
 - LeaveValidationMode, 181
 - Result, 182
 - ScreenBasedCalibrationValidation, 180
 - StartCollectingData, 182
 - State, 182
 - ToString, 182
 - ValidationState, 180
- Tobii.Research.Addons.Utility, 31
- TobiiApplicationPath
 - GazeUtilityLibrary.ConfigItem, 100
- TobiiCalibrate, 31
 - GazeUtilityLibrary.ConfigItem, 100
- TobiiCalibrate.App, 39
- TobiiCalibrateArguments
 - GazeUtilityLibrary.ConfigItem, 100
- TopLeft
 - GazeUtilityLibrary.ConfigScreenArea, 104
 - GazeUtilityLibrary.ScreenArea, 177
- TopRight
 - GazeUtilityLibrary.ConfigScreenArea, 104
 - GazeUtilityLibrary.ScreenArea, 178
- ToString
 - Tobii.Research.Addons.CalibrationValidationPoint, 84
 - Tobii.Research.Addons.CalibrationValidationResult, 87
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 182
- TrackerDevice
 - GazeUtilityLibrary.ConfigItem, 100
- TrackerDisabled
 - GazeUtilityLibrary.Tracker.BaseTracker, 54
- TrackerEnabled
 - GazeUtilityLibrary.Tracker.BaseTracker, 55
- TrackerLogger
 - GazeUtilityLibrary.TrackerLogger, 189
- trackerMessageBox
 - GazeUtilityLibrary.Tracker.BaseTracker, 53
- TrialId
 - GazeControlLibrary.PipeCommand, 167
 - GazeToMouse.App, 38
- Update
 - GazeUtilityLibrary.DriftCompensation, 108
- UpdateGazePoint
 - CustomCalibrationLibrary.Models.CalibrationModel, 69
- UserPosition
 - CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel, 198
- UserPositionData
 - GazeUtilityLibrary.DataStructs.UserPositionData, 193
- UserPositionDataHandler
 - GazeUtilityLibrary.Tracker.BaseTracker, 52
- UserPositionDataReceived
 - GazeUtilityLibrary.Tracker.BaseTracker, 55
- UserPositionGuide
 - CustomCalibrationLibrary.Models.CalibrationModel, 70
 - CustomCalibrationLibrary.Views.UserPositionGuide, 196
- UserPositionGuideChanged
 - CustomCalibrationLibrary.Models.CalibrationModel, 71
- UserPositionGuideViewModel
 - CustomCalibrationLibrary.ViewModels.UserPositionGuideViewModel, 197
- V1
 - GazeUtilityLibrary.ScreenTriangle, 188
- V2
 - GazeUtilityLibrary.ScreenTriangle, 188
- V3
 - GazeUtilityLibrary.ScreenTriangle, 188
- ValidationCloseCommand
 - CustomCalibrationLibrary.ViewModels.ValidationResultViewModel, 200

- ValidationData
 - CustomCalibrationLibrary.Models.CalibrationModel, 70
 - CustomCalibrationLibrary.ViewModels.ValidationResultViewModel, 201
- ValidationDurationThreshold
 - GazeUtilityLibrary.ConfigItem, 101
- ValidationLogColumnOrder
 - GazeUtilityLibrary.ConfigItem, 101
- ValidationLogColumnName
 - GazeUtilityLibrary.ConfigItem, 101
- ValidationLogWriteOutput
 - GazeUtilityLibrary.ConfigItem, 101
- ValidationOutputValue
 - GazeUtilityLibrary.DataStructs, 29
- ValidationPoints
 - GazeUtilityLibrary.ConfigItem, 101
- ValidationRestartCommand
 - CustomCalibrationLibrary.ViewModels.ValidationResultViewModel, 201
- ValidationResult
 - CustomCalibrationLibrary.Views.ValidationResult, 199
- ValidationResultViewModel
 - CustomCalibrationLibrary.ViewModels.ValidationResultViewModel, 200
- ValidationState
 - Tobii.Research.Addons.ScreenBasedCalibrationValidation, 180
- ValidationTimer
 - GazeUtilityLibrary.ConfigItem, 101
- ValidityLeft
 - GazeUtilityLibrary.DataStructs.GazeCalibrationData, 124
- ValidityRight
 - GazeUtilityLibrary.DataStructs.GazeCalibrationData, 124
- Visibility
 - GazeUtilityLibrary.DataStructs.LiveGazePoint, 156
- Warning
 - GazeUtilityLibrary.TrackerLogger, 190
- Width
 - GazeUtilityLibrary.ConfigScreenArea, 104
 - GazeUtilityLibrary.ScreenArea, 178
- WriteToCalibrationOutput
 - GazeUtilityLibrary.GazeConfiguration, 131
- WriteToGazeOutput
 - GazeUtilityLibrary.GazeConfiguration, 132
- WriteToValidationOutput
 - GazeUtilityLibrary.GazeConfiguration, 132
- X
 - GazeUtilityLibrary.DataStructs.LiveGazePoint, 156
- XCoord
 - GazeUtilityLibrary.DataStructs.GazeCalibrationData, 124
- XCoordLeft
 - GazeUtilityLibrary.DataStructs.GazeCalibrationData, 124
 - GazeUtilityLibrary.DataStructs.UserPositionData, 194
- XCoordRight
 - GazeUtilityLibrary.DataStructs.GazeCalibrationData, 125
 - GazeUtilityLibrary.DataStructs.UserPositionData, 194
- Y
 - GazeUtilityLibrary.DataStructs.LiveGazePoint, 156
- YCoord
 - GazeUtilityLibrary.DataStructs.GazeCalibrationData, 125
- YCoordLeft
 - GazeUtilityLibrary.DataStructs.GazeCalibrationData, 125
 - GazeUtilityLibrary.DataStructs.UserPositionData, 194
- YCoordRight
 - GazeUtilityLibrary.DataStructs.GazeCalibrationData, 125
 - GazeUtilityLibrary.DataStructs.UserPositionData, 194
- ZCoordLeft
 - GazeUtilityLibrary.DataStructs.UserPositionData, 194
- ZCoordRight
 - GazeUtilityLibrary.DataStructs.UserPositionData, 195